



# **Technical Publications**

**Direction DOC0660253**

**Revision 15**

## **NM General Purpose 600/800 Series \* DICOM CONFORMANCE STATEMENT**

\* Refer to Section 1 for a list of products in NM General Purpose 600/800 Series family to which this Conformance Statement applies.

**Copyright© 2018 by General Electric Co.**

**Do not duplicate**



## LIST OF REVISIONS

REV	DATE	DESCRIPTION	PAGES	APPR.
1	May 2010	Discovery NM/CT 670 Release	All	T. Cohen
2	Mar 2011	Discovery NM 630 Release Update Private Dictionary and content of NM Images created	1.1 2.1 2.2.1 2.7 3.3 3.4 3.5 4.3 4.4 4.5 6.1	M. Mesh
3	Feb 2012	Brivo NM 615 and Optima NM/CT 640 Release	Conformance statement overview 2.1, 2.2.1, 2.3.1, 2.3.1.2.1.2 3.1, 4.1, 5.1 Added new sections: Section 8 Section 9	M. Mesh
4	May 2012	Update after external review	Conformance statement overview 8.1 8.3 8.5.4 9.7.1.4, 9.7.1.5, 9.7.1.6, 9.7.1.7	M. Mesh
5	July 2012	Modify list of attributes , update attributes descriptions	3.4.1.1 3.4.2.1 3.4.2.4 3.4.3.1 3.4.3.4	M. Mesh



			3.4.6.7 3.4.6.9 3.4.6.15 3.4.6.16 3.5.2 4.4.1.1 4.4.2.1 4.4.2.2 – removed 4.4.5.5 4.5.1 4.5.3 8.4.3.1 8.4.3.5 8.4.3.6 8.4.3.9 8.5.3 5.2.1.2 9.7.1.5.1	
6	April 2013	Update CT Radiation Dose report IOD attributes            MPPS Support            New Private Tag	9.3 Table 9-2 9.4.1.1 9.4.1.2 9.4.1.3 9.4.3.3 9.4.3.4 9.4.3.5 9.5.1 9.5.2 9.5.3 Conformance Statement Overview 1.1 2.1 2.2.2. 2.2.3 2.2.4 2.3.1 2.3.1.2 2.3.1.2.5 – new 2.3.1.1.4 7 3.4.3.4 Table 3-11 3.5.2 Table 3-35	M. Mesh
7	June 2014	Discovery NM/CT 670-Pro product release	Conformance Statement Overview 1.1 1.6 2.1 2.2.1 2.3 (Note)	M.Mesh



			2.3.1.2.1 (Note) 2.3.1.2.2 (Note) 2.3.1.2.3 (Note) 2.3.1.2.4 (Note) 3.1 3.4.4.1 3.4.5.1 3.4.6.10 3.4.6.15 3.4.6.16 4.1 5.1 6.1 7.1 7.4.4 8.1 (Note)	
8	September 2014	Optima NM/CT 640 product release	Conformance Statement Overview 2.2.3 2.2.4.4 2.3.1.2.2.2.1 2.3.1.2.5.1 2.3.1.2.5.2.1 3.4.6.7 3.4.6.9 3.4.6.10 9.3.2.1 9.3.3.1	M. Mesh
9	August 2015	Discovery NM/CT 670 ES product Release	Conformance statement overview 1.1 2.1 2.2.1 2.3.1.2.1 (Note) 3.1 (Note) 3.4.4.1 3.4.5.1 3.4.6.10 3.4.6.15 3.4.6.16 4.1 (Note) 5.1 6.1 (Note) 7.1 7.4.4 8.1 (Note)	M. Mesh
10	November 2015	Discovery NM/CT 670 CZT product Release	Conformance statement overview 1.1 2.1 2.2.1 2.3.1.2.1 (Note)	M.Mesh



			<ul style="list-style-type: none"> <li>3.1 (Note)</li> <li>3.4.4.1</li> <li>3.4.5.1</li> <li>3.4.6.10</li> <li>3.4.6.15</li> <li>3.4.6.16</li> <li>4.1 (Note)</li> <li>5.1</li> <li>6.1 (Note)</li> <li>7.1</li> <li>7.4.4</li> <li>8.1 (Note)</li> </ul>	
11	August 2016	<p>New release of NM 600/800 Series Software</p> <p>New types of Secondary Capture</p> <p>Radiopharmaceutical Radiation Dose SR</p>	<p>Conformance statement overview</p> <ul style="list-style-type: none"> <li>2.3.1</li> <li>2.3.1.2.1.1</li> <li>3.4.3.4</li> <li>3.4.5.1</li> <li>3.4.6.6</li> <li>4</li> <li>10 (New)</li> </ul>	M.Mesh
12	January 2017	Discovery 670 DR product Release	<p>Conformance statement overview</p> <ul style="list-style-type: none"> <li>1.1</li> <li>3.1 Note</li> <li>4.3.5.1</li> <li>7.1</li> <li>10.3.2.1</li> </ul>	M.Mesh
13	November 2017	<p>800 Series products release</p> <p>New private Tags</p>	<p>Conformance statement overview</p> <ul style="list-style-type: none"> <li>1.1</li> <li>1.6</li> <li>2.1</li> <li>2.2.1</li> <li>2.3</li> <li>2.3.1</li> <li>2.3.1.2.1</li> <li>2.3.1.2.2</li> <li>2.3.1.2.3</li> <li>2.3.1.2.4</li> <li>3.1</li> <li>3.4.6</li> <li>3.4.6.10</li> <li>3.4.6.16</li> <li>3.4.6.17</li> <li>3.5.1</li> <li>4.1</li> <li>5.1</li> <li>7.1</li> <li>7.4.4</li> <li>8.1</li> <li>10.1</li> </ul>	M.Mesh
14	December 2017	New release of NM 600/800 Series	Conformance Statement Overview	M.Mesh



		Software New private Tags	1.1 1.2 2.3.1 2.3.1.2.1.2 2.5.3 (New) 2.5.4 (New) 3.3 3.4.3.13.4.6.10 3.4.6.16 3.5.2 5.4 5.4.1.2 5.4.2.1 8.5.4 9.2.1 9.3.1.1 9.3.3.5 10.2.1 10.3.1.1 10.3.3.5 10.6.1.2.1	
15	May 2018	New release of NM 600/800 Series Software	Conformance Statement Overview 1.1 1.2 2.3.1 2.3.1.2.1.2 2.3.1.2.2.2 2.5.3 2.5.4 2.8.1 3.3 3.4.3.1 3.4.6.8 3.4.6.9 3.4.6.10 3.4.6.16 3.5.2 5.4 5.4.1.2 5.4.2.1 8.5.4 9.2.1 9.3.1.1 9.3.3.5 10.2.1 10.3.1.1 10.3.2.1 Table 10-6 10.3.3.5 10.6.1.2.1	M.Mesh

**THIS PAGE LEFT INTENTIONALLY BLANK**

## CONFORMANCE STATEMENT OVERVIEW

Throughout this document the term "NM General Purpose 600/800 Series" refers to the following products:

- **Discovery NM/CT 670**
- **Discovery NM/CT 670 Pro**
- **Discovery NM/CT 670 ES**
- **Discovery NM/CT 670 CZT**
- **Discovery 670 DR**
- **Discovery NM 630**
- **Optima NM/CT 640**
- **Brivo NM 615**
- **NM/CT 850**
- **NM/CT 860**
- **NM 830**
- **NM/CT 870 DR**
- **NM/CT 870 CZT**

The Discovery NM/CT 670 system is a combination of an NM camera and a BrightSpeed Elite 16 slices CT Scanner.

The Discovery NM/CT 670 Pro, Discovery NM/CT 670 CZT ,Discovery 670 DR ,NM/CT 870 DR and NM/CT 870 CZT systems are a combination of an NM camera and an Optima CT540 16 Slices CT scanner.

The Discovery NM/CT 670 ES system is a combination of an NM camera and an Optima CT540 8 slices CT scanner.

NM/CT 850 system is combination of NM Camera and CT 850 8 Slices CT Scanner.

NM/CT 860 system is combination of NM Camera and CT 860 8 Slices CT Scanner.

**Discovery NM/CT 670\*** term is used when functionality described is the same for all Hybrid 600 Series configurations (Discovery NM/CT 670, Discovery NM/CT 670 Pro, Discovery NM/CT 670 ES, Discovery 670 DR and Discovery NM/CT 670 CZT).

**NM/CT 800\*** term is when functionality described is the same for all Hybrid 800 Series configurations (NM/CT 850, NM/CT 860, NM/CT 870 DR and NM/CT 870 CZT).

**NM/CT 870\*** term is when functionality described is the same for Hybrid800 Series configurations (NM/CT 850, NM/CT 860, NM/CT 870 DR and NM/CT 870 CZT).

Discovery NM/CT 670\* and NM/CT 800\* systems work in three different modes:

1. NM standalone - Generates NM images.
2. NM/CT hybrid - A single NM scan (Tomo, Gated Tomo only) partnered with the related Hybrid CT protocol with one or more CT acquisition(s), NM and CT series will be included in the same study.
3. CT standalone exams scheduled via CT UI. Generates CT images only.

The Discovery NM/CT 670\* DICOM implementations allow the user to send Nuclear Medicine image data, which may be partnered with CT image data acquired through the front-end acquisition system, and Secondary Capture Objects, created as reports of several Quality Control (QC) operations on the NM camera, Radiopharmaceutical Radiation Dose Structured Reports (RRDR) and RRDR Secondary Capture Objects to another DICOM station. Note that CT images and secondary capture are sent by the CT scanner independently (See [Reference A](#) in Section 1.6).

The NM/CT 800\* DICOM implementations allow the user to send Nuclear Medicine image data, which may be partnered with CT image data acquired through the front-end acquisition system, and Secondary Capture Objects,



DIR DOC0660253 REV 15

created as reports of several Quality Control (QC) operations on the NM camera, Radiopharmaceutical Radiation Dose Structured Reports (RRDR) and RRDR Secondary Capture Objects to another DICOM station. Note that CT images, CT Dose Reports and secondary captures are sent by the CT 800\* scanner independently (See [Reference B](#) in Section 1.6).

The Discovery NM 630 system is NM Standalone camera, which generates only NM Images, Radiopharmaceutical Radiation Dose Structured Reports (RRDR) and RRDR Secondary Capture Object Its DICOM implementation is the same as the Discovery NM/CT 670\* systems working in NM standalone mode.

The NM 830 system is NM Standalone camera, which generates only NM Images, Radiopharmaceutical Radiation Dose Structured Reports (RRDR) and RRDR Secondary Capture Object. Its DICOM implementation is the same as the NM/CT 800\* systems working in NM standalone mode.

The Brivo NM 615 is single-detector NM Standalone camera derived from Discovery NM 630, which generates only NM Images It's DICOM implementation is the same as Discovery NM 630.

The NM 830, Discovery NM 630 and Brivo NM 615 DICOM implementations allows the user to send Nuclear Medicine image data acquired through the front-end acquisition system, and Secondary Capture Objects, created as reports of several Quality Control (QC) operations on the NM camera, to another DICOM station.

The Optima NM/CT 640 system is NM/CT Hybrid camera, which generates both NM and CT images .It works in 2 modes:

1. NM standalone - Generates NM images.
2. NM/CT hybrid - A single NM scan (Tomo, Gated Tomo only) partnered with the related Hybrid CT protocol with one or more CT acquisition(s), NM and CT series will be included in the same study

The Optima NM/CT 640 DICOM implementation allows the user to send Nuclear Medicine image data and Computerized Tomography image data acquired through the front-end acquisition system, CT Radiation Dose reports of patient exposure during Hybrid acquisition and Secondary Capture Objects, created as reports of several Quality Control (QC) operations on the, to another DICOM station.

**NM General Purpose 600/800 Series** DICOM implementations support storage commitment for the already transferred data. This guarantees the user that the acquired data is safely archived for future use.

**NM General Purpose 600/800 Series** DICOM implementations also support receiving Worklist information from a remote AE. They support receiving more than one Scheduled procedure step per study instance. Similarly, **NM General Purpose 600/800 Series** products support scheduling locally more than one protocol to be performed for a study.

**NM General Purpose 600/800 Series** DICOM implementation also provide a verification mechanism by which a remote AE can verify application-level communication **NM General Purpose 600/800 Series** DICOM Server. Also provided is a mechanism by **NM General Purpose 600/800 Series** user can verify application-level communication with a remote AE.

**NM General Purpose 600/800 Series** DICOM implementation creates and updates Modality Performed Procedure Step instances managed by a remote AE in association with image acquisition. Completion or discontinuation of the MPPS is performed as the result of an operator action.

**Note:** Modality Performed Procedure Step SOP Class is supported by Discovery NM/CT 670 system in CT Standalone mode via BrightSpeed Elite 16 slices CT Scanner, Discovery NM/CT 670 ES system in CT Standalone mode via Optima CT540 8 Slices CT scanner, by Discovery NM/CT 670 Pro, Discovery 670 DR, Discovery NM/CT 670 CZT and NM/CT 870\* systems in CT Standalone mode via Optima CT540 16 slices CT scanner, by NM/CT 850 system in standalone mode by CT 850 Scanner and by NM/CT 860 system in standalone mode by CT 860 Scanner. Refer to appropriate CT Scanner DICOM Conformance Statements (Section 1.6).

DIR DOC0660253 REV 15

Table 0.1 provides an overview of the network services supported by NM General Purpose 600/800 Series

Table 0.1 – NETWORK SERVICES

SOP Classes	User of Service (SCU)	Provider of Service (SCP)
<b>Transfer</b>		
Secondary Capture Image Storage	Yes	No
Nuclear Medicine Image Storage	Yes	No
Computerized Tomography Image Storage	Yes (*)	No
X-Ray Radiation Dose SR	Yes (****)	No
Enhanced SR	Yes (**)	No
Radiopharmaceutical Radiation Dose SR	Yes (***)	No
<b>Workflow Management</b>		
Storage Commitment Push Model SOP Class	Yes	No
Modality Worklist Information Model – FIND SOP Class	Yes	No
Modality Performed Procedure Step SOP Class	Yes	No

**Note:** (\*)– Discovery NM/CT 670\*,NM/CT 800\* and Optima NM/CT 640 DICOM implementations only.  
 (\*\*) – Discovery NM/CT 670\* and NM/CT 800\* implementations (for CT Radiation Dose SR objects created by CT Scanner) and Optima NM/CT 640 DICOM implementation (for transfer CT Radiation Dose SR Objects to systems, which do not support X-Ray Radiation Dose SR SOP Class)  
 (\*\*\*) – Except Brivo NM 615 DICOM implementation  
 (\*\*\*\*) - Optima NM/CT 640 DICOM implementations only

<b>1. INTRODUCTION .....</b>	<b>18</b>
1.1 Overview .....	18
1.2 Overall DICOM Conformance Statement Document Structure.....	21
1.3 Intended Audience .....	22
1.4 Scope and Field of Application .....	22
1.5 Important Remarks .....	22
1.6 References.....	23
1.7 Definitions.....	23
1.8 Symbols and Abbreviations.....	25
<b>2. NETWORK CONFORMANCE STATEMENT .....</b>	<b>26</b>
2.1 Introduction.....	26
2.2 Implementation Model.....	27
2.2.1 Overview.....	27
2.2.2 Application Data Flow Diagram .....	28
2.2.3 Functional Definition of AE's.....	29
2.2.4 Sequencing of Real-World Activities.....	29
2.2.4.1 Default Protocol Workflow.....	29
2.2.4.2 Local Scheduling.....	30
2.2.4.3 Modality Worklist (MWL).....	30
2.2.4.4 Report of Performed Procedure Step progress .....	30
2.2.4.5 Send Images (Manual) .....	30
2.3 AE Specifications .....	31
2.3.1 NM Camera AE Specification.....	31
2.3.1.1 Association Establishment Policies.....	31
2.3.1.2 Association Initiation Policy .....	32
2.3.1.3 Association Acceptance Policy .....	40
2.4 Communication Profiles .....	43
2.4.1 Supported Communication Stacks .....	43
2.4.2 Physical Media Support .....	43
2.5 Extensions / Specializations/ Privatizations .....	43
2.5.1 Standard Extended / Specialized / Private SOP Classes .....	43
2.5.1.1 Standard Extended SOP Classes .....	43
2.5.1.2 Private SOP Classes .....	43
2.5.2 Private Transfer Syntaxes .....	43
2.5.3 Additional Protocols .....	43

DIR DOC0660253 REV 15

2.5.4	IPv4 and IPv6 Support .....	43
<b>2.6</b>	<b>Configuration .....</b>	<b>43</b>
2.6.1	AE Title/Presentation Address Mapping .....	43
2.6.2	Configurable Parameters .....	43
<b>2.7</b>	<b>Support of Extended Character Sets .....</b>	<b>44</b>
<b>2.8</b>	<b>Codes and Controlled Terminology .....</b>	<b>44</b>
2.8.1	Fixed Coded Terminology .....	44
2.8.2	Mapped Coded Terminology .....	45
2.8.3	Configurable Coded Terminology .....	45
<b>2.9</b>	<b>Security Profiles .....</b>	<b>45</b>
<b>3.</b>	<b>NM INFORMATION OBJECT IMPLEMENTATION .....</b>	<b>46</b>
<b>3.1</b>	<b>Introduction .....</b>	<b>46</b>
<b>3.2</b>	<b>NM Camera Mapping of DICOM Entities .....</b>	<b>46</b>
<b>3.3</b>	<b>IOD Module Table .....</b>	<b>46</b>
<b>3.4</b>	<b>Information Module Definitions .....</b>	<b>48</b>
3.4.1	Patient Entity Modules .....	48
3.4.1.1	Patient Module .....	48
3.4.2	Study Entity Modules .....	49
3.4.2.1	General Study Module .....	49
3.4.2.2	Patient Study Module .....	51
3.4.2.3	Standard Extended Study Module .....	51
3.4.2.4	Private Study Module .....	52
3.4.3	Series Entity Modules .....	53
3.4.3.1	General Series Module .....	53
3.4.3.2	Standard Extended Series Module .....	56
3.4.3.3	NM/PET Patient Orientation Module .....	56
3.4.3.4	Private Series Module .....	58
3.4.4	Frame Of Reference Entity Modules .....	60
3.4.4.1	Frame Of Reference Module .....	60
3.4.5	Equipment Entity Modules .....	61
3.4.5.1	General Equipment Module .....	61
3.4.6	Image Entity Modules .....	63
3.4.6.1	General Image Module .....	63
3.4.6.2	Image Pixel Module .....	64
3.4.6.3	Acquisition Context Module .....	64
3.4.6.4	NM Image Pixel Module .....	65
3.4.6.5	Multi-Frame Module .....	66
3.4.6.6	NM Multi-frame Module .....	66
3.4.6.7	NM Image Module .....	68
3.4.6.8	NM Isotope Module .....	71
3.4.6.9	NM Detector Module .....	72
3.4.6.10	NM Tomo Acquisition Module .....	74
3.4.6.11	NM Multi-gated Acquisition Module .....	75

DIR DOC0660253 REV 15

3.4.6.12	NM Phase Module .....	76
3.4.6.13	VOI LUT Module .....	77
3.4.6.14	SOP Common Module .....	77
3.4.6.15	Private Image Module .....	77
3.4.6.16	Private NM Image Module.....	78
3.4.6.17	Private Image Tomo Module .....	81
3.4.6.18	Private Image Multi-Gated Module .....	82
3.4.6.19	Private Image GSPECT Module .....	82
<b>3.5</b>	<b>Standard Extended and Private Data Attributes .....</b>	<b>82</b>
3.5.1	Standard Extended Attributes .....	82
3.5.2	Private Group QUASAR_INTERNAL_USE.....	83
3.5.3	Private Group GEMS_GENIE_1 .....	85
3.5.4	Private Group GEMS_XELPRV_01.....	86
3.5.5	Private Group APEX_PRIVATE.....	87
<b>4.</b>	<b>SECONDARY CAPTURE INFORMATION OBJECT IMPLEMENTATION .....</b>	<b>88</b>
<b>4.1</b>	<b>Introduction.....</b>	<b>88</b>
<b>4.2</b>	<b>NM Camera Mapping of DICOM Entities .....</b>	<b>88</b>
<b>4.3</b>	<b>IOD Module Table .....</b>	<b>89</b>
<b>4.4</b>	<b>Information Module Definitions .....</b>	<b>90</b>
4.4.1	Patient Entity Modules.....	90
4.4.1.1	Patient Module .....	90
4.4.2	Study Entity Modules.....	91
4.4.2.1	General Study Module .....	91
4.4.3	Series Entity Modules .....	91
4.4.3.1	General Series Module.....	91
4.4.3.2	Private SC Series.....	92
4.4.4	Equipment Entity Modules.....	93
4.4.4.1	SC Equipment Module.....	93
4.4.5	Image Entity Modules .....	93
4.4.5.1	General Image Module.....	93
4.4.5.2	Image Pixel Module .....	94
4.4.5.3	SC Image Module .....	95
4.4.5.4	SOP Common Module .....	95
4.4.5.5	Standard Extended Image Module .....	96
4.4.5.6	Private SC Image Module .....	96
<b>4.5</b>	<b>Standard Extended and Private Data Attributes .....</b>	<b>96</b>
4.5.1	Standard Extended Attributes .....	96
4.5.2	Private Group QUASAR_INTERNAL_USE.....	96
4.5.3	Private Group GEMS_GENIE_1 .....	97
<b>4.6</b>	<b>Standard Extended and Private Context Groups .....</b>	<b>97</b>
<b>5.</b>	<b>MODALITY WORKLIST QUERY IMPLEMENTATION .....</b>	<b>98</b>

DIR DOC0660253 REV 15

<b>5.1</b>	<b>Introduction.....</b>	<b>98</b>
<b>5.2</b>	<b>NM Camera Mapping of DICOM Entities .....</b>	<b>98</b>
<b>5.3</b>	<b>Worklist Query Module Table.....</b>	<b>99</b>
<b>5.4</b>	<b>Worklist Query Module Definitions .....</b>	<b>99</b>
5.4.1	Common Scheduled Procedure Step Entity Modules.....	99
5.4.1.1	SOP Common Module .....	99
5.4.1.2	Scheduled Procedure Step Module .....	100
5.4.2	Common Requested Procedure Entity Modules.....	101
5.4.2.1	Requested Procedure Module .....	101
5.4.3	Common Imaging Service Request Entity Modules.....	103
5.4.3.1	Imaging Service Request Module .....	103
5.4.4	Common visit Entity Modules.....	103
5.4.4.1	Visit Identification .....	103
5.4.4.2	Visit Status .....	103
5.4.4.3	Visit Relationship.....	104
5.4.4.4	Visit Admission.....	104
5.4.5	Common Patient Entity Modules .....	104
5.4.5.1	Patient Relationship .....	104
5.4.5.2	Patient Identification .....	105
5.4.5.3	Patient Demographic.....	105
5.4.5.4	Patient Medical .....	106
<b>6.</b>	<b>STORAGE COMMITMENT PUSH MODEL IMPLEMENTATION.....</b>	<b>107</b>
<b>6.1</b>	<b>Storage Commitment Push Model Information Object Definition.....</b>	<b>107</b>
6.1.1	STORAGE COMMITMENT MODULE FOR N-ACTION.....	107
6.1.2	STORAGE COMMITMENT MODULE FOR N-EVENT-REPORT .....	108
6.1.2.1	Processing of Failure Reason when received in a N-Event-Report.....	108
<b>7.</b>	<b>MODALITY PERFORMED PROCEDURE STEP IMPLEMENTATION .....</b>	<b>109</b>
<b>7.1</b>	<b>Introduction.....</b>	<b>109</b>
<b>7.2</b>	<b>Relationship between Scheduled and Performed Procedure Steps .....</b>	<b>109</b>
<b>7.3</b>	<b>Modality Performed Procedure Step Module Table.....</b>	<b>109</b>
<b>7.4</b>	<b>Modality Performed Procedure Step Module Definitions.....</b>	<b>110</b>
7.4.1	SOP Common Module .....	110
7.4.2	Performed Procedure Step Relationship Module .....	110
7.4.3	Performed Procedure Step Information Module .....	112
7.4.4	Image Acquisition Results Module .....	114
<b>7.5</b>	<b>Standard Extended and Private Data Attributes .....</b>	<b>116</b>
<b>7.6</b>	<b>Standard Extended and Private Context Groups .....</b>	<b>116</b>

**8. CT INFORMATION OBJECT IMPLEMENTATION .....117**

**8.1 Introduction..... 117**

**8.2 NM Camera Mapping of DICOM Entities ..... 117**

**8.3 IOD Module Table ..... 117**

**8.4 Information Module Definitions ..... 118**

8.4.1 Series Entity Modules ..... 119

8.4.1.1 Private CT Series Module ..... 119

8.4.2 Frame Of Reference Entity Modules ..... 119

8.4.2.1 Frame Of Reference Module ..... 119

8.4.3 Image Entity Modules ..... 119

8.4.3.1 General Image Module ..... 119

8.4.3.2 Image Plane Module ..... 120

8.4.3.3 Image Pixel Module ..... 120

8.4.3.4 Contrast/Bolus Module ..... 121

8.4.3.5 Standard Extended Image Module ..... 121

8.4.3.6 CT Image Module ..... 122

8.4.3.7 VOI LUT module ..... 128

8.4.3.8 SOP Common Module ..... 129

8.4.3.9 Private CT Image Module ..... 129

**8.5 Standard Extended and Private Data Attributes ..... 131**

8.5.1 Standard Attributes ..... 131

8.5.2 Private Group GEMS\_GENIE\_1 ..... 131

8.5.3 Private Group QUASAR\_INTERNAL\_USE ..... 132

8.5.4 Private Group GEHC\_HYBRID\_01 ..... 132

**9. CT RADIATION DOSE STRUCTURED REPORT INFORMATION OBJECT IMPLEMENTATION ..... 133**

**9.1 Introduction..... 133**

**9.2 NM Camera Mapping of DICOM Entities ..... 133**

9.2.1 IOD Module Table ..... 133

**9.3 Information Module Definitions ..... 134**

9.3.1 Series Entity Modules ..... 134

9.3.1.1 SR Document Series Module ..... 134

9.3.1.2 Standard Extended Series Module ..... 135

9.3.1.3 Private Series Module ..... 136

9.3.2 Equipment Entity Modules ..... 137

9.3.2.1 Enhanced General Equipment Module ..... 137

9.3.3 Document Entity Modules ..... 138

9.3.3.1 SR Document General Module ..... 138

9.3.3.2 SR Document Content Module ..... 138

9.3.3.3 Standard Extended SR Document Module ..... 140

9.3.3.4 Private SR Document Module ..... 140

9.3.3.5 SOP Common Module ..... 140

DIR DOC0660253 REV 15

<b>9.4</b>	<b>Standard Extended and Private Data Attributes .....</b>	<b>141</b>
9.4.1	Standard Attributes .....	141
9.4.2	Private Group QUASAR_INTERNAL_USE.....	141
9.4.3	Private Group GEMS_GENIE_1 .....	142
<b>9.5</b>	<b>Standard Extended and Private Context Groups .....</b>	<b>142</b>
<b>9.6</b>	<b>Standard, Standard Extended and Private Templates .....</b>	<b>142</b>
9.6.1	Standard Templates.....	142
9.6.1.1	Template ID 10011 CT Radiation Dose.....	142
9.6.1.2	TID 1002 Observer Context.....	144
9.6.1.3	TID 1004 Device Observer Identifying Attributes.....	144
9.6.1.4	TID 10012 CT Accumulated Dose Data .....	145
9.6.1.5	TID 10013 CT Irradiation Event Data .....	146
9.6.1.6	TID 10014 Scanning Length.....	149
9.6.1.7	TID 10015 CT Dose Check Details .....	149
9.6.2	Standard Extended Templates.....	152
9.6.3	Private Templates.....	152
<b>10.</b>	<b>RADIOPHARMACEUTICAL RADIATION DOSE STRUCTURED REPORT INFORMATION OBJECT IMPLEMENTATION .....</b>	<b>153</b>
<b>10.1</b>	<b>Introduction.....</b>	<b>153</b>
<b>10.2</b>	<b>NM Camera Mapping of DICOM Entities .....</b>	<b>153</b>
10.2.1	IOD Module Table.....	153
<b>10.3</b>	<b>Information Module Definitions .....</b>	<b>154</b>
10.3.1	Series Entity Modules .....	154
10.3.1.1	SR Document Series Module .....	154
10.3.1.2	Standard Extended Series Module .....	155
10.3.1.3	Private Series Module .....	155
10.3.2	Equipment Entity Modules.....	156
10.3.2.1	Enhanced General Equipment Module.....	156
10.3.3	Document Entity Modules .....	156
10.3.3.1	SR Document General Module .....	156
10.3.3.2	SR Document Content Module .....	157
10.3.3.3	Standard Extended SR Document Module.....	158
10.3.3.4	Private SR Document Module.....	158
10.3.3.5	SOP Common Module .....	158
<b>10.4</b>	<b>Standard Extended and Private Data Attributes .....</b>	<b>159</b>
10.4.1	Standard Attributes .....	159
10.4.2	Private Group QUASAR_INTERNAL_USE.....	159
10.4.3	Private Group GEMS_GENIE_1 .....	160
<b>10.5</b>	<b>Standard Extended and Private Context Groups .....</b>	<b>160</b>
<b>10.6</b>	<b>Standard, Standard Extended and Private Templates .....</b>	<b>160</b>
10.6.1	Standard Templates.....	160
10.6.1.1	TID 10021 Radiopharmaceutical Radiation Dose.....	160
10.6.1.2	TID 10022 Radiopharmaceutical Administration Event Data.....	161



DIR DOC0660253 REV 15

10.6.1.3	TID 10024 Radiopharmaceutical Administration Patient Characteristics.....	165
10.6.2	Standard Extended Templates.....	165
10.6.3	Private Templates.....	166

## 1. INTRODUCTION

### 1.1 OVERVIEW

Throughout this document the term "NM General Purpose 600/800 Series" refers to the following products:

- **Discovery NM/CT 670**
- **Discovery NM/CT 670 Pro**
- **Discovery NM/CT 670 CZT**
- **Discovery NM/CT 670 ES**
- **Discovery 670 DR**
- **Discovery NM 630**
- **Optima NM/CT 640**
- **Brivo NM 615**
- **NM/CT 850**
- **NM/CT 860**
- **NM 830**
- **NM/CT 870 DR**
- **NM/CT 870 CZT**

**Discovery NM/CT 670\*** term is used when functionality described is the same for all Hybrid 600 Series configurations (Discovery NM/CT 670, Discovery NM/CT 670 Pro, Discovery NM/CT 670 ES, Discovery 670 DR and Discovery NM/CT 670 CZT).

**NM/CT 800\*** term is when functionality described is the same for all Hybrid 800 Series configurations (NM/CT 850, NM/CT 860, NM/CT 870 DR and NM/CT 870 CZT).

**NM/CT 870\*** term is when functionality described is the same for Hybrid800 Series configurations (NM/CT 850, NM/CT 860, NM/CT 870 DR and NM/CT 870 CZT).

The Discovery NM/CT 670 system is a combination of an NM camera and a BrightSpeed Elite 16 slices CT Scanner.

The Discovery NM/CT 670 Pro, Discovery NM/CT 670 CZT ,Discovery 670 DR ,NM/CT 870 DR and NM/CT 870 CZT systems are a combination of an NM camera and an Optima CT540 16 Slices CT scanner.

The Discovery NM/CT 670 ES system is a combination of an NM camera and an Optima CT540 8 slices CT scanner.

NM/CT 850 system is combination of NM Camera and CT 850 8 Slices CT Scanner.

NM/CT 860 system is combination of NM Camera and CT 860 8 Slices CT Scanner.

Discovery NM/CT 670\* and NM/CT 800\* systems work in three different modes:

- 1) NM standalone - Generates NM images.
- 2) NM/CT hybrid - A single NM scan (Tomo, Gated Tomo only) partnered with the related Hybrid CT protocol with one or more CT acquisition(s), NM and CT series will be included in the same study.
- 3) CT standalone exams scheduled via CT UI. Generates CT images only.

The Discovery NM/CT 670\* DICOM implementations allow the user to send Nuclear Medicine image data, which may be partnered with CT image data acquired through the front-end acquisition system, and Secondary Capture Objects, created as reports of several Quality Control (QC) operations on the NM camera, Radiopharmaceutical Radiation Dose Structured Reports (RRDR) and RRDR Secondary Capture Objects to another DICOM station.

The NM/CT 800\* DICOM implementations allow the user to send Nuclear Medicine image data, which may be partnered with CT image data acquired through the front-end acquisition system, and Secondary Capture Objects, created as reports of several Quality Control (QC) operations on the NM camera, Radiopharmaceutical Radiation Dose Structured Reports (RRDR) and RRDR Secondary Capture Objects to another DICOM station. Note that CT images, CT Dose Reports and secondary captures are sent by the CT 800\* scanner independently. Discovery NM/CT 670\* and NM/CT 870\* CT scanner DICOM Conformance Statement is described in a separate document ([Reference A](#) in Section 1.6).

NM/CT 850 and NM/CT 860 CT Scanner DICOM Conformance Statement is described in a separate document ([Reference B](#) in Section 1.6).

The Discovery NM 630 system is NM Standalone camera, which generates only NM Images, Radiopharmaceutical Radiation Dose Structured Reports (RRDR) and RRDR Secondary Capture Object. Its DICOM implementation is the same as the Discovery NM/CT 670\* systems working in NM standalone mode.

The NM 830 system is NM Standalone camera, which generates only NM Images, Radiopharmaceutical Radiation Dose Structured Reports (RRDR) and RRDR Secondary Capture Object. Its DICOM implementation is the same as the NM/CT 800\* systems working in NM standalone mode.

The Optima NM/CT 640 system is NM/CT Hybrid camera, which generates both NM and CT images. It works in 2 modes:

- 1) NM standalone - Generates NM images.
- 2) NM/CT hybrid - A single NM scan (Tomo, Gated Tomo only) partnered with the related Hybrid CT protocol with one or more CT acquisition(s), NM and CT series will be included in the same study.

The Brivo NM 615 is single-detector NM Standalone camera derived from Discovery NM 630. Its DICOM implementation is the same as Discovery NM 630.

The DICOM Conformance related to the NM standalone functionality as well as DICOM behavior of the NM/CT hybrid scanning is described in this document.

This DICOM Conformance Statement is divided into Sections as described below:

**Section 1 (Introduction)**, which describes the overall structure, intent, and references for this Conformance Statement

**Section 2 (Network Conformance Statement)**, which specifies the GEHC equipment compliance to the DICOM requirements for the implementation of Networking features.

**Section 3 (NM Image Information Object Implementation)**, which specifies the GEHC equipment compliance to DICOM requirements for the implementation of a NM Image Information Object.

**Section 4 (Secondary Capture Image Information Object Implementation)**, which specifies the GEHC equipment compliance to DICOM requirements for the implementation of a SC Image Information Object.

**Section 5 (Modality Worklist Information Model)**, which specifies the GEHC equipment compliance to DICOM requirements for the implementation of the Modality Worklist service.

**Section 6 (Storage Commitment PUSH Model Implementation)**, which specifies the GEHC equipment compliance to DICOM requirements for the implementation of the Storage Commitment service.

DIR DOC0660253 REV 15

**Section 7 (Modality Performed Procedure Step Model Implementation)**, which specifies the GEHC equipment compliance to DICOM requirements for the implementation of the Modality Performed Procedure Step service.


**Section 8 (CT Image Information Object Implementation)**, which specifies the GEHC equipment compliance to DICOM requirements for the Optima NM/CT 640 implementation of a CT Image Information Object.

**Section 9 (CT Radiation Dose Structured Report Information Object Implementation)**, which specifies the GEHC equipment compliance to DICOM requirements for the Optima NM/CT 640 implementation of a CT Radiation Dose SR Information Object.

**Section 10 (Radiopharmaceutical Radiation Dose Structured Report Information Object Implementation)**, which specifies the GEHC equipment compliance to DICOM requirements for the implementation of a Radiopharmaceutical Radiation Dose SR Information Object.

1.2 OVERALL DICOM CONFORMANCE STATEMENT DOCUMENT STRUCTURE

The Documentation Structure of the GEHC DICOM Conformance Statements is shown in the Illustration below.

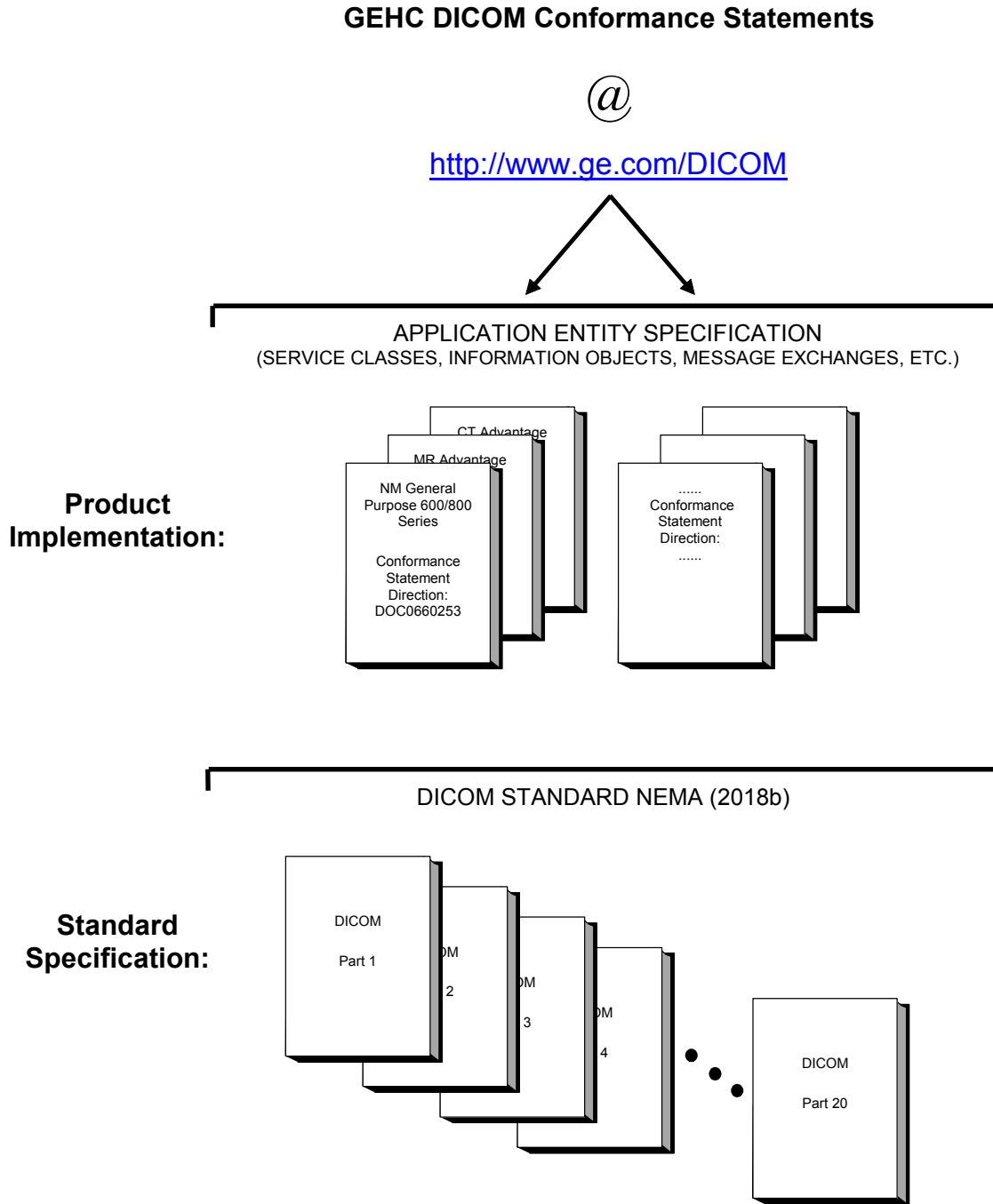

DIR DOC0660253 REV 15

This document specifies the DICOM implementation. It is entitled:

**NM General Purpose 600/800 Series**  
 Conformance Statement for DICOM  
 Direction **DOC0660253**

This DICOM Conformance Statement documents the DICOM Conformance Statement and Technical Specification required interoperating with the GEHC network interface.

The GEHC Conformance Statement, contained in this document, also specifies the Lower Layer communications which it supports (e.g., TCP/IP). However, the Technical Specifications are defined in the DICOM Part 8 standard.

For more information regarding DICOM, copies of the Standard may be obtained on the Internet at <http://medical.nema.org>. Comments on the Standard may be addressed to:

DICOM Secretariat  
 NEMA  
 1300 N. 17<sup>th</sup> Street, Suite 1752  
 Rosslyn, VA 22209  
 USA  
 Phone: +1.703.841.3200

### 1.3 INTENDED AUDIENCE

The reader of this document is concerned with software design and/or system integration issues. It is assumed that the reader of this document is familiar with the DICOM Standard and with the terminology and concepts which are used in that Standard.

### 1.4 SCOPE AND FIELD OF APPLICATION

It is the intent of this document to provide an unambiguous specification for GEHC implementations. This specification, called a Conformance Statement, includes a DICOM Conformance Statement and is necessary to ensure proper processing and interpretation of GEHC medical data exchanged using DICOM. The GEHC Conformance Statements are available to the public.

The reader of this DICOM Conformance Statement should be aware that different GEHC devices are capable of using different Information Object Definitions. For example, a GEHC CT Scanner may send images using the CT Information Object, MR Information Object, Secondary Capture Object, etc.

Included in this DICOM Conformance Statement are the Module Definitions which define all data elements used by this GEHC implementation. If the user encounters unspecified private data elements while parsing a GEHC Data Set, the user is well advised to ignore those data elements (per the DICOM standard). Unspecified private data element information is subject to change without notice. If, however, the device is acting as a "full fidelity storage device", it should retain and re-transmit all of the private data elements which are sent by GEHC devices.

### 1.5 IMPORTANT REMARKS

The use of these DICOM Conformance Statements, in conjunction with the DICOM Standards, is intended to facilitate communication with GE imaging equipment. However, **by itself, it is not sufficient to ensure that inter-operation will be successful**. The **user (or user's agent)** needs to proceed with caution and address at least four issues:

- **Integration** - The integration of any device into an overall system of interconnected devices goes beyond the scope of standards (DICOM v3.0), and of this introduction and associated DICOM Conformance Statements when interoperability with non-GE equipment is desired. The responsibility to analyze the applications requirements and to design a solution that integrates GE imaging equipment with

DIR DOC0660253 REV 15

non-GE systems is the **user's** responsibility and should not be underestimated. The **user** is strongly advised to ensure that such an integration analysis is correctly performed.

- **Validation** - Testing the complete range of possible interactions between any GE device and non-GE devices, before the connection is declared operational, should not be overlooked. Therefore, the **user** should ensure that any non-GE provider accepts full responsibility for all validation required for their connection with GE devices. This includes the accuracy of the image data once it has crossed the interface between the GE imaging equipment and the non-GE device and the stability of the image data for the intended applications.

Such a validation is required before any clinical use (diagnosis and/or treatment) is performed. It applies when images acquired on GE imaging equipment are processed/displayed on a non-GE device, as well as when images acquired on non-GE equipment is processed/displayed on a GE console or workstation.

- **Future Evolution** - GE understands that the DICOM Standard will evolve to meet the user's growing requirements. GE is actively involved in the development of the DICOM Standard. DICOM will incorporate new features and technologies and GE may follow the evolution of the Standard. The GEHC protocol is based on DICOM as specified in each DICOM Conformance Statement. Evolution of the Standard may require changes to devices which have implemented DICOM. **In addition, GE reserves the right to discontinue or make changes to the support of communications features (on its products) described by these DICOM Conformance Statements.** The **user** should ensure that any non-GE provider, which connects with GE devices, also plans for the future evolution of the DICOM Standard. Failure to do so will likely result in the loss of function and/or connectivity as the DICOM Standard changes and GE Products are enhanced to support these changes.

- **Interaction** - It is the sole responsibility of the **non-GE provider** to ensure that communication with the interfaced equipment does not cause degradation of GE imaging equipment performance and/or function.

## 1.6 REFERENCES

- |             |  |
|-------------|--|
| NEMA PS3    | Digital Imaging and Communications in Medicine (DICOM) Standard, available free at <a href="http://medical.nema.org/">http://medical.nema.org/</a> |
| Reference A | DOC0636565 ,“LightSpeed*, Discovery*, BrightSpeed* and Optima* Dicom Conformance Statement”  |
| Reference B | DOC2074801 NM-CT 850/860 Dicom Conformance Statement   |

## 1.7 DEFINITIONS

Informal definitions are provided for the following terms used in this Conformance Statement. The DICOM Standard is the authoritative source for formal definitions of these terms.

**Abstract Syntax** – the information agreed to be exchanged between applications, generally equivalent to a Service/Object Pair (SOP) Class. Examples: Verification SOP Class, Modality Worklist Information Model Find SOP Class, Computed Radiography Image Storage SOP Class.

**Application Entity (AE)** – an end point of a DICOM information exchange, including the DICOM network or media interface software; i.e., the software that sends or receives DICOM information objects or messages. A single device may have multiple Application Entities.

**Application Entity Title** – the externally known name of an *Application Entity*, used to identify a DICOM application to other DICOM applications on the network.

**Application Context** – the specification of the type of communication used between *Application Entities*. Example: DICOM network protocol.

DIR DOC0660253 REV 15

**Association** – a network communication channel set up between *Application Entities*.

**Attribute** – a unit of information in an object definition; a data element identified by a *tag*. The information may be a complex data structure (Sequence), itself composed of lower level data elements. Examples: Patient ID (0010,0020), Accession Number (0008,0050), Photometric Interpretation (0028,0004), Procedure Code Sequence (0008,1032).

**Information Object Definition (IOD)** – the specified set of *Attributes* that comprise a type of data object; does not represent a specific instance of the data object, but rather a class of similar data objects that have the same properties. The *Attributes* may be specified as Mandatory (Type 1), Required but possibly unknown (Type 2), or Optional (Type 3), and there may be conditions associated with the use of an Attribute (Types 1C and 2C). Examples: MR Image IOD, CT Image IOD, Print Job IOD.

**Joint Photographic Experts Group (JPEG)** – a set of standardized image compression techniques, available for use by DICOM applications.

**Module** – a set of *Attributes* within an *Information Object Definition* that are logically related to each other. Example: Patient Module includes Patient Name, Patient ID, Patient Birth Date, and Patient Sex.

**Negotiation** – first phase of *Association* establishment that allows *Application Entities* to agree on the types of data to be exchanged and how that data will be encoded.

**Presentation Context** – the set of DICOM network services used over an *Association*, as negotiated between *Application Entities*; includes *Abstract Syntaxes* and *Transfer Syntaxes*.

**Protocol Data Unit (PDU)** – a packet (piece) of a DICOM message sent across the network. Devices must specify the maximum size packet they can receive for DICOM messages.

**Security Profile** – a set of mechanisms, such as encryption, user authentication, or digital signatures, used by an *Application Entity* to ensure confidentiality, integrity, and/or availability of exchanged DICOM data

**Service Class Provider (SCP)** – role of an *Application Entity* that provides a DICOM network service; typically, a server that performs operations requested by another *Application Entity* (*Service Class User*). Examples: Picture Archiving and Communication System (image storage SCP, and image query/retrieve SCP), Radiology Information System (Modality Worklist SCP).

**Service Class User (SCU)** – role of an *Application Entity* that uses a DICOM network service; typically, a client. Examples: imaging modality (image storage SCU, and modality Worklist SCU), imaging workstation (image query/retrieve SCU)

**Service/Object Pair (SOP) Class** – the specification of the network or media transfer (service) of a particular type of data (object); the fundamental unit of DICOM interoperability specification. Examples: Ultrasound Image Storage Service, Basic Grayscale Print Management.

**Service/Object Pair (SOP) Instance** – an information object; a specific occurrence of information exchanged in a *SOP Class*. Examples: a specific x-ray image.

**Tag** – a 32-bit identifier for a data element, represented as a pair of four digit hexadecimal numbers, the “group” and the “element”. If the “group” number is odd, the tag is for a private (manufacturer-specific) data element. Examples: (0010,0020) [Patient ID], (07FE,0010) [Pixel Data], (0019,0210) [private data element]

**Transfer Syntax** – the encoding used for exchange of DICOM information objects and messages. Examples: *JPEG* compressed (images), Little Endian Explicit value representation.

**Unique Identifier (UID)** – a globally unique “dotted decimal” string that identifies a specific object or a class of objects; an ISO-8824 Object Identifier. Examples: Study Instance UID, SOP Class UID, SOP Instance UID.

**Value Representation (VR)** – the format type of an individual DICOM data element, such as text, an integer, a person’s name, or a code. DICOM information objects can be transmitted with either explicit identification of the type of each data element (Explicit VR), or without explicit identification (Implicit VR); with Implicit VR, the receiving application must use a DICOM data dictionary to look up the format of each data element.



1.8 SYMBOLS AND ABBREVIATIONS

AE	Application Entity
AET	Application Entity Title
CT	Computerized Tomography
DHCP	Dynamic Host Configuration Protocol
DICOM	Digital Imaging and Communications in Medicine
DNS	Domain Name System
HIS	Hospital Information System
IOD	Information Object Definition
ISO	International Organization for Standards
JPEG	Joint Photographic Experts Group
LUT	Look-up Table
MWL	Modality Worklist
MPPS	Modality Performed Procedure Step
NM	Nuclear Medicine
O	Optional (Key Attribute)
OSI	Open Systems Interconnection
PACS	Picture Archiving and Communication System
PDU	Protocol Data Unit
PPS	Performed Procedure Step
R	Required (Key Attribute)
RIS	Radiology Information System
RRDR	Radiopharmaceutical Radiation Dose Report
SC	Secondary Capture
SR	Structured Report
SCP	Service Class Provider
SCU	Service Class User
SDO	Series Data Object
SOP	Service-Object Pair
SPS	Scheduled Procedure Step
TCP/IP	Transmission Control Protocol/Internet Protocol
U	Unique (Key Attribute)
UI	User Interface
UL	Upper Layer
VM	Value Multiplicity
VR	Value Representation

## 2. NETWORK CONFORMANCE STATEMENT

### 2.1 INTRODUCTION

This section of the DICOM Conformance Statement specifies the NM General Purpose 600/800 Series compliance to DICOM requirements for networking features. Refer to Section 1 for detailed description of products of NM General Purpose 600/800 Series.

This section details the roles and DICOM Service Classes supported by the NM General Purpose 600/800 Series.

The Discovery NM/CT 670\* and NM/CT 800\* DICOM implementations allow the user to acquire, store Nuclear Medicine Image data, Radiopharmaceutical Radiation Dose Structured Reports (RRDR), RRDR Secondary Capture Objects, CT Image data, RRDR, CT Radiation Dose Reports and Secondary Capture Objects, created as reports of several Quality Control (QC) operations, and send them to another DICOM station.

The Discovery NM 630, NM 830 and Brivo NM 615 DICOM implementations allow the user to acquire, store Nuclear Medicine Image data, Radiopharmaceutical Radiation Dose Structured Reports (RRDR), RRDR Secondary Capture Objects and Secondary Capture Objects, created as reports of several Quality Control (QC) operations, and send them to another DICOM station.

The Brivo NM 615 DICOM implementation allows the user to acquire, store Nuclear Medicine Image data and Secondary Capture Objects, created as reports of several Quality Control (QC) operations, and send them to another DICOM station.

The Optima NM/CT 640 DICOM implementation allows the user to send acquired Nuclear Medicine image data and CT image data, Radiopharmaceutical Radiation Dose Structured Reports (RRDR), RRDR Secondary Capture Objects, CT Radiation Dose reports of patient exposure during Hybrid acquisition and Secondary Capture Objects In this situation

NM General Purpose 600/800 Series products provide the DICOM C-STORE service as a service class user (SCU).

The NM General Purpose 600/800 Series DICOM implementations support storage commitment for the already transferred data. This guarantees the user that the acquired image data are safely archived for future use. In this situation NM General Purpose 600/800 Series products provide the DICOM Storage Commitment Service as Service Class User (SCU).

The NM General Purpose 600/800 Series DICOM implementations support receiving Worklist information from a remote AE. For NM standalone exams (all products) and for hybrid exams (Discovery NM/CT 670\*, NM/CT 800\* and Optima NM/CT 640) DICOM implementations support receiving more than one Scheduled Procedure Step per study instance, enabling the acquisition of their matching number of protocols for this study. Similarly, NM General Purpose 600/800 Series products support scheduling locally more than one protocol to be performed for a study. Note that each Scheduled Procedure Step is performed independently.

For CT standalone exams the Discovery NM/CT 670\* and NM/CT 800\* DICOM implementations support scheduling a single Scheduled Procedure Step per study instance (via Worklist or by local scheduling) enabling the acquisition of a single protocol per study.

The NM General Purpose 600/800 Series DICOM implementations creates and updates Modality Performed Procedure Step (MPPS) instances managed by a remote AE in association with image acquisition. Completion or discontinuation of the MPPS is performed as the result of an operator action.

DIR DOC0660253 REV 15

The Discovery NM/CT 670\* and NM/CT 800\* DICOM Implementations use two different AE titles, one for the NM Camera and one for the CT Scanner. A verification mechanism is provided, by which a remote application entity (AE) can verify application-level communication with the NM Camera DICOM Server and CT Scanner DICOM server. Also provided is a mechanism by which a Discovery NM/CT 670\* or NM/CT 800\* user can verify application-level communication with a remote DICOM AE. In these situations, Discovery NM/CT 670\* and NM/CT 800\* provide the DICOM C-ECHO service as both a SCP and SCU, respectively, from NM and CT.

The other NM General Purpose 600/800 Series DICOM implementations provide a verification mechanism by which a remote application entity (AE) can verify application-level communication with the NM General Purpose 600/800 Series DICOM Server. Also provided is a mechanism by which a NM General Purpose 600/800 Series product user can verify application-level communication with a remote DICOM AE. In these situations, NM General Purpose 600/800 Series products provide the DICOM C-ECHO service as both a SCP and SCU, respectively.

## 2.2 IMPLEMENTATION MODEL

### 2.2.1 Overview

Two different Application Entities (AE) logically provide all DICOM functionality:

- **CT Scanner DICOM Server AE** for Discovery NM/CT 670\* and NM/CT 800\* CT standalone mode.
- **NM Camera AE** for all other modes of NM General Purpose 600/800 Series products The NM Camera AE is commanded to perform DICOM services through the use of the NM Camera user interface.

The both AEs also listen on pre-defined ports for incoming connections from remote DICOM AEs.

**Note:** In this section term **NM Camera AE DICOM implementation** refers to **DICOM implementation of all NM General Purpose 600/800 Series products, except Discovery NM/CT 670\* and NM/CT 800\* working in CT Standalone mode ( CT standalone exams scheduled via CT UI, CT images only are generated).**

**This section describes the only NM Camera Application Entity Implementation Model.**


For Discovery NM/CT 670\*, NM/CT 870 DR and NM/CT 870 CZT CT Scanner DICOM Server Application Entity Implementation model, refer to CT Scanner DICOM Conformance Statement ([Reference A](#) in Section 1.6).

For NM/CT 850 and NM/CT 860 CT Scanner DICOM Server Application Entity Implementation model, refer to CT Scanner DICOM Conformance Statement ([Reference B](#) in Section 1.6).

2.2.2 Application Data Flow Diagram

The network application model for the NM Camera DICOM Implementation is shown in the following Illustration:

ILLUSTRATION 2-1  
NM CAMERA NETWORK APPLICATION MODEL AND DATA FLOW DIAGRAM

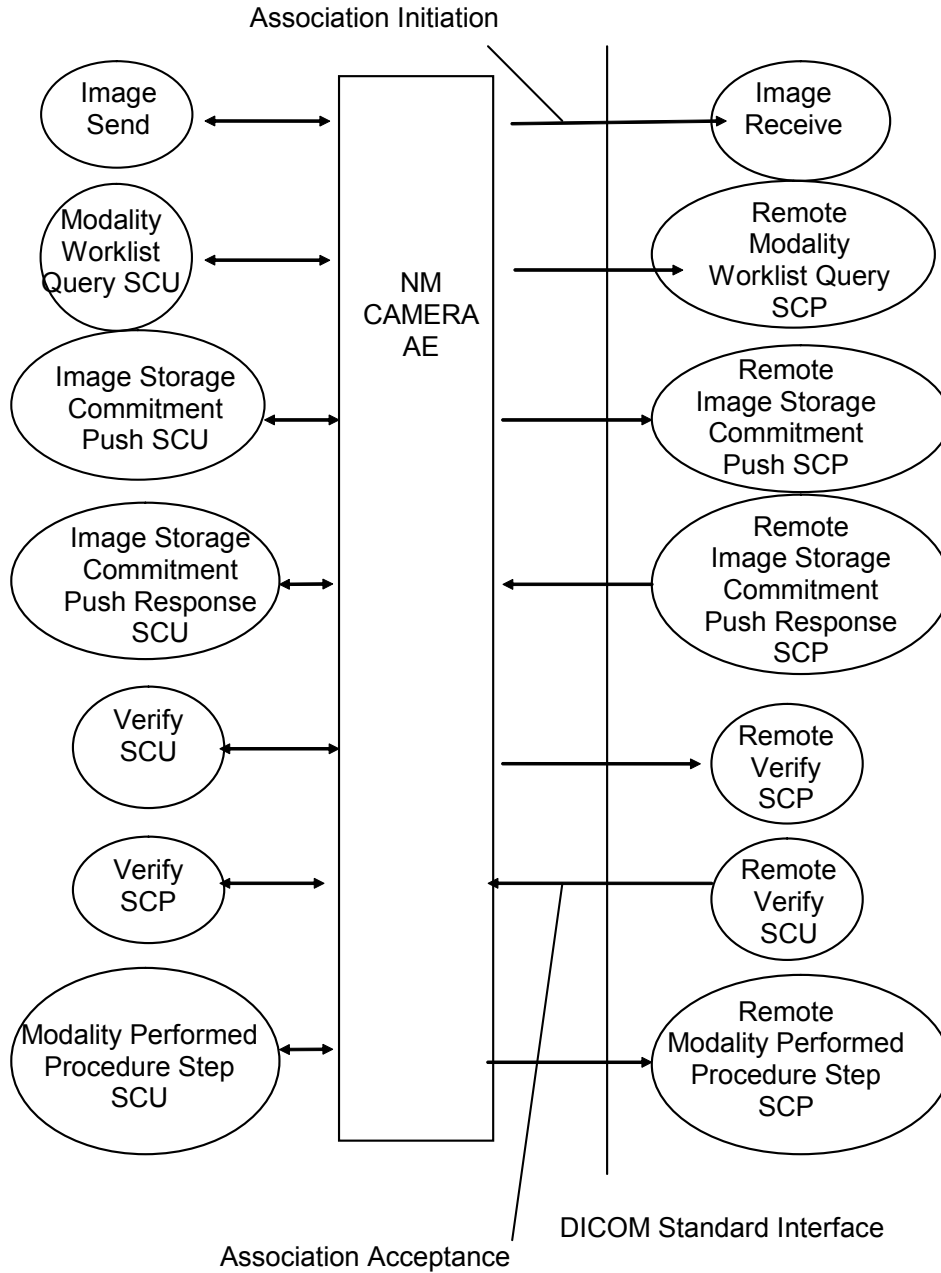

DIR DOC0660253 REV 15

### 2.2.3 Functional Definition of AE's

The NM Camera Application Entity (AE) initiates the following functions:

- *Store*: Initiates a DICOM association in order to send images to a remote AE. If the remote AE accepts a presentation context applicable to the image(s) being sent, the NM Camera AE will send the images via the C-STORE service.
- *Storage commitment*: Initiates a DICOM association in order to request a storage commitment from a remote AE. If the remote AE supports storage commitment the NM Camera AE will request a storage commitment for the image(s) previously sent successfully via the N-ACTION-RQ.
- *Verify*: Initiates a DICOM association in order to send a verification message to a remote AE via a C-ECHO-RQ message.
- *Modality Work List (MWL)*: Initiates a DICOM association in order to query the work list from a remote AE. If the remote AE accepts a presentation context applicable to the modality work list request being sent, the NM Camera AE will receive appropriate MWL responses via the C-FIND service.
- *Modality Performed Procedure Step (MPPS)*: Initiates a DICOM association in order to report the progress of the procedure step for each of the following operations:
  - *Start PPS*: in order to create a DICOM Modality Performed Procedure Step SOP instance in the remote AE when first image of the protocol is stored in database .If the remote AE accepts a presentation context applicable to Modality performed Procedure Step, the PPS DICOM SCU will issue a request to create the SOP instance in the remote AE via the N-CREATE service.
  - *Complete PPS*: in order to update a DICOM Modality Performed Step instance that is already created with the remote AE when all scan acquisitions are completed for the protocol or the operator marks this protocol as completed. If the remote AE accepts a presentation context applicable to Modality performed Procedure Step, the PPS DICOM SCU will issue a request to update the SOP instance in the remote AE via the N-SET service. The PPS Status is set to 'COMPLETED'.
  - *Discontinue PPS*: Initiates a DICOM association in order to update a DICOM Modality Performed Step instance that is already created with the remote AE when the operator marks this protocol discontinued. If the remote AE accepts a presentation context applicable to Modality performed Procedure Step, the PPS DICOM SCU will issue a request to update the SOP instance in the remote AE via the N-SET service. The PPS Status is set to 'DISCONTINUED'.

The NM Camera AE responds to the following functions:

- *Verify*: Responds to incoming C-ECHO-RQ messages by returning a C-ECHO-RSP message with a status of "Success."
- *Storage Commitment Response*: Responds to incoming N-EVENT-REPORT messages arriving from Remote AE with the status of storage commitment for images previously requested by NM Camera AE.

### 2.2.4 Sequencing of Real-World Activities

#### 2.2.4.1 Default Protocol Workflow

The NM Camera AE queries the remote station for the Modality Worklist; performs acquisition according to local schedules or by Worklist procedures; reports the start of procedure using MPPS ;sends images to auto-transfer destinations at the end of each scan of the protocol and then requests Storage Commitment for previously sent images, reports the progress of procedure using MPPS.

DIR DOC0660253 REV 15

#### 2.2.4.2 Local Scheduling

The system either does not have a Modality Worklist Server AE installed or a Modality Worklist Server AE installed but no Worklist information is obtained from HIS/RIS system for the current procedure that is being performed. The information required for performing the procedure is supplied through the user interface of the system. The use of this information in MPPS SOP instances is described later in this document (see section 7).

#### 2.2.4.3 Modality Worklist (MWL)

The system has a Modality Worklist Server AE installed. Work-List information is obtained from HIS/RIS system through the use of Basic Worklist Management Service. The user or the system initiates a MWL query (as a MWL SCU) to the MWL SCP with a given set of query parameters. The MWL SCP returns responses, which match the query parameters.

Items from the returned Worklist responses are presented to the user within “To-Do List” browser.

A subset of attributes corresponding to operator selected returned Worklist responses will be included in acquired NM, CT and CT Radiation Dose SR SOP instances related to the responses (see sections 3.4, 8.4 , 9.3 accordingly).

#### 2.2.4.4 Report of Performed Procedure Step progress

The system initiates a ‘Start PPS’ after the first acquired image is of the protocol is installed to database.

NM Camera AE initiates a MPPS N-CREATE request to the remote AE (MPPS SCP), in-order to create a MPPS SOP instance at the remote AE. The MPPS SCP returns response indicating the success/failure of the request execution. The PPS state information is updated on the system based on the response data.

When all scans of the protocol are successfully completed, system initiates ‘Complete PPS’. The remote AE returns response indicating the success/failure of the request execution. The PPS state information is updated in the system based on the response data.

The user may decide to finish the protocol even the protocol is not completed (not all scans are performed or acquisition of any scan(s) failed). In this case the system initiates ‘Complete PPS’ or ‘Discontinue PPS’ based on the choice selected by the user using the user interface provided. The remote AE returns response indicating the success/failure of the request execution. The PPS state information is updated in the system based on the response data.

#### 2.2.4.5 Send Images (Manual)

The user has to select data from Local or QC database browser on the “Data Management” panel and click on the transfer destination at the bottom the “Data Management” panel.

DIR DOC0660253 REV 15

2.3 AE SPECIFICATIONS

**Note:** For Discovery NM/CT 670\*, NM/CT 870\* CT DICOM Server AE specifications, refer to CT Scanner DICOM Conformance Statement ([Reference A](#) in Section 1.6.)

For NM/CT 850 and NM/CT 860\* CT Scanner DICOM Server Application Entity Implementation model, refer to CT Scanner DICOM Conformance Statement ([Reference B](#) in Section 1.6).

2.3.1 NM Camera AE Specification

The NM Camera Application Entity provides Standard Conformance to the following DICOM SOP Classes as an SCU and/or as an SCP:

SOP Class Name	SOP Class UID	SCU	SCP
Verification SOP Class	1.2.840.10008.1.1	Yes	Yes
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7	Yes	No
Nuclear Medicine Image Storage	1.2.840.10008.5.1.4.1.1.20	Yes	No
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	Yes (*)	No
X-Ray Radiation Dose SR	1.2.840.10008.5.1.4.1.1.88.67	Yes (*)	No
Enhanced SR	1.2.840.10008.5.1.4.1.1.88.22	Yes (***)	No
Radiopharmaceutical Radiation Dose SR	1.2.840.10008.5.1.4.1.1.88.68	Yes (**)	No
Modality Worklist Information Model - FIND	1.2.840.10008.5.1.4.31	Yes	No
Modality Performed Procedure Step SOP Class	1.2.840.10008.3.1.2.3.3	Yes	No
Storage Commitment Push Model	1.2.840.10008.1.20.1	Yes	No

**Note:** (\*) –Optima NM/CT 640 DICOM implementations only  
 (\*\*) – Except Brivo NM 615 DICOM implementations.  
 (\*\*\*) –Optima NM/CT 640 DICOM implementations only (for transfer CT Radiation Dose SR objects to systems which do not support X-Ray Radiation Dose SR SOP Class as SCP).

2.3.1.1 Association Establishment Policies

2.3.1.1.1 General

The DICOM Application Context Name (ACN), which is always proposed, is:

Application Context Name	1.2.840.10008.3.1.1.1
--------------------------	-----------------------

The maximum length PDU receive size for the NM Camera AE is:

Maximum Length PDU	64234 (Not Configurable)
--------------------	--------------------------

The SOP Class Extended Negotiation is not supported.

DIR DOC0660253 REV 15

**2.3.1.1.2 Number of Associations**

The NM Camera AE (SCU) will initiate a single DICOM association to perform a single selected object (study/series/image ) send to a remote node. One association is opened per image both in manual send and in auto-send. Multiple Send operations can be performed. The Image Storage Commitment Request (SCU) initiates a new single association for all the images that were successfully stored on the remote AE. The NM Camera AE (SCU) will initiate a single DICOM association for each PPS message sent to remote AE. NM Camera AE can initiate a maximum of 5 simultaneous associations to remote nodes.

The NM Camera AE (SCP) can accept multiple DICOM associations opened simultaneously to service verifications. The NM Camera AE can support a maximum of 5 simultaneous associations initiated by remote nodes.

**2.3.1.1.3 Asynchronous Nature**

Asynchronous mode is not supported. All operations will be performed synchronously.

**2.3.1.1.4 Implementation Identifying Information**

The Implementation UID for this DICOM Implementation is:

<b>NM Camera Implementation UID</b>	<b>1.2.840.113619.6.280</b>
<b>NM Camera Implementation Version Name</b>	<b>NUEVO_MC3_401</b>

**2.3.1.2 Association Initiation Policy**

When the NM Camera Application Entity initiates an Association for any Real-World Activity, it will propose the Presentation Contexts for all Real-World Activities; i.e., there is only a single, comprehensive Presentation Context Negotiation proposed for the AE.

The NM Camera AE proposes only a single Transfer Syntax in each Presentation Context; i.e., for each Abstract Syntax in the following Presentation Context Tables, the AE proposes one Presentation Context for each specified Transfer Syntax.

The NM Camera AE initiates a new association in the following cases:

- Due to an image send operation being initiated from the NM Camera user interface, or by auto send option.
- Due to a storage commitment request operation being initiated from the NM Camera user interface upon successful image transfer or by auto send option
- Due to a Verify operation initiated to determine whether the remote DICOM station is operational.
- Due to Modality Worklist request being initiated from the NM Camera user interface
- For every PPS operation initiated.

**2.3.1.2.1 Real-World Activity Image Send**

**Note:** For hybrid scanning in Discovery NM/CT 670\* and NM/CT 800\* DICOM implementations selecting CT image(s) to be transferred and selecting the send destination is done from the CT scanner. For CT Image Send (Push) Real World activity description please refer to Discovery NM/CT 670\* and NM//CT 870\* CT Scanner DICOM Conformance Statement ([Reference A](#) in Section 1.6) and NM/CT 850 /NM/CT 860 CT Scanner DICOM Conformance Statement ([Reference B](#) in Section 1.6).

**Note:** For all other NM General Purpose 600/800 Series products sending acquired data (manually or automatically) is performed via the NM Camera Operator Console screen.



DIR DOC0660253 REV 15

**2.3.1.2.1.1 Associated Real-World Activity**

There are two ways to send NM data: manual and automatic. In the manual way, in order to send data, the operator must both select objects (Study/Series/ Image(s)) to be transferred from the Data Management and select a destination by pressing the destination button.

Once these selections have been made, the operator pushes the Transfer Destination button to initiate a send operation. In the automatic way, when any single acquisition (scan) has been completed (in a NM-only context as well as in the hybrid context) the data is automatically sent to the destination(s) which are predefined within system configuration; the NM Camera AE will then initiate an association with the remote AE in order to send the selected object(s) – one object per association – and will accept and interpret responses received from the remote AE.

The UI will indicate the status of the object (Study, Series, Image) being transferred. The status can be one of IN-PROGRESS, SUCCESS, or FAILURE. The associated error messages due to a failed status can be found in system log.

**2.3.1.2.1.2 Proposed Presentation Context Table**

Presentation Context Table – Proposed by NM Camera AE for “Image Send “Activity					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Nuclear Medicine Image Storage	1.2.840.10008.5.1.4.1.1.20	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
Secondary Image Capture Storage	1.2.840.10008.5.1.4.1.1.7	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
CT Image Storage(*)	1.2.840.10008.5.1.4.1.1.2	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
X-Ray Radiation Dose SR (*)	1.2.840.10008.5.1.4.1.1.88.67	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
Enhanced SR (***)	1.2.840.10008.5.1.4.1.1.88.22	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
Radiopharmaceutical Radiation Dose SR (***)	1.2.840.10008.5.1.4.1.1.88.68	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		

**Note:** (\*) –Optima NM/CT 640 DICOM implementation only  
 (\*\*) – Except Brivo NM 615 DICOM implementations.  
 (\*\*\*) –Optima NM/CT 640 DICOM implementation only (for transfer CT Radiation Dose SR objects to systems which do not support X-Ray Radiation Dose SR SOP Class as SCP).

**2.3.1.2.1.2.1 SOP Specific DICOM Conformance Statement for All Storage SOP Classes**

The NM Camera AE implementation includes optional data elements in the SOP Instances as described in Sections 3, 4 and 8.

This implementation can perform a single C-STORE operation over a single association.

All the operations used by this SOP class support an **Association Timer**. This timer is started when the association request is sent, and is stopped when the respective response is received. The default time-out value is 60 seconds and is not configurable.

DIR DOC0660253 REV 15

Each C-STORE operation also supports an **Operation Inactivity Timer**. This time-out starts once the first C-STORE request has been issued (on association) or received and is reset each time a C-STORE response has been received or when subsequent C-STores are sent. This time-out is 5 minutes. It is non-configurable.

If any of the two timers mentioned above expires, the connection is closed and the operation in progress is considered failed.

Upon receiving a C-STORE confirmation containing a status other than Successful or Warning, this implementation will consider the current request to be a failure but will continue to attempt to send any remaining images in the request on a different association.

When NM Camera AE initiates an association to issue a C-STORE, the image will be transmitted by the NM Camera AE with the same elements as was created locally.

Transfer log shows one of these statuses for store request: JOB\_SUCCEEDED, JOB\_FAILED, JOB\_IN\_PROGRESS. The specific error codes can be observed in the log.

Following are the status codes that are more specifically processed when receiving messages from a **Storage SCP** equipment :

Service Status	Status Code	Further Meaning	Application Behavior When Receiving Status Code
Failure	A700-A7FF	Refused: Out of resources	The message " JOB FAILED: error sending image" is displayed in Transfer Log.
	A900-A9FF	Error: Data Set does not match SOP Class	The message " JOB FAILED: error sending image" is displayed in Transfer Log
	C000-CFFF	Error: Cannot Understand	The message " JOB FAILED: error sending image" is displayed in Transfer Log
	0122	SOP Class Not Supported	The message " JOB FAILED: error sending image" is displayed in Transfer Log
Warning	B000	Coercion of Data Elements	The message "JOB SUCCEEDED" posted to the Transfer Log.
	B006	Elements Discarded	The message "JOB SUCCEEDED" posted to the Transfer Log.
	B007	Data Set does not match SOP Class	The message "JOB SUCCEEDED" posted to the Transfer Log.
Success	0000		The message "JOB SUCCEEDED" posted to the Transfer Log.
*	*	Any other status code.	The message " JOB FAILED: error sending image" is displayed in Transfer Log.

**2.3.1.2.2 Real-World Activity Modality Worklist Query SCU**

**Note:** For NM standalone exams and NM/CT hybrid exams the Worklist query is performed via the NM Camera Operator Console screen. Worklist for CT standalone scanning is performed via CT Scanner Operator Console. For Real World activity: Worklist Query for CT standalone please refer to Discovery NM/CT 670\* and NM/CT 870\* CT Scanner DICOM Conformance Statement

DIR DOC0660253 REV 15

([Reference A](#) in Section 1.6) and to NM/CT 850 and NM/CT 860 CT Scanner DICOM Conformance Statement ([Reference B](#) in Section 1.6).

**2.3.1.2.2.1 Associated Real-World Activity**

NM Camera AE queries the remote AE for a Modality Worklist (MWL) in the following cases:

- When NM Camera application is started, MWL query is automatically performed for updating entries displayed in the NM Camera “To Do List”.
- User opens “Filter...” button in NM Camera UI. MWL Filter dialog is opened, user defines relevant MWL matching keys and presses on “Query” button.
- Users requires MWL query using latest MWL matching keys defined by pressing on “Refresh” button in the NM camera UI.
- Users requires MWL query using latest MWL matching keys defined to map of the MWL Requested or Scheduled Procedures to the NM Camera acquisition protocols .

**2.3.1.2.2.2 Proposed Presentation Context Table**

Presentation Context Table – Proposed by NM Camera AE for Modality Worklist Query SCU Activity					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Modality Worklist Information Model - FIND	1.2.840.10008.5.1.4.31	Implicit VR Little Endian	1.2.840.10008.1.2	SCU	None
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

**2.3.1.2.2.2.1 SOP Specific DICOM Conformance Statement for the Modality Worklist Information Model - FIND SOP Class**

The NM Camera includes matching keys in the Modality Worklist queries as described in Section 5

If Modality Worklist query failed, the user receives a notification message.

All the operations used by this SOP class support an **Association Timer**. This timer is started when the association request is sent, and is stopped when the respective response is received. The default time-out value is 60 seconds, and is not configurable.

All the operations used by this SOP class support a **Session Timer**. This timer is started when the association is established, and stopped when association is ended. The default time-out value is 180 seconds.

If any of the two timers mentioned above expires, the connection is closed and the operation in progress is considered failed

Sending C-FIND CANCEL is not supported by the NM Camera AE.

If NM Camera AE receives a success response with no matching response, a pop-up will show saying that the Worklist query has failed.

When the C-FIND response received from the Worklist SCP does not include one of the tags defined as type 1 (e.g. Patient name, Patient ID, Study Instance UID) the NM Camera AE will show a pop-up message listing which tags are missing, the Worklist item will be rejected.

DIR DOC0660253 REV 15

Following are the status codes that are more specifically processed when receiving messages from a Modality Worklist SCP equipment :

Service Status	Status Code	Further Meaning	Application Behavior When Receiving Status Code
Failure	A700	Refused: Out of resources	User receives a notification message. Only locally scheduled studies or studies which acquisition is in progress are displayed in the NM Camera "To Do List"
	A900	Error: Identifier does not match SOP Class	User receives a notification message. Only locally scheduled studies or studies which acquisition is in progress are displayed in the NM Camera "To Do List"
	C000-CFFF	Error: Unable to process	User receives a notification message. Only locally scheduled studies or studies which acquisition is in progress are displayed in the NM Camera "To Do List"
	0122	SOP Class Not Supported	User receives a notification message. Only locally scheduled studies or studies which acquisition is in progress are displayed in the NM Camera "To Do List"
Cancel	FE00	Matching terminated due to cancel	Not Applicable
Success	0000	Matching is complete - No final identifier is supplied	All accepted MWL entries, received from MWL SCP system, along with locally scheduled studies or studies which acquisition is in progress are displayed in the NM Camera "To Do" List
Pending	FF00	Matches are continuing - Current Match is supplied and any Optional Keys were supported in the same manner as Required Keys.	"MWL query is in progress" notification is displayed along with MWL query progress indicator.
*	FF01	Matches are continuing - Warning that one or more Optional Keys were not supported for existence for this Identifier	"MWL query is in progress" notification is displayed along with MWL query progress indicator.
	*	Any other status code.	User receives a notification message. Only locally scheduled studies or studies which acquisition is in progress are displayed in the NM Camera "To Do List"

**2.3.1.2.3 Real-World Activity Image Storage Commitment Push SCU**

**Note:** For Real World activity Image Storage Commitment Push SCU of CT Scanner AE please refer to please refer to Discovery NM/CT 670\* and NM/CT 870\* CT Scanner DICOM Conformance Statement ([Reference A](#) in Section 1.6) and to

DIR DOC0660253 REV 15

NM/CT 850 and NM/CT 860 CT Scanner DICOM Conformance Statement  
 ([Reference B](#) in Section 1.6) Associated Real-World Activity.

The operator must both select objects (Study/Series/Image(s)) to be transferred from the Data Transfer panel, and select a destination from the list of previously defined destinations. Once these selections have been made, the operator pushes the "Destination" button to initiate an image send operation.

If the destination is configured as storage commitment capable or the destination is configured to use other storage commitment capable devices, the NM Camera AE initiates the following operations:

- Negotiates and establishes association with remote Storage Commitment Provider.
- Sends the selected images to the remote DICOM AE.
- Closes the association.
- If all the images are transferred without failures the following steps will be executed. If there are any failures the job will be marked as failed and the Storage Commitment request will not be sent.
- Establishes a new association for sending the commitment request for all successfully transferred images.
- Receives the response on same association or on a different association.
- Updates the archive flag information for successful instances.

The Transfer Log shows the status of the storage commitment request progress. The status can be either JOB\_SUCCEEDED, JOB\_FAILED, or JOB\_IN\_PROGRESS.

**2.3.1.2.3.1 Proposed Presentation Context Table**

Presentation Context Table – Proposed by NM Camera AE for Image Storage Commitment Push SCU Activity					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Storage Commitment Push Model	1.2.840.10008.1.20.1	Implicit VR Little Endian	1.2.840.10008.1.2	SCU	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		

**2.3.1.2.3.1.1 SOP Specific DICOM Conformance Statement for the Storage Commitment Push Model SOP Class SCU**

The storage commitment request (N-ACTION) can perform a storage commitment request for multiple images over a single association. A new association is initiated for the storage commitment request for every group of successfully transferred image(s).

Upon receiving a N-ACTION confirmation containing a “Successful” status, the next N-ACTION-RQ operation is performed for the new association.

Upon receiving a N-ACTION confirmation containing a “Refused” status, the association is terminated. The reason for termination is recorded in the system log file.

DIR DOC0660253 REV 15

Upon receiving a N-ACTION confirmation containing a status other than the DICOM standard defined values, the current request is considered to be a failure and will terminate the association. The reason for termination is recorded in the system log file.

The NM Camera AE uses DICOM network storage services to transfer SOP Instances which are to be committed.

The NM Camera AE may request Storage Commitment for Instances of any of the Composite SOP Classes it supports as an SCU (see Section 2.3.1).

The Storage Commitment Information Object is described in Section 6.

All the operations used by this SOP class support an **Association Timer**. This timer is started when the association request is sent, and is stopped when the respective response is received. The default time-out value is 60 seconds, and is not configurable.

All the operations used by this SOP class support a **Session Timer**. This timer is started when the association is established, and stopped when the association is ended. The default time-out value is 10 minutes and is not configurable.

If any of the two timers mentioned above expires, the connection is closed and the operation in progress is considered failed.

Following are the status codes that are more specifically processed when receiving N-Action responses from a **Storage Commitment** SCP equipment:

Service Status	Status Code	Further Meaning	Application Behavior When Receiving Status Code
Failure	0119	Class-instance conflict	The message " JOB FAILED: error storage commitment" is displayed in Transfer Log; Storage commitment will be considered as failed;
	0112	No such SOP Instance	The message " JOB FAILED: error storage commitment" is displayed in Transfer Log; Storage commitment will be considered as failed;
	0110	Processing failure	The message " JOB FAILED: error storage commitment" is displayed in Transfer Log; Storage commitment will be considered as failed;
	0213	Resource limitation	The message " JOB FAILED: error storage commitment" is displayed in Transfer Log; Storage commitment will be considered as failed;
	0122	Referenced SOP Class Not Supported	The message " JOB FAILED: error storage commitment" is displayed in Transfer Log; Storage commitment will be considered as failed;
	0131	Duplicate Transaction UID	The message " JOB FAILED: error storage commitment" is displayed in Transfer Log; Storage commitment will be considered as failed;

DIR DOC0660253 REV 15

Success	0000		The message " JOB IN PROGRESS: waiting for commitment" is displayed in Transfer Log. The request for storage comment is considered successfully sent. A timer is started which will expire if no N-EVENT-REPORT for the Transaction UID is received within a configurable timeout period.
*	*	Any other status code.	The message " JOB FAILED: error storage commitment" is displayed in Transfer Log; Storage commitment will be considered as failed;

As part of the storage commitment implementation, Remote AE (SCP) will initiate an association to this implementation and will send an N-EVENT-REPORT. The attribute of the N-EVENT-REPORT message will include an indication on all images for which a commitment has succeeded and those for which it has failed.

The receipt of a N-EVENT-REPORT on an association that NM Camera has initiated is not supported. The Remote AE (SCP) must initiate a new association in order to send the new N-EVENT-REPORT (see Section 2.3.1.3.2.2.1).

**2.3.1.2.4 Real-World Activity Verify SCU**

**Note:** For CT scanner AE Real Time Activity: Verify refer to Discovery NM/CT 670\* and NM/CT 870\* CT Scanner DICOM Conformance Statement ([Reference A](#) in Section 1.6) and NM/CT 850 and NM/CT 860\* CT Scanner DICOM Conformance Statement ([Reference B](#) in Section 1.6).

**2.3.1.2.4.1 Associated Real-World Activity**

Service personnel invoke the DICOM Station Configuration Utility from NM Camera UI. The user selects any of defined remote DICOM AE and presses on "Refresh status" button. The NM Camera AE initiates an association with the remote DICOM AE in order to verify communication at the application level. The success or failure of the verification process is displayed to the user.

**2.3.1.2.4.2 Proposed Presentation Context Table**

Presentation Context Table – Proposed by NM Camera AE for Verify SCU Activity					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Verification SOP Class	1.2.840.10008.1.1	Implicit VR Little Endian	1.2.840.10008.1.2	SCU	None

**2.3.1.2.4.2.1 SOP Specific DICOM Conformance Statement for Verification SOP Class**

The NM Camera AE provides standard conformance to the DICOM Verification Service Class. All the operations used by this SOP class support an Association Timer. This timer is started when the association request is sent, and is stopped when the respective response is received. The default time-out value is 15 seconds and is not configurable.

**2.3.1.2.5 Real-World Activity Performed Procedure Step creation and update**

DIR DOC0660253 REV 15

**2.3.1.2.5.1 Associated Real-World Activity**

The real-world activities are mentioned in section 2.2.4.4 describing PPS Sequencing of Real-World Activities. Each of the real world activity results in either creating a new Performed procedure Step SOP instance at the remote SCP or updating an already created Performed Procedure Step SOP instance as per the DICOM standard.

NM Camera AE generates an association establishment request upon the following events:

- In order to issue an N-CREATE message with Performed Procedure Step status = IN-PROGRESS when the first image of the first scan in the protocol is stored in database.
- In order to send a final intermediate N-SET message with Performed Procedure Step status =COMPLETED in the following events:
  - When the acquisition of the last scan in a protocol is completed leaving no non-acquired scans in the protocol, the protocol is automatically considered completed.
  - When the operator presses the "Protocol Completed" button for a selected protocol.
  - When the operator presses the "Protocol Completed" button where a study instance is selected, the N-SET message is sent for all protocols in a study that has acquired data.
- In order to send a final N-SET message with Performed Procedure Step status = DISCONTINUED when the operator presses the "Discontinue Protocol" button for a selected protocol or for all of the protocols in a study when study is selected. The message will be set with the reason for discontinuing the protocol as selected by the operator (see Table 7-5 for list of possible discontinue reasons).

**2.3.1.2.5.2 Proposed Presentation Context Table**

Presentation Context Table – Proposed by NM Camera AE for Performed Procedure Step creation and update Activity					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Modality Performed Procedure Step	1.2.840.10008.3.1.2.3.3	Implicit VR Little Endian	1.2.840.10008.1.2	SCU	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		

**2.3.1.2.5.2.1 SOP Specific DICOM Conformance Statement for Modality Performed Procedure Step SOP Class**

When sending a PPS message, the system looks for stations configured to act as PPS manager.

If the remote AE does not support the proposed Presentation context, an appropriate error message is logged. Only one association is opened at a time.

The NM Camera AE includes attributes in the Modality Performed Procedure Step N-CREATE and N-SET messages as described in Section 7.4.

NM Camera AE does not define any extended error codes. The standard error codes are handled. On a response with status ‘Success’ for the N-CREATE or N-SET request, the system updates the state and indicates the same on the user interface. On a response with status other than ‘Success’ the operation is deemed ‘Failed’.

If the operation is ‘Failed’, detailed message is logged into system log-file and system automatically retries the failed operation.

If the N-CREATE operation fails, the system marks the condition to enable sending N-CREATE again upon next scan's acquisition start.

This is an update that will be reported again by sending of the final N-SET COMPLETE or DISCONTINUE message.

The PPS N-SET request can inform of completion of multiple images/SR objects over a single association.

**2.3.1.3 Association Acceptance Policy**



DIR DOC0660253 REV 15

The NM Camera AE places no limitation on whom may connect to it. The maximum number of associations accepted in parallel is limited to 5.

Any remote AE can open an association to the NM Camera AE for the purpose of application level communication verification.

As part of the storage commitment implementation, the NM Camera AE Server responds to N-EVENT-REPORT received from the remote AE.

**2.3.1.3.1 Real-World Activity Verify SCP**

**2.3.1.3.1.1 Associated Real-World Activity**

The NM Camera AE is always listening for associations. No operator action is required to respond to a Verification request.

The real-world activity associated with the Verification request is to send a C-ECHO-RSP message with a status of “Success” to the requesting AE.

**2.3.1.3.1.2 Accepted Presentation Context Table**

Presentation Context Table - Accepted by NM Camera AE for Verify SCP Activity					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Verification SOP Class	1.2.840.10008.1.1	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		

**2.3.1.3.1.2.1 SOP Specific DICOM Conformance Statement for Verification SOP Class**

The NM Camera AE provides standard conformance to the DICOM verification service class.

**2.3.1.3.2 Real-World Activity Image Storage Commitment Push Response SCU**

**2.3.1.3.2.1 Associated Real-World Activity**

As part of the storage commitment implementation, Remote AE (SCP) initiates an association to this implementation and sends an N-EVENT-REPORT. The attribute of the N-EVENT-REPORT message includes an indication on all images for which a commitment has succeeded and those for which it has failed.

The receipt of an N-EVENT-REPORT on an association that NM Camera AE has initiated is not supported. The Remote AE (SCP) must initiate a new association in order to send the new N-EVENT-REPORT.

On reception of a successful N-EVENT-REPORT-RQ notification from the Storage Commitment Provider, the images are flagged as committed in the database.

**2.3.1.3.2.2 Accepted Presentation Context Table**

Presentation Context Table - Accepted by NM Camera AE for Image Storage Commitment Push Response SCU Activity					
Abstract Syntax		Transfer Syntax		Role	Extended

DIR DOC0660253 REV 15

Name	UID	Name List	UID List		Negotiation
Storage Commitment Push Model	1.2.840.10008.1.20.1	Implicit VR Little Endian	1.2.840.10008.1.2	SCU	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		

**2.3.1.3.2.2.1 SOP Specific DICOM Conformance Statement for the Storage Commitment Push Model SOP Class SCU**

The NM Camera AE will only accept the SCU role (which must be proposed via SCP/SCU Role Selection Negotiation) within a Presentation Context for the Storage Commitment Push Model SOP Class.

Upon receiving a Storage Commitment N-EVENT-REPORT (Storage Commitment Result), the NM Camera AE will validate the Transaction UID against its list of outstanding Storage Commitment Request Transaction UIDs. If it matches an outstanding Request, the AE will mark all SOP Instances for which a success status is indicated with an Archived flag, shown on the user interface as “Archive” Icon. User shall explicitly request deleting of Instances , which are not marked as “Archived”

If the Storage Commitment Result indicates any failure status the error will be written to the **/home/ctuser/neuvo/logfiles/nwscp.log** error log.

Any retry must be manually reinitiated as a new Storage request following by Storage Commitment Request (see Section 2.3.1.2.3). The list of specific Failure Reason Codes that this AE will be able to process is described in Section 6.1.2.1

Following are the status codes the Application may send back in the N-EVENT-REPORT response command to the Storage Commitment SCP Equipment that sent the N-EVENT-REPORT request:

Service Status	Status Code	Further Meaning	Status Code Explanation	Related Fields Sent Back to the SCU
Failure	0110	Processing failure	This status code is sent if SCANNER AE failed to understand N_EVENT_RESPONSE message	None
Success	0000		Send in case that SCANNER AE successfully processed N_EVENT_RESPONSE message	None

**2.3.1.3.2.3 Presentation Context Acceptance Criterion**

The NM Camera AE evaluates each Presentation Context independently, and accepts any Presentation Context that matches an Abstract Syntax for any Real-World Activity.

**2.3.1.3.2.4 Transfer Syntax Selection Policies**

Within each Presentation Context, the NM Camera AE will select Transfer Syntaxes according to the following priority (highest priority first):

1. Explicit VR Little Endian
2. Implicit VR Little Endian

DIR DOC0660253 REV 15

## 2.4 COMMUNICATION PROFILES

### 2.4.1 Supported Communication Stacks

The DICOM Upper Layer Protocol is supported using TCP/IP, as specified in DICOM PS3.8.

The TCP/IP stack is inherited from the Linux Operating System.

### 2.4.2 Physical Media Support

Ethernet 802.3 provides the physical network layer for this product.

## 2.5 EXTENSIONS / SPECIALIZATIONS/ PRIVATIZATIONS

### 2.5.1 Standard Extended / Specialized / Private SOP Classes

#### 2.5.1.1 Standard Extended SOP Classes

The product provides Standard Extended Conformance to all supported SOP Classes, through the inclusion of additional Type 3 Standard Elements and Private Data Elements. The extensions are defined in Sections 3.5, 4.5 , 8.5, 9.4 and 10.4.

#### 2.5.1.2 Private SOP Classes

The NM General Purpose 600/800 Series systems do not implement any Private SOP Class.

### 2.5.2 Private Transfer Syntaxes

The NM General Purpose 600/800 Series systems do not implement any Private transfer syntax.

### 2.5.3 Additional Protocols

The NM General Purpose 600/800 Series systems support DHCP Protocol.

### 2.5.4 IPv4 and IPv6 Support

The NM General Purpose 600/800 Series systems support IPv4 only.

## 2.6 CONFIGURATION

The NM General Purpose 600/800 Series systems are configured by GE Healthcare Field Service Engineers. The DICOM configuration items below are configurable or re-configurable by a Field Service Engineer but are not accessible through the implementation user interface.

### 2.6.1 AE Title/Presentation Address Mapping

The NM General Purpose 600/800 Series systems allow for the configuration of the mapping of remote AE titles to IP addresses and ports. The IP address of a remote AE may be in a different subnet (using routing). A router is configurable to ensure communication from one sub-net to another. This configuration is performed by GE Healthcare Field Service Engineers.

### 2.6.2 Configurable Parameters

The following fields are configurable for NM Camera AE (local):

- Local AE Title

DIR DOC0660253 REV 15

**Note:** Listening Port Number (default value is 4006), PDU length and any time-outs are not configurable for NM Camera AE. The configuration of IP routers and subnet mask is available on a OS level.

The following fields are configurable for every remote DICOM AE:

- Remote AE Title
- Remote IP Address
- Listening TCP/IP Port Number
- Remote AE functionality flags:
  - PPS provider
  - Send destination
  - Auto-send destination
  - Modality Worklist Provider
  - Storage Commit Server only
  - Storage Commitment on (AE Title of one of previously defined Storage Commitment Servers)
  - Auto-Processing destination (shall be used for Xeleris Workstations only)
  - Auto-Reconstruction destination (shall be used for Smart Console© systems only)

**Note:** All configurations must be performed by a GE Field Engineer. The DICOM configuration items are configurable or re-configurable by a Field Service Engineer but are not accessible through the NM Camera user interface.

## 2.7 SUPPORT OF EXTENDED CHARACTER SETS

The NM Camera AE supports only a single single-byte extended character set ISO\_IR 100 (Latin alphabet Number 1 supplementary set).

The NM Camera user interface will allow the user to enter characters from the console keyboard that is within ASCII or the configured extended character set. The product specifies ISO\_IR 100 (Latin alphabet Number 1) extended character set in Specific Character Set (0008,0005) whether any such extended characters are included in SOP Instances or not.

The product will accept, as a Modality Worklist SCU, Scheduled Procedure Step Identifiers with any value of Specific Character Set (0008,0005) defined by DICOM standard. It will map that Specific Character Set value without change into the images created pursuant to that Scheduled Procedure Step. Text attributes of the Scheduled Procedure Step Identifier, including Patient and Physician names that include extended characters will be displayed as described above. If Modality Worklist entries do not contain Specific Character Set (0008,0005) value, the NM Camera AE adds ISO\_IR 100 (Latin alphabet Number 1) extended character set in Specific Character Set (0008,0005) to the images created pursuant to that Scheduled Procedure Step.

## 2.8 CODES AND CONTROLLED TERMINOLOGY

### 2.8.1 Fixed Coded Terminology

The NM Camera AE uses the fixed (non-configurable, non-extensible) coded terminology in NM Image SOP Instance as described in Section 3 and in Modality Performed Procedure attributes as described in Section 7

The NM Camera DICOM implementation is capable of supporting arbitrary coding schemes for Procedure and Protocol Codes. During installation, a service technician will establish a mapping between the site-specific codes and the Protocol Names used internally to identify acquisition protocols. A remote AE station configured to act as Worklist provider is configured to map according to one of the DICOM tags:

- (0040,0007) - Scheduled Procedure Step Description

DIR DOC0660253 REV 15

- (0032,1060) - Requested Procedure Description
- (0040,0008) - Scheduled Protocol Code Sequence

NM Camera Implementation extends the following sets of coded terminology:

Context Group	Default Value Set	Use
Radionuclides	CID 18	Selected value from this group is used in NM Image Radiopharmaceutical Information Sequence (0054,0016) > Radionuclide Code Sequence (0054,0300)  (GEHC-NUC, 20001, ^166^Holmium)
Radiopharmaceuticals	CID 25	Selected value from this group is used in NM Image Radiopharmaceutical Information Sequence (0054,0016) > Radiopharmaceutical Code Sequence (0054,0304)  (SRT,C-B10C4, Technetium^99m Ethyl cysteinate dimer) (SRT,C-B1184, Strontium-89 chloride) (SRT,C-B10A5, Yttrium-90 microspheres) (GEHC-NUC,10001, Iodine^123^ DaTSCAN Ioflupane) (GEHC-NUC,10004, Lutetium^177) (GEHC-NUC,10005, Xeneisol Xe^133^) (GEHC-NUC,10006, Technetium^99m TRODAT-1) (GEHC-NUC,10007, Radium Ra 223 dichlorid) (GEHC-NUC,10008, Lutetium^177 DOTATATE) (GEHC-NUC,10009, Lutetium^177 NAAG n-acetylaspartylglutamate)

**2.8.2 Mapped Coded Terminology**

The product does not use any mapped coded terminology

**2.8.3 Configurable Coded Terminology**

The product does not use any configurable coded terminology

**2.9 SECURITY PROFILES**

The product does not conform to any defined DICOM Security Profiles.

It is assumed that the product is used within a secured environment. It is assumed that a secured environment includes at a minimum:

1. Firewall or router protections to ensure that only approved external hosts have network access to the product.
2. Firewall or router protections to ensure that the product only has network access to approved external hosts and services.
3. Any communications with external hosts and services outside the locally secured environment use appropriate secure network channels (such as a Virtual Private Network (VPN))

### 3. NM INFORMATION OBJECT IMPLEMENTATION

#### 3.1 INTRODUCTION

**Note:** In this section term NM Camera AE DICOM implementation refers to DICOM implementation of NM General Purpose 600/800 Series products except Discovery NM/CT 670\* and NM/CT 800\* in CT standalone mode.

This section specifies the use of the DICOM NM Image IOD to represent the information included in NM Images produced by this implementation. Corresponding attributes are conveyed using the module construct.

#### 3.2 NM CAMERA MAPPING OF DICOM ENTITIES

The NM Camera AE maps DICOM Information Entities to local Information Entities in the product’s database and user interface.

**TABLE 3-1  
MAPPING OF DICOM ENTITIES TO NM CAMERA ENTITIES**

DICOM IE	NM Camera Entity
Patient	Patient
Study	Study
Series	Series
Image	Dataset

#### 3.3 IOD MODULE TABLE

The Nuclear Medicine Information Object Definition comprises the modules of the following table, plus Standard Extended and Private attributes. Standard Extended and Private Attributes are described in Section 3.5.

**TABLE 3-2  
NM IMAGE IOD MODULES**

Entity Name	Module Name	Usage	Reference
Patient	Patient	Used	3.4.1.1
	Clinical Trial Subject	Not Used	N/A
Study	General Study	Used	3.4.2.1
	Patient Study	Used	3.4.2.2
	Standard Extended Study	Used	3.4.2.3
	Private Study	Used	3.4.2.4
	Clinical Trial Study	Not Used	N/A
Series	General Series	Used	3.4.3.1

DIR DOC0660253 REV 15

	Clinical Trial Series	Not Used	N/A
	Standard Extended Series	Used	3.4.3.2
	NM/PET Patient Orientation	Used	3.4.3.3
	Private Series	Used	3.4.3.4
Frame of Reference	Frame of Reference	Used for images where Image Type (0008,0008) Value 3 is TOMO or GATED TOMO	3.4.4.1
	Synchronization	Not Used	N/A
Equipment	General Equipment	Used	3.4.5.1
Image	General Image	Used	3.4.6.1
	Image Pixel	Used	3.4.6.2
	Acquisition Context	Used for Cardiac SPECT images only	3.4.6.3
	Device	Not Used	N/A
	NM Image Pixel	Used	3.4.6.4
	Specimen	Not Used	N/A
	Multi-frame	Used	3.4.6.5
	NM Multi-frame	Used	3.4.6.6
	NM Image	Used	3.4.6.7
	NM Isotope	Used	3.4.6.8
	NM Detector	Used	3.4.6.9
	NM Tomo-Acquisition	Used for images where Image Type (0008,0008) Value 3 is TOMO or GATED TOMO	3.4.6.10
	NM Multi-gated Acquisition	Used for images where Image Type (0008,0008) Value 3 is GATED or GATED TOMO	3.4.6.11
	NM Phase	Used for images where Image Type (0008,0008) Value 3 is DYNAMIC	3.4.6.12
	NM Reconstruction	Not Used	N/A
	Overlay Plane	Not Used	N/A
	Multi-frame Overlay	Not Used	N/A
	VOI LUT	Used	3.4.6.13
	SOP Common	Used	3.4.6.14
	Frame Extraction	Not Used	N/A
Private Image	Used	3.4.6.15	
Private NM Image	Used	3.4.6.16	

Private Image Tomo	Used for images where Image Type (0008,0008) Value 3 is TOMO or GATED TOMO	3.4.6.17
Private Image Multi-Gated	Used for images where Image Type (0008,0008) Value 3 is GATED or GATED TOMO	3.4.6.18
Private Image GSPECT	Used for images where Image Type (0008,0008) Value 3 is GATED TOMO	3.4.6.19

**3.4 INFORMATION MODULE DEFINITIONS**

Please refer to DICOM Part 3 (Information Object Definitions) for a description of each of the entities, modules, and attributes contained within the NM Information Object.

The following modules are included to convey Enumerated Values, Defined Terms, and Optional Attributes supported Type 1 & Type 2 Attributes are also included for completeness and to define what values they may take and where these values are obtained from when generating the instance. It should be noted that they are the same ones as defined in the DICOM Standard Part 3 (Information Object Definitions). Also note that Attributes which are not present in tables are not supported.

NM Camera private attributes are defined in private modules, each of which follows the related Standard module. Private data element tags are assigned following the rules given in Part 5 of the DICOM v3.0 Standard, and are identified using the (gggg, xxe) format, where xx represents a reserved block of element numbers within the group gggg.

Note that any element not listed in table(s) means that it is not supported (not stored in the created images).

**3.4.1 Patient Entity Modules**

**3.4.1.1 Patient Module**

**TABLE 3-3  
PATIENT MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Patient's Name	(0010,0010)	2	Patient's full name. (*) Compound from user Last Name and First Name for locally scheduled protocols(**) (***)
Patient ID	(0010,0020)	2	Primary hospital identification number or code for the patient. (*) (**) (***)
Issuer of Patient ID	(0010,0021)	3	Not Used
Issuer of Patient ID Qualifiers Sequence	(0010,0024)	3	Not Used
Patient's Birth Date	(0010,0030)	2	Birth date of the patient. (*) (**) (***) Sent as ZERO LENGTH value if value is not received either from MWL or from user input.



Patient's Sex	(0010,0040)	2	Sex of the named patient. (*) (**)(***) Enumerated Values: M = male F = female O = other Sent as ZERO LENGTH value if value is not received either from MWL or from user input.
Other Patient IDs	(0010,1000)	3	Not Used
Other Patient IDs Sequence	(0010,1002)	3	Not Used

**Note 1 :** (\*) - Attributes copied from the Worklist if the study source was actually copied from a Worklist query result.

**Note 2 :** (\*\*) - Attributes copied from the user input for Locally scheduled Protocols

**Note 3 :** (\*\*\*) - Cannot be modified by user if meaningful value is received from MWL

3.4.2 Study Entity Modules

3.4.2.1 General Study Module

TABLE 3-4  
GENERAL STUDY MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Study Instance UID	(0020,000D)	1	Unique identifier for the Study.(*)(***) Generated by the system for Locally Scheduled protocols
Study Date	(0008,0020)	2	Date the Study started. Set to start date of the first protocol in the study.
Study Time	(0008,0030)	2	Time the Study started. Set to start time of the first protocol in the study.
Accession Number	(0008,0050)	2	A RIS generated number that identifies the order for the Study. (*) (**) Sent as ZERO LENGTH value if value is not received either from MWL or from user input. User can modify value received from MWL
Referring Physician's Name	(0008,0090)	2	Name of the patient's referring physician (*) (**) Only Last and First Names received from MWL are displayed in UI Sent as ZERO LENGTH value if value is not received either from MWL or from user input. User can modify value received from MWL

Study ID	(0020,0010)	2	User or equipment generated Study identifier. Automatically assigned to the short name of the first Protocol in the study (**) May be updated by user Always sent as non-empty value.
Study Description	(0008,1030)	3	Study Description. (*) (**) Automatically assigned to the short name of the first Protocol in the study for the Locally Scheduled Protocols. Copied from Requested Procedure Description (0032,1060) of the first SPS in the study when read from Worklist Sent as ZERO LENGTH value if value is not received neither from MWL nor from user input User can modify value received from MWL
Name of Physician(s) Reading Study	(0008,1060)	3	Names of the physician(s) reading the Study. (**) First Name and/or Last Name are copied from user input if entered, otherwise sent as ZERO LENGTH.
Procedure Code Sequence.	(0008,1032)	3	Procedure Code Sequence.(*) Mapped without change from Modality Worklist Requested Procedure Code Sequence (0032,1064). May present in locally protocols appended to study with already scheduled MWL protocol.
> Include 'Code Sequence Macro'			
Referenced Study Sequence	(0008,1110)	3	Referenced Study Sequence.(*) Only a single item is permitted in this sequence, if sent. Copied from 1 <sup>st</sup> valid item of Referenced Study Sequence sent in Worklist. If no valid items exist, not sent. May present in locally protocols appended to study with already scheduled MWL protocol.
>Include 'SOP Instance Reference Macro'			

**Note 1 :** (\*) - Attributes copied from the Worklist if the study source was actually copied from a Worklist query result.

**Note 2 :** (\*\*) - Attributes copied from the user input for Locally scheduled Protocols

**Note 3 :** (\*\*\*) - Cannot be modified by user if received from MWL

3.4.2.2 Patient Study Module

**TABLE 3-5  
PATIENT STUDY MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Admitting Diagnoses Description	(0008,1080)	3	Description of the admitting diagnosis (diagnoses) (**) Sent as ZERO LENGTH value if value is not received from user input for locally scheduled protocols
Patient's Age	(0010,1010)	3	Age of the Patient Calculated from Patient Birth Date if Patient Birth Date is not empty. Not sent is Patient Birth Date is not defined. Cannot be updated if Patient Birth Date is entered from MWL.
Patient's Size	(0010,1020)	3	Length or size of the Patient, in meters. (*) (**) Sent as ZERO LENGTH value if value is not received either from MWL or from user input. User can modify value received from MWL
Patient's Weight	(0010,1030)	3	Weight of the Patient, in kilograms.(*) (**) Sent as ZERO LENGTH value if value is not received either from MWL or from user input. User can modify value received from MWL

**Note 1 :** (\*) - Attributes copied from the Worklist if the study source was actually copied from a Worklist query result.

**Note 2 :** (\*\*) - Attributes copied from the user input for Locally scheduled Protocols

3.4.2.3 Standard Extended Study Module

**TABLE 3-6  
STANDARD EXTENDED STUDY MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Requested Procedure Comments	(0040,1400)	3	User-defined Study notes Sent as ZERO LENGTH value if value is not received from user input
Allergies	(0010,2110)	3	Description of prior reaction to contrast agents, or other patient allergies or adverse reactions (*) (**) Sent as ZERO LENGTH value if value is not received either from MWL or from user input. User can modify value received from MWL

Pregnancy Status	(0010,21C0)	3	<p>Describes pregnancy state of patient. (*)(**)</p> <p>Enumerated Values:</p> <p>0001 = not pregnant</p> <p>0002 = possibly pregnant</p> <p>0003 = definitely pregnant</p> <p>0004 = unknown</p> <p>Sent as “unknown” if value is not received either from MWL or from user input and if tag (0010,1040) - Patient’s Sex – is other than “F” (Female)</p> <p>User can modify value received from MWL if tag (0010,1040) - Patient’s Sex – is “F” (Female) only</p>
------------------	-------------	---	---

**Note 1 :** (\*) – Attributes copied from the Worklist if the study source was actually copied from a Worklist query result.

**Note 2 :** (\*\*) – Attributes copied from the user input for Locally scheduled Protocols

3.4.2.4 Private Study Module

**TABLE 3-7  
PRIVATE STUDY MODULE ATTRIBUTES**

Attribute Name	Tag	Private Creator ID	Attribute Description
Auto-Processing Application	(0009,xx1E)	QUASAR_INTERNAL_USE	Auto-Processing Application Description. Not sent if Auto-processing is not defined for protocol.
Acquisitions Count	(0009,xx1F)	QUASAR_INTERNAL_USE	Number of acquisitions per study. Not sent if information is not provided by protocol
Completed Acquisitions Count	(0009,xx20)	QUASAR_INTERNAL_USE	Number of completed acquisitions per study. Not sent if information is not provided by protocol
Acquisition flag	(0009,xx42)	QUASAR_INTERNAL_USE	Used for indicating if the study is acquired. Always sent. Default Value “acquired”
Patient Unique Key	(0009,xx39)	QUASAR_INTERNAL_USE	Patient unique key. Always sent with non-empty value.

DIR DOC0660253 REV 15

3.4.3 Series Entity Modules

3.4.3.1 General Series Module

TABLE 3-8  
GENERAL SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Modality	(0008,0060)	1	Type of equipment that originally acquired the data used to create the images in this Series.  Defined Terms used for data created on this system: NM = Nuclear Medicine CT = Computed Tomography SR = Structured Report
Series Instance UID	(0020,000E)	1	Internally generated unique identifier of the Series.
Series Number	(0020,0011)	2	A number that identifies this Series.  Internally generated
Laterality	(0020,0060)	2C	Laterality of (paired) body part examined. Copied from the user input.  Not sent if value is not received from user input.  Enumerated Values: R = right L = left
Performing Physicians' Name	(0008,1050)	3	Name of the physician(s) administering this Series. (* (**) Sent as empty string if no user input provided.
Protocol Name	(0018,1030)	3	User-defined description of the conditions under which the Series was performed.  The full path of the performed protocol name. E.g. Factory&Service&Static  Always sent with non-empty value.
Series Description	(0008,103E)	3	Description of the Series  Defined by the acquired Scan name.  Always sent.
Operators' Name	(0008,1070)	3	Name(s) of the operator(s) supporting the Series Copied from user input.  Sent as empty string if no user input provided.
Referenced Performed Procedure Step Sequence	(0008,1111)	3	Uniquely identifies the Performed Procedure Step SOP Instance to which the Series is related. The sequence has exactly 1 item. Sequence is added to all image(s) created by system.  Always sent.
>Referenced SOP Class UID	(0008,1150)	1	Set to "1.2.840.10008.3.1.2.3.3"

DIR DOC0660253 REV 15

>Referenced SOP Instance UID	(0008,1155)	1	Uniquely identifies the referenced SOP Instance. Internally generated.
Body Part Examined	(0018,0015)	3	Text description of the part of the body examined. Sent as empty string if no user input provided. Defined Terms used on this system: ABDOMEN ANKLE ARM BREAST CHEST CLAVICLE COCCYX CSPINE ELBOW EXTREMITY FOOT HAND HEAD HEART HIP JAW KNEE LEG LSPINE NECK PELVIS SHOULDER SKULL SSPINE TSPINE WHOLEBODY NECKCHESTABDPERV
Patient Position	(0018,5100)	2C	Patient position descriptor relative to the Equipment. For NM Images sent in Standard Extended Series Module 3.4.3.2 Sent for CT Images. The Defined Terms used are: HFP = Head First-Prone HFS = Head First-Supine FFP = Feet First-Prone FFS = Feet First-Supine Attribute is copied from the user input.
Request Attributes Sequence	(0040,0275)	3	Sequence that contains attributes from the Imaging Service Request. The sequence has exactly 1 item. Used when Series as created as result of MWL request or protocol, appended to MWL request, not sent otherwise.

DIR DOC0660253 REV 15

>Requested Procedure ID	(0040,1001)	1C	Identifier that identifies the Requested Procedure in the Imaging Service Request.(*)
>Accession Number	(0008,0050)	3	A RIS generated number that identifies the order for the Study(*) May be updated by user manually.
>Study Instance UID	(0020,000D)	3	Unique identifier for the Study. (*) (***)
>Referenced Study Sequence	(0008,1110)	3	Uniquely identifies the Study SOP Instances associated with this SOP Instance. The sequence has exactly 1 item. (*) (***)
>> Include 'SOP Instance Reference Macro'			
>Requested Procedure Description	(0032,1060)	3	Institution-generated administrative description or classification of Requested Procedure. (*) (***)
>Requested Procedure Code Sequence	(0032,1064)	3	Not Used
>Scheduled Procedure Step ID	(0040,0009)	1C	Identifier that identifies the Scheduled Procedure Step. (*) (***)
>Scheduled Procedure Step Description	(0040,0007)	3	Institution-generated description or classification of the Scheduled Procedure Step to be performed. (*) (***)
>Scheduled Protocol Code Sequence	(0040,0008)	3	Sequence describing the Scheduled Protocol following a specific coding scheme. (*) (***) The sequence has exactly 1 item. If MWL request contains more than 1 item, only the first valid item is copied.
>>Include 'Code Sequence Macro'			
Comments on the Performed Procedure Step	(0040,0280)	3	User-defined comments on the Performed Procedure Step (**). Sent as empty string if no user input provided.
Performed Procedure Step ID	(0040,0253)	3	Equipment generated identifier of the protocol carried out within this step. The PPS ID is unique within a study. For MWL scheduled protocols set with "W_" + <SPS ID>. For locally scheduled protocols set with "L_" + numbered id starting from 1 (L_1, L_2). Sent as empty string for QC protocols. Always sent.
Performed Procedure Step Start Date	(0040,0244)	3	The date that the protocol (SPS) acquisition actually started (doesn't matter if the protocol originated from MWL or was locally scheduled). A locally scheduled protocol is an SPS that is created/added in the camera. Always sent.

DIR DOC0660253 REV 15

Performed Procedure Step Start Time	(0040,0245)	3	The time that the protocol (SPS) acquisition actually started (doesn't matter if the protocol originated from MWL or was locally scheduled). A locally scheduled protocol is an SPS that is created/added in the camera. Always sent
Performed Procedure Step Description	(0040,0254)	3	The full path of the performed protocol name. E.g. Factory&Cardiology&One Day Always sent
Performed Protocol Code Sequence	(0040,0260)	3	Not Used Assisted protocol setting is not supported

**Note 1 : (\*)** – Attributes copied from the Worklist if the study source was actually copied from a Worklist query result (if available).

**Note 2 : (\*\*)** – Attributes copied from the user input for Locally scheduled Protocols

**Note 3 : (\*\*\*)** – Cannot be modified by user if received from MWL

3.4.3.2 Standard Extended Series Module

TABLE 3-9  
STANDARD EXTENDED SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Patient Position	(0018,5100)	3	Patient position descriptor relative to the Equipment: The Defined Terms are: HFP = Head First-Prone HFS = Head First-Supine HFDR = Head First-Decubitus Right HFDL = Head First-Decubitus Left FFDR = Feet First-Decubitus Right FFDL = Feet First-Decubitus Left FFP = Feet First-Prone FFS = Feet First-Supine Attribute is copied from the user input. Not sent, if another Patient Position is used for acquisition.

3.4.3.3 NM/PET Patient Orientation Module

TABLE 3-10  
NM/PET PATIENT ORIENTATION MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
----------------	-----	------	-----------------------



DIR DOC0660253 REV 15

Patient Orientation Code Sequence	(0054,0410)	2	<p>Describes the orientation of the patient with respect to gravity.</p> <p>Contains 1 items if Patient Position is one of the following:</p> <p>HFP = Head First-Prone                      HFS = Head First-Supine                      HFDR = Head First-Decubitus Right                      HFDL = Head First-Decubitus Left                      FFDR = Feet First-Decubitus Right                      FFDL = Feet First-Decubitus Left                      FFP = Feet First-Prone                      FFS = Feet First-Supine</p> <p>Otherwise contain 0 items.</p>
>Include Code Sequence Macro			Always (F-10450,99SDM, "recumbent") is sent
> Patient Orientation Modifier Code Sequence	(0054,0412)	2C	Patient Orientation Modifier. Required if needed to fully specify the orientation of the patient with respect to gravity. Contains exactly 1 item, if sent
>>Include 'Code Sequence Macro'			<p>The following codes are supported:</p> <p>(F-10310, 99SDM, "prone")                      (F-10340, 99SDM, "supine")                      (F-10317, 99SDM, "right lateral decubitus")                      (F-10319, 99SDM, "left lateral decubitus")</p>
Patient Gantry Relationship Code Sequence	(0054,0414)	2	<p>Describes the orientation of the patient with respect to the gantry. Contains 1 item if Patient Position is one of the following:</p> <p>HFP = Head First-Prone                      HFS = Head First-Supine                      HFDR = Head First-Decubitus Right                      HFDL = Head First-Decubitus Left                      FFDR = Feet First-Decubitus Right                      FFDL = Feet First-Decubitus Left                      FFP = Feet First-Prone                      FFS = Feet First-Supine</p> <p>Otherwise contain 0 items.</p>
>Include Code Sequence Macro			<p>The following codes are supported:</p> <p>(F-10470, 99SDM, "headfirst")                      (F-10480, 99SDM, "feet-first")</p>

DIR DOC0660253 REV 15

3.4.3.4 Private Series Module

TABLE 3-11  
PRIVATE SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Sequence Type	(0009,xx13)	QUASAR_INTERNAL_USE	Acquired Sequence Type. Always sent as non-empty value
Sequence Name	(0009,xx14)	QUASAR_INTERNAL_USE	Acquired Sequence Name. Always sent as non-empty value
Administrations Plan	(0009, xx25)	QUASAR_INTERNAL_USE	List of Scan Planned Administrations
Protocol Scheduled Date	(0009,xx40)	QUASAR_INTERNAL_USE	Protocol Scheduled Date (*). Taken from the SPS Start date of the first SPS in the study – Tag (0040,0002).  For QC scans – sent as ZERO Length Always sent
Protocol Scheduled Time	(0009,xx41)	QUASAR_INTERNAL_USE	Protocol Scheduled Time (*) Taken from the SPS Start time of the first SPS in the study – Tag (0040,0003).  For QC scans – sent as ZERO Length. Always sent
Matched protocol	(0009,xx43)	QUASAR_INTERNAL_USE	For Worklist items. The originally matched protocol vs. protocol name which is the protocol actually acquired.  Sent if originally matched protocol is modified by user, otherwise not sent  Not sent for locally scheduled protocols.
Private SPS ID	(0009,xx44)	QUASAR_INTERNAL_USE	Defines protocols that were appended to the original MWL protocol.  Contains SPS ID of original MWL Protocol.  Sent empty for images of locally scheduled protocols and appended to originally locally scheduled protocols.
Pre-Medication	(0009,xx45)	QUASAR_INTERNAL_USE	Keeps the Pre-Medication as appears in the “To-Do” list (*) (**)  Sent as empty string if not received from user input.
Anatomic Reference	(0009,xx48)	QUASAR_INTERNAL_USE	Keeps the anatomic reference for the specific scan (*). Not sent if not received from user input.

Acquisition Mode	(0009,xx59)	QUASAR_INTERNAL_USE	Acquisition mode. Possible Values: NORMAL_ACQ = 0 GATED_ACQ = 1 TOMO_ACQ = 2 ENERGY_ACQ = 3 WB_CONT_ACQ = 4 DYNAMIC_TOMO_ACQ = 5 TOMO_ACQ_DURING_STEP_ACQ = 6 GATED_TOMO_ACQ_DURING_STEP_ACQ = 7
Series Data Sequence	(0033,xx70)	GEMS_XELPRV_01	Sequence of item contains information about acquisition parameters. May contain from 1 or more items. Each item describes specific parameters set. Sent for NM data only
>Object Type	(0033,xx08)	GEMS_XELPRV_01	Object Type. Contains string "SERIES DATA "
>Modified Flag	(0033,xx10)	GEMS_XELPRV_01	Default value = 0 (Not Modified)
>Name	(0033,xx11)	GEMS_XELPRV_01	SDO Name
>Database Object Unique ID	(0033,xx16)	GEMS_XELPRV_01	Database UID of SDO; contains value of SDO UID tag (0033,xx72) generated at time of object creation.
>Date	(0033,xx17)	GEMS_XELPRV_01	SDO Creation date
>Time	(0033,xx18)	GEMS_XELPRV_01	SDO Creation time
>Series Data Flags	(0033,xx19)	GEMS_XELPRV_01	SDO Flags. Default value = 0
>Protocol Name	(0033,xx1A)	GEMS_XELPRV_01	Name of Protocol created SDO
>Relevant Data UID	(0033,xx1B)	GEMS_XELPRV_01	UID(s) of SOP Instance(s) relative to SDO
>Bulk Data	(0033,xx1C)	GEMS_XELPRV_01	SDO parameter(s) stored as binary buffer(s)
>Int Data	(0033,xx1D)	GEMS_XELPRV_01	List of SDO parameters stored as integers
>Double Data	(0033,xx1E)	GEMS_XELPRV_01	List of SDO parameters stored as doubles
>String Data	(0033,xx1F)	GEMS_XELPRV_01	List of SDO parameters stored as list of strings
>Bulk Data Format	(0033,xx20)	GEMS_XELPRV_01	Format of bulk parameters; contains information about name and size of bulk buffers
>Int Data Format	(0033,xx21)	GEMS_XELPRV_01	Format of integer parameters; contains information about name and number of integers in list
>Double Data Format	(0033,xx22)	GEMS_XELPRV_01	Format of double parameters; contains information about name and number of doubles in list
>String Data Format	(0033,xx23)	GEMS_XELPRV_01	Format of string parameters; contains information about name and number of strings in list
>Description	(0033,xx24)	GEMS_XELPRV_01	User or equipment generated SDO description
>SDO Private SOP Class UID	(0033,xx71)	GEMS_XELPRV_01	SDO Private SOP Class UID- "1.2.840.113619.4.17"
>SDO Instance UID	(0033,xx72)	GEMS_XELPRV_01	SDO Instance UID; Internally generated

DIR DOC0660253 REV 15

**Note 1 : (\*)** – *Attributes copied from the Worklist if the study source was actually copied from a Worklist query result (if available).*

**Note 2 : (\*\*)** – *Attributes copied from the user input for Locally scheduled Protocols*

### 3.4.4 Frame Of Reference Entity Modules

#### 3.4.4.1 Frame Of Reference Module

This section specifies the Attributes necessary to uniquely identify a Frame Of Reference which insures the spatial relationship of Images within a Series. It also allows Images across multiple Series to share the same Frame Of Reference. This Frame Of Reference (or coordinate system) shall be constant for all Images related to a specific Frame Of Reference.

A hybrid CT/NM scan is composed of a single NM scan partnered with one or more CT scans. The two modalities share the same imaging space and the body imaged by the two modalities is represented, in most of the cases, by spatially aligned images. There are situations for which optimal NM imaging and optimal CT imaging impose changing the table height during the hybrid scan. In this case, the imaging space of both modalities remains the same, but the NM and CT images of the body are no longer spatially aligned. In order to prevent accidental fusion of such images, the same Frame Of Reference UID value shared by two series of different modalities will show that the images are spatially related and that the imaged body was scanned spatially aligned between the two images.

With this approach the handling of Frame Of Reference UID is done as follows:

1. In all combinations, images acquired with the same landmark and same table height will have the same Frame of Reference UID value. The generated images are “inherently” aligned.
2. Discovery NM/CT 670\* and Optima NM/CT 640 uses a private tag named “Acquisition Parent UID”:
  - a. In CT images (Discovery NM/CT 670\* and Optima NM/CT 640 implementation), it is (0031, xx02) Private Creator “GEHC\_HYBRID\_01”.
  - b. In NM images (Discovery NM/CT 670\* and Optima NM/CT 640 implementations) and CT Images (Optima NM/CT 640 implementation), it is (0011,xx31) Private Creator “GEMS\_GENIE\_1”.

All the images (NM projections and CT slices obtained by on-the-fly reconstruction or retrospective reconstruction) resulting from scans defined as part of a single hybrid procedure and using the same landmark shall have the same Acquisition Parent UID value. Hence, when the operator sets a new CT landmark, overriding the landmark set on NM persistence or during a previous CT scan belonging to the same hybrid procedure, a new value is assigned to the Acquisition Parent UID attribute. The Acquisition Parent UID attribute allows Xeleris workstation, or other Image Displayer implementations, to check that CT and NM series were acquired as part of the same hybrid procedure and are using the same landmark.

3. When NM is acquired first, and the CT series is acquired at a different table height than used by NM, but with same landmark, then the CT series will get a different Frame Of Reference UID value than that used by NM, but the same Acquisition Parent UID value as NM. This shows that the imaged body is no longer aligned between the two modalities, although they have been acquired by a single hybrid scan (using the same imaging space). In this case, image registration must be applied in order to register the two modalities. The information stored by each modality in Table Height (0018,1130) can be used in this registration process.
4. When CT is acquired first, and the NM series is acquired at a different table height than used by CT, but with same landmark, the NM series will get the same Frame Of Reference UID value than that used by CT, the same Acquisition Parent UID value as CT and the image position of the NM will be modified to compensate for the table height differences existing between the two modalities. Therefore, in this case, the imaged body in

DIR DOC0660253 REV 15

aligned between the two modalities and the same Frame Of Reference UID value shows it.

5. CT images acquired with a different landmark than the NM scan or than previous CT series of the same hybrid scan will get a Frame of Reference UID value different of the NM scan or of the previous CT series, showing that the images are not spatially related. The Acquisition Parent UID value of these images will be also changed, as indicated above at 2. In this way, Xeleris (and potentially other image displayers) will be informed that these images are not spatially aligned and image registration must be applied in order to achieve it.

Note that for CT standalone exams (Discovery NM/CT 670\* implementations), the same Frame of Reference UID (0020,0052) value does not imply imaged body alignment, but just image space alignment. This attribute must be used in conjunction with the Table Height (0018,1130) to determine if two imaged bodies are spatially aligned. See CT Scanner DICOM Conformance Statement. ([Reference A](#) in Section 1.6).

The Frame of Reference Module Attributes appear for TOMO and GATED TOMO scan types. They are no available for STATIC, WHOLE BODY, GATED and DYNAMIC scans.

**TABLE 3-12**  
**FRAME OF REFERENCE MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Frame of Reference UID	(0020,0052)	1	Uniquely identifies the frame of reference for a Series. See explanation above
Position Reference Indicator	(0020,1040)	2	Sent as ZERO LENGTH value if it is not received from user input.

### 3.4.5 Equipment Entity Modules

#### 3.4.5.1 General Equipment Module

**TABLE 3-13**  
**GENERAL EQUIPMENT MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Manufacturer	(0008,0070)	2	Manufacturer of the equipment that produced the composite instances. Default Value "GE MEDICAL SYSTEMS" Always sent
Institution Name	(0008,0080)	3	Institution where the equipment that produced the composite instances is located. Always taken from system configuration. Sent as empty string if it is not received from user input. Always sent.
Institution Address	(0008,0081)	3	Mailing address of the institution where the equipment that produced the composite instances is located. Always taken from system configuration Sent as empty string if it is not received from user input. Always sent.

DIR DOC0660253 REV 15

Institutional Department Name	(0008,1040)	3	<p>Department in the institution where the equipment that produced the composite instances is located.</p> <p>Always taken from system configuration</p> <p>Sent as empty string if it is not received from user input.</p> <p>Always sent.</p>
Manufacturer's Model Name	(0008,1090)	3	<p>Manufacturer's model name of the equipment that produced the composite instances.</p> <p>Set to "Tandem_Discovery_670" for Discovery NM/CT 670 implementation</p> <p>Set to "Tandem_Discovery_670_ES" for Discovery NM/CT 670 ES implementation</p> <p>Set to "Tandem_Discovery_670_Pro" for Discovery NM/CT 670 Pro and Discovery 670 DR implementations.</p> <p>Set to "Tandem_Apollo" for Discovery NM/CT 670 CZT implementation</p> <p>Set to "Tandem_Discovery_630" for Discovery NM 630 implementation</p> <p>Set to "Tandem_Optima_640" for Optima NM/CT 640 implementation</p> <p>Set to "Tandem_Brivo_615" for Brivo NM 615 implementation</p> <p>Set to "Tandem_830" for NM 830 implementation</p> <p>Set to "Tandem_850" for NM/CT 850 implementation</p> <p>Set to "Tandem_860" for NM/CT 860 implementation</p> <p>Set to "Tandem_870_DR" for NM/CT 870DR implementation</p> <p>Set to "Tandem_870_CZT" for NM/CT 870CZT implementation</p> <p>Always sent</p>
Device Serial Number	(0018,1000)	3	<p>Manufacturer's serial number of the equipment that produced the composite instances.</p> <p>Always taken from system configuration</p> <p>Sent as empty string if it is not received from user input.</p> <p>Always sent</p>
Station Name	(0008,1010)	3	<p>User defined name identifying the machine that produced the composite instances.</p> <p>Camera name taken from configuration is used</p> <p>Sent as empty string if it is not received from user input.</p> <p>Always sent</p>
Software Versions	(0018,1020)	3	<p>Manufacturer's designation of software version of the equipment that produced the composite instances</p> <p>Software/Hardware versions of current release</p> <p>e.g. "1.003.317.1\HARDWARE_VERSION_1"</p> <p>Always sent as non-empty value.</p>

DIR DOC0660253 REV 15

Spatial Resolution	(0018,1050)	3	The inherent limiting resolution in mm of the acquisition equipment for high contrast objects for the data gathering and reconstruction technique chosen. Always taken from system configuration. Default value is 5. Always sent
Date of Last Calibration	(0018,1200)	3	Date when the image acquisition device calibration was last changed in any way. Always taken from system configuration Always Sent
Time of Last Calibration	(0018,1201)	3	Time when the image acquisition device calibration was last changed in any way. Always taken from system configuration Always Sent

### 3.4.6 Image Entity Modules

#### 3.4.6.1 General Image Module

TABLE 3-14  
GENERAL IMAGE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Instance Number	(0020,0013)	2	A number that identifies this image. Sent as ZERO Length value.
Patient Orientation	(0020,0020)	2C	Not sent for NM (not required )
Content Date	(0008,0023)	2C	The date the image pixel data creation started. Always sent if image is part of a series in which the images are temporally related. May be sent otherwise.
Content Time	(0008,0033)	2C	The time the image pixel data creation started. Always sent if image is part of a series in which the images are temporally related. May be sent otherwise.
Image Type	(0008,0008)	3	See 3.4.6.7.1
Acquisition Date	(0008,0022)	3	The date the acquisition of data that resulted in this image started Always sent with non-empty value
Acquisition Time	(0008,0032)	3	The time the acquisition of data that resulted in this image started Always sent with non-empty value
Image Comments	(0020,4000)	3	Contains additional information about image. Always sent with non-empty value

Quality Control Image	(0028,0300)	3	Indicates whether or not this image is a quality control or phantom image. Enumerated Values: YES NO Always sent
-----------------------	-------------	---	---

3.4.6.2 Image Pixel Module

TABLE 3-15  
IMAGE PIXEL MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Samples per Pixel	(0028,0002)	1	See 3.4.6.4 for NM Images
Photometric Interpretation	(0028,0004)	1	See 3.4.6.4 for NM Images
Rows	(0028,0010)	1	Number of rows in the image.
Columns	(0028,0011)	1	Number of columns in the image
Bits Allocated	(0028,0100)	1	See 3.4.6.4 for NM Images
Bits Stored	(0028,0101)	1	See 3.4.6.4 for NM Images
High Bit	(0028,0102)	1	See 3.4.6.4 for NM Images
Pixel Representation	(0028,0103)	1	Data representation of the pixel samples. Each sample shall have the same pixel representation. Enumerated Values used: 0000H = unsigned integer.
Pixel Data	(7FE0,0010)	1	A data stream of the pixel samples that comprise the Image.
Planar Configuration	(0028,0006)	1C	Not Used (number of Samples per Pixel is always 1)
Pixel Aspect Ratio	(0028,0034)	1C	Not Used
Smallest Image Pixel Value	(0028,0106)	3	The minimum actual pixel value encountered in this image.
Largest Image Pixel Value	(0028,0107)	3	The maximum actual pixel value encountered in this image.

3.4.6.3 Acquisition Context Module

This section specifies Attributes for the description of the conditions present during data acquisition.

TABLE 3-16  
ACQUISITION CONTEXT MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Acquisition Context Sequence	(0040,0555)	2	A sequence of Items that describes the conditions present during the acquisition of the data of the SOP Instance.  The Acquisition context sequence contains 0 items when acquisition context in scan is left "UNKNOWN", otherwise contains 1 item.



DIR DOC0660253 REV 15

>Concept Name Code Sequence	(0040,A043)	1	A concept that constrains the meaning of (i.e. defines the role of) the Observation Value. This sequence contains 1 item
>>Include 'Code Sequence Macro'			(109054, DCM, "Patient State") is supported as defined in TID 3470
>Concept Code Sequence	(0040,A168)	1C	This is the Value component of a Name/Value pair when the Concept implied by Concept Name Code Sequence (0040,A043) is a Coded Value. This sequence contains 1 item
>>Include 'Code Sequence Macro'			DCID (3101) NM Procedural State Values is supported as defined in TID 3470: The following values are used: <ul style="list-style-type: none"> <li>• (F-01604 ,SRT , "Resting State")</li> <li>• (F-05019 ,SRT , "Cardiac Stress State")</li> <li>• (109092 ,DCM , "Reinjection State")</li> <li>• (109093 ,DCM , "Redistribution State")</li> <li>• (109094 ,DCM , "Delayed Redistribution State")</li> </ul>

3.4.6.4 NM Image Pixel Module

This section specifies the Attributes that describe the pixel data of a NM image.

TABLE 3-17  
NM IMAGE PIXEL MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Samples per Pixel	(0028,0002)	1	Number of samples (planes) in this image. The value always set to 1.
Photometric Interpretation	(0028,0004)	1	Specifies the intended interpretation of the pixel data Enumerated Values supported : MONOCHROME2
Bits Allocated	(0028,0100)	1	Number of bits allocated for each pixel sample. Each sample shall have the same number of bits allocated. Enumerated Values supported : 16.
Bits Stored	(0028,0101)	1	Number of bits stored for each pixel sample. Value equal to Bit Allocated (0028,0100)
High Bit	(0028,0102)	1	Most significant bit for pixel sample data. Value equal to Bit Stored (0028,0101) – 1

Pixel Spacing	(0028,0030)	2	Physical distance in the patient between the center of each pixel, specified by a numeric pair – adjacent row spacing (delimiter) adjacent column spacing, in mm.
---------------	-------------	---	---

3.4.6.5 Multi-Frame Module

TABLE 3-18  
MULTI-FRAME MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Number of Frames	(0028,0008)	1	Number of frames in a Multi-frame Image.
Frame Increment Pointer	(0028,0009)	1	See 3.4.6.6.1 for further specialization.

3.4.6.6 NM Multi-frame Module

TABLE 3-19  
NM MULTI-FRAME MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Frame Increment Pointer	(0028,0009)	1	See 3.4.6.6.1 for further specialization.
Energy Window Vector	(0054,0010)	1C	Defines energy set window to which each frame belongs. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Energy Window Vector (0054,0010).
Number of Energy Windows	(0054,0011)	1	Number of energy set windows in SOP Instance. Possible values: 1 to 16
Detector Vector	(0054,0020)	1C	Defines detector to which each frame belongs. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Detector Vector (0054,0020).
Number of Detectors	(0054,0021)	1	Number of detectors in SOP Instance. Possible values: 1 or 2. For Brivo NM 615 DICOM implementation the possible value is 1.
Phase Vector	(0054,0030)	1C	Defines phase to which each frame belongs. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Phase Vector (0054,0030).
Number of Phases	(0054,0031)	1C	Number of phases in SOP Instance. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Phase Vector (0054,0030). Supported values: 1, 2, 3, 4, 5
Rotation Vector	(0054,0050)	1C	Defines rotation to which each frame belongs. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Rotation Vector (0054,0050).

DIR DOC0660253 REV 15

Number of Rotations	(0054,0051)	1C	Number of Rotations in SOP Instance. Always set to 1. Sent if Image Type (0008,0008), Value 3 is TOMO and GATED TOMO.
R-R Interval Vector	(0054,0060)	1C	Defines R-R Interval to which each frame belongs. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for R-R Interval Vector (0054,0060).
Number of R-R Intervals	(0054,0061)	1C	Number of R-R Intervals in SOP Instance. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for R-R Interval Vector (0054,0060).
Time Slot Vector	(0054,0070)	1C	Defines time slot, within cardiac cycle, to which each frame belongs. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Time Slot Vector (0054,0070).
Number of Time Slots	(0054,0071)	1C	Number of time slots in SOP Instance. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Time Slot Vector (0054,0070).
Angular View Vector	(0054,0090)	1C	Defines angular view number to which each frame belongs. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Angular View Vector (0054,0090).
Time Slice Vector	(0054,0100)	1C	Defines frame numbers within each phase. Sent if the value of the Frame Increment Pointer (0028,0009) includes the Tag for Time Slice Vector (0054,0100).

**3.4.6.6.1 Frame Increment Pointer**

The Frame Increment Pointer (0028,0009) defines which frame index vectors are present in the NM Image instance. The Frame Increment Pointer is supported per the DICOM specification for all image types defined in Table 3-20.

**TABLE 3-20  
ENUMERATED VALUES FOR FRAME INCREMENT POINTER**

<b>Image Type (0008,0008), Value 3</b>	<b>Frame Increment Pointer (0028,0009)</b>
STATIC or WHOLE BODY	0054H 0010H \ 0054H 0020H Sequencing is by Energy Window Vector (0054,0010), Detector Vector (0054,0020).
DYNAMIC	0054H 0010H \ 0054H 0020H \ 0054H 0030H \ 0054H 0100H Sequencing is by Energy Window Vector (0054,0010), Detector Vector (0054,0020), Phase Vector (0054,0030), Time Slice Vector (0054,0100)
GATED	0054H 0010H \ 0054H 0020H \ 0054H 0060H \ 0054H 0070H Sequencing is by Energy Window Vector (0054,0010), Detector Vector (0054,0020), R-R Interval Vector (0054,0060), Time Slot Vector (0054,0070)

DIR DOC0660253 REV 15

TOMO	0054H 0010H \ 0054H 0020H \ 0054H 0050H \ 0054H 0090H Sequencing is by Energy Window Vector (0054,0010), Detector Vector (0054,0020), Rotation Vector (0054,0050), Angular View Vector (0054,0090)
GATED TOMO	0054H 0010H \ 0054H 0020H \ 0054H 0050H \ 0054H 0060H \ 0054H 0070H \ 0054H 0090H Sequencing is by Energy Window Vector (0054,0010), Detector Vector (0054,0020), Rotation Vector (0054,0050), R-R Interval Vector (0054,0060), Time Slot Vector (0054,0070), Angular View Vector (0054,0090).

3.4.6.7 NM Image Module

TABLE 3-21  
NM IMAGE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Image Type	(0008,0008)	1	See 3.4.6.7.1 for specialization.
Image ID	(0054,0400)	3	User or equipment generated Image identifier. Taken from Scan name
Counts Accumulated	(0018,0070)	2	Sum of all gamma events for all frames in the image.
Acquisition Termination Condition	(0018,0071)	3	Description of how the data collection was stopped. (*) Defined Terms are used: CNTS = counts DENS = density, count limit reached within ROI MANU = manual TIME = time TRIG = physiological trigger Always sent as non-empty value
Actual Frame Duration	(0018,1242)	1C	Elapsed time for one frame acquisition in msec.. Sent when the Image Type (0008,0008), Value 3, is equal to STATIC or WHOLE BODY.
Count Rate	(0018,1243)	3	Maximum count rate achieved during the acquisition in counts/sec Always sent as non-empty value

DIR DOC0660253 REV 15

Corrected Image	(0028, 0051)	3	<p>Corrections have been applied to the image.</p> <p>Defined Terms are used:</p> <ul style="list-style-type: none"> <li>UNIF = flood corrected</li> <li>COR = center of rotation corrected</li> <li>ATTN = attenuation corrected</li> <li>SCAT = scatter corrected</li> <li>NRGY = energy corrected</li> <li>LIN = linearity corrected</li> <li>CLN = count loss normalization</li> <li>BADP = bad Pixels correction (for Discovery NM/CT 670 CZT only)</li> </ul> <p>Always sent</p>
Table Height	(0018,1130)	3	<p>The height of the patient table in mm.</p> <p>Not sent when the Image Type (0008,0008), Value 3, is equal to TOMO or GATED TOMO.</p> <p>Otherwise sent as non-empty value</p>
Table Traverse	(0018,1131)	3	<p>Location of the patient table in mm. Not sent when the Image Type (0008,0008), Value 3, is equal to TOMO or GATED TOMO.</p> <p>Otherwise sent as non-empty value</p>
Whole Body Technique	(0018,1301)	3	<p>The type of scan performed. Sent if Image Type (0008,0008), Value 3, contains the value WHOLE BODY.</p> <p>Enumerated Values used:</p> <ul style="list-style-type: none"> <li>1PS = one pass - for Continues Acquisition.</li> <li>MSP= multiple static frames collected into a whole body frame – for Step And Shoot acquisition.</li> </ul>
Scan Velocity	(0018,1300)	2C	<p>The speed of the camera motion over the body in mm/sec. Sent as ZERO Length value if Image Type (0008,0008) Value 3 contains the value WHOLE BODY and the WHOLEBODY is acquired in the Step-and-Shoot technique, where the table does not move during acquisition.</p>
Scan Length	(0018,1302)	2C	<p>Size of the imaged area in the direction of scanning motion, in mm.</p> <p>Sent as ZERO Length value if Image Type (0008,0008) Value 3 contains the value WHOLE BODY and the WHOLEBODY is acquired in the Step-and-Shoot technique, where the table does not move during acquisition.</p>

**Note1:** (\*) – Attribute value is taken from user input

DIR DOC0660253 REV 15

**3.4.6.7.1 Image Type**

The following values of Image Type (0008,0008) are be sent :

Value 1 shall have the following Enumerated Values:

- ORIGINAL identifies an Original Image
- DERIVED identifies a Derived Image

Value 2 shall have the following Enumerated Value:

- PRIMARY identifies a Primary Image

The following Enumerated Values of Value 3 are created:

- STATIC - Identifies a Static Image
- DYNAMIC - Identifies a Dynamic Image
- GATED - Identifies a Multi-Gated Image
- WHOLE BODY - Identifies a Whole Body Image
- TOMO - Identifies a Tomographic (SPECT) Image
- GATED TOMO - Identifies a Multi-gated Tomographic Image

The following Enumerated Values of Value 4 are created:

- EMISSION - Transmission source is NOT active during image acquisition

3.4.6.8 NM Isotope Module

This section contains Attributes that describe the isotope administered for the acquisition.

TABLE 3-22  
NM ISOTOPE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Energy Window Information Sequence	(0054,0012)	2	Sequence of Items that describe the energy window groups used. Contain from 1 to 4 Items. The number of items shall be equal to Number of Energy Windows (0054,0011)
> Energy Window Name	(0054,0018)	3	A user defined name which describes this Energy Window.
>Energy Window Range Sequence	(0054,0013)	3	Sequence describing window energy limits. Contains from 1 to 16 items.
>> Energy Window Lower Limit	(0054,0014)	3	The lower limit of the energy window in KeV.
>> Energy Window Upper Limit	(0054,0015)	3	The upper limit of the energy window in KeV.
Radiopharmaceutical Information Sequence	(0054,0016)	2	Information on radiopharmaceutical(s) used. May contain from 0 to 3 items
> Radionuclide Code Sequence	(0054,0300)	2	Sequence that identifies the radionuclide. Always contains 0 items for Brivo NM 615 implementation. Contains 1 item for all other DICOM implementations.
>> <i>Include 'Code Sequence Macro'</i>			Not Used for Brivo NM 615 implementation. For all other implementations: Baseline CID is 18"Isotopes in Radiopharmaceuticals". Extended values are used (see 2.8.1) The Coding Scheme Designator (0008,0102) has value "99SDM"
> Radiopharmaceutical Route	(0018,1070)	3	Route of injection. (*)
> Administration Route Code Sequence	(0054,0302)	3	Sequence that identifies the administration route for the radiopharmaceutical. Not Used for Brivo NM 615 implementation. Contains 1 item for all other DICOM implementations, if Radiopharmaceutical Route(0018,1070) is sent with non-empty value. Code Meaning (0008,0104) is sent.
>> <i>Include 'Code Sequence Macro'</i>			Not Used for Brivo NM 615 implementation. For all other implementations: Baseline CID 11 "Route of Administration"
> Radiopharmaceutical Volume	(0018,1071)	3	Volume of injection in cubic cm. (*)

DIR DOC0660253 REV 15

> Radiopharmaceutical Start Time	(0018,1072)	3	Time of start of injection. (*)
> Radiopharmaceutical Stop Time	(0018,1073)	3	Time of end of injection. (*) Used in Brivo NM 615 implementations only.
> Radionuclide Total Dose	(0018,1074)	3	Total amount of radionuclide injected in MBq (*)
> Radiopharmaceutical Administration Event UID	(0008,3012)	3	Unique identification of the administration of the radiopharmaceutical to the patient.
> Radiopharmaceutical	(0018,0031)	3	Name of the radiopharmaceutical. (*)
> Radiopharmaceutical Code Sequence	(0054,0304)	3	Sequence that identifies the radiopharmaceutical. Not Used for Brivo NM 615 implementation. Contains 1 item for all other DICOM implementations, if Radiopharmaceutical (0018,0031) is sent with non-empty value.
>> <i>Include 'Code Sequence Macro'</i>			Not Used for Brivo NM 615 implementation. For all other implementations: Baseline CID 25 "Radiopharmaceuticals". Extended values are used (see 2.8.1) Code Meaning (0008,0104) is sent.
Intervention Drug Information Sequence	(0018,0026)	3	Sequence of Items that describes the intervention drugs used. May contain from 1 to 3 Items. Not sent if no user input exists.
>Intervention Drug Name	(0018,0034)	3	Name of intervention drug. (*)
>Intervention Drug Code Sequence	(0018,0029)	3	Not Used
>Administration Route Code Sequence	(0054,0302)	3	Not Used
>Intervention Drug Start Time	(0018,0035)	3	Time of administration of the intervention drug, using the same time base as for the Acquisition Start Time (0008,0032). (*)
>Intervention Drug Stop Time	(0018,0027)	3	Time of completion of administration of the intervention drug, using the same time base as for the Acquisition Start Time (0008,0032). (*)
>Intervention Drug Dose	(0018,0028)	3	Intervention drug dose, in mg. (*)

**Note1:** (\*) – Attribute value is taken from user input if it's not empty, otherwise is not sent

3.4.6.9 NM Detector Module

TABLE 3-23  
NM DETECTOR MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Detector Information Sequence	(0054,0022)	2	Sequence of Items that describe the detectors used. Contains 1 or 2 Items.



DIR DOC0660253 REV 15

> Collimator/Grid Name	(0018,1180)	3	Label describing the collimator used, e.g. LEGP LEHR, etc. Always sent as non-empty value.
> Collimator Type	(0018,1181)	2	Collimator type. Defined Terms: PARA = Parallel (default) NONE = No collimator UNKN = Unknown Always sent as non-empty value
> Focal Distance	(0018,1182)	2	Focal distance, in mm. Default value is 0.
> Zoom Center	(0028,0032)	3	The amount of offset from (0, 0) applied to each pixel in the image before application of the zoom factor, specified by a numeric pair (in mm).(*) Always sent as non-empty value
> Zoom Factor	(0028,0031)	3	The amount of magnification applied to each pixel in the image.(*) Typical Range 1.0 to 4.0 Always sent as non-empty value
> Distance Source to Detector	(0018,1110)	2C	Not sent
> Start Angle	(0054,0200)	3	Position of the detector about the patient for the start of the acquisition, in degrees. (*) Sent if Image Type (0008,0008), Value 3, is other than TOMO, GATED TOMO. If Image Type (0008,0008), Value 3, is WHOLE BODY, sent when Whole Body Technique (0018,1301) is "MSP".
> Radial Position	(0018,1142)	3	Not Sent
> Image Orientation (Patient)	(0020,0037)	2	The direction cosines of the first row and the first column with respect to the patient. Set for first frame in dataset
> Image Position (Patient)	(0020,0032)	2	The x, y, and z coordinates of the upper left hand corner (center of the first voxel transmitted) of the image, in mm. Set for first frame in dataset.
> View Code Sequence	(0054,0220)	3	Sequence that describes the projection of the anatomic region of interest on the image receptor. May contain 0 or 1 item. Sent when Image Type (0008,0008) value3 contains WHOLE BODY, otherwise is not sent. Used in all configurations except Brivo NM 615.
>>Include 'Code Sequence Macro'			CID 26 Nuclear Medicine Projections is supported The following values are used: (SRT, G-5215, Anterior projection) (SRT, G-5216, Posterior projection)

**Note1:** (\*) – Attribute value is taken from user input

**3.4.6.10 NM Tomo Acquisition Module**

This module is present when the Image Type (0008,0008) Value 3, is equal to TOMO or GATED TOMO.

**TABLE 3-24  
NM TOMO ACQUISITION MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Rotation Information Sequence	(0054,0052)	2	Sequence of Items that describe TOMO rotational groups. Contain only 1 item.
> Start Angle	(0054,0200)	1	Position of the detector about the patient for the start of the acquisition, in degrees. (*)
> Angular Step	(0018,1144)	1	The angular scan arc step between views of the TOMO acquisition, in degrees (*)
> Rotation Direction	(0018,1140)	1	Direction of rotation of the detector about the patient. (*) Enumerated Values: CW = clockwise (decreasing angle) CC = counter-clockwise (increasing angle).
> Scan Arc	(0018,1143)	1	The effective angular range of the scan data in degrees. (*) The value is always positive.
> Actual Frame Duration	(0018,1242)	1	Nominal acquisition time per angular position, in msec.
> Radial Position	(0018,1142)	3	Radial distance of the detector from the center of rotation, in mm. Sent as list – one value per angular view.
> Distance Source to Detector	(0018,1110)	2C	Not sent
> Number of Frames in Rotation	(0054,0053)	1	Number of angular views in this rotation.(*)
> Table Traverse	(0018,1131)	3	Table longitudinal position at acquisition start .in mm. Always sent.
> Table Height	(0018,1130)	3	The distance in mm of the top of the patient table to the center of rotation.  Sent as non-empty value in Discovery NM/CT 670* and Optima NM/CT 640 implementations.  In a hybrid scan it is expressed relative to CT iso center, and in non-hybrid scan it is expressed as in normal NM only – measured from floor and up.

DIR DOC0660253 REV 15

Type of Detector Motion	(0054,0202)	3	Describes the detector motion during acquisition.(*). Enumerated Values: STEP AND SHOOT = Interrupted motion, acquire only while stationary. CONTINUOUS = Gantry motion and acquisition are simultaneous and continuous. ACQ DURING STEP= Interrupted motion, acquire while both stationary and motion Always sent as non-empty value
-------------------------	-------------	---	---

**Note1:** (\*) – Attribute value is taken from user input.

**3.4.6.11 NM Multi-gated Acquisition Module**

Describe the conditions under which this module is present in this implementation Module applies to a GATED Multi-frame Image. This module is present when the Image Type (0008,0008) Value 3, is equal to GATED or GATED TOMO

**TABLE 3-25  
NM MULTI-GATED ACQUISITION MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Beat Rejection Flag	(0018,1080)	3	Heart beat duration sorting has been applied. Enumerated Values: Y = yes N = no Always sent.
PVC Rejection	(0018,1085)	3	Description of type of arrhythmic beat rejection criteria used. Always sent as “Reject beats out of pvc window”
Skip Beats	(0018,1086)	3	Number of beats skipped after a detected arrhythmia Always sent.
Heart Rate	(0018,1088)	3	Average number of heart beats per minute for the collection period for these frames Always sent.
Gated Information Sequence	(0054,0062)	2C	Sequence of Items that describe R-R intervals. Sent if the Frame Increment Pointer (0028,0009) contains the Tag for R-R Interval Vector (0054,0060) Contains only 1 item if presents
> Cardiac Framing Type	(0018,1064)	3	Description of type of framing performed. Always sent.
> Data Information Sequence	(0054,0063)	2	Sequence of Items that describe gating criteria. Contains only 1 item.
>> Frame Time	(0018,1063)	1	Nominal time per individual frame in msec

DIR DOC0660253 REV 15

>> Low R-R Value	(0018,1081)	3	R-R interval lower limit for beat rejection, in msec Always sent.
>> High R-R Value	(0018,1082)	3	R-R interval upper limit for beat rejection, in msec Always sent.
>> Intervals Acquired	(0018,1083)	3	Number of heartbeats that fall within Low R-R Value (0018,1081) and High R-R Value (0018,1082), and were therefore accepted and contribute gamma events to this R-R Interval. Always sent.
>> Intervals Rejected	(0018,1084)	3	Number of heartbeats that fall outside Low R-R (0018,1081) and High R-R Value (0018,1082), and do not contribute gamma events to this R-R Interval. Always sent.
>> Time Slot Information Sequence	(0054,0072)	2C	Sequence of Items that describe Time Slot Information. Sent if the Frame Increment Pointer (0028,0009) contains the Tag for Time Slot vector (0054,0070) Contains 1 or more items if it presents, the number of items is equal to Number of Time Slots (0054,0071).
>>> Time Slot Time	(0054,0073)	3	Not Used

3.4.6.12 NM Phase Module

This section contains Attributes that describe dynamic phases of a dynamic acquisition image performed on the patient. This module is present when the Image Type (0008,0008) Value 3, is equal to DYNAMIC. New phase may be defined by user from Scan setup. Default number of phases is 1.

TABLE 3-26  
NM PHASE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Phase Information Sequence	(0054,0032)	2C	Sequence of Items that describes each dynamic phase. Sent if the Frame Increment Pointer (0028,0009) contains the Tag for Phase Vector (0054,0030). Contains 1 to 5 items, if presents.
> Phase Delay	(0054,0036)	1	Time paused between the last frame of the previous phase and the first frame of this phase, in msec. (*) Set to 0 for 1 <sup>st</sup> phase.
> Actual Frame Duration	(0018,1242)	1	Nominal time of acquisition per individual frame, in msec.(*)
> Pause Between Frames	(0054,0038)	1	Time paused between each frame of this phase (in msec).(*)
> Number of Frames in Phase	(0054,0033)	1	Number of frames in this phase. (*)

**Note1:** (\*) – Attribute value is taken from user input

**3.4.6.13 VOI LUT Module**

**TABLE 3-27  
VOI LUT MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Window Center	(0028,1050)	1C	Window Center for display. Only single value is present. Calculated from actually acquired maximal and minimal pixel values. Always sent
Window Width	(0028,1051)	1C	Window Width for display. Only single value is present. Calculated from actually acquired maximal and minimal pixel values. Always sent.

**3.4.6.14 SOP Common Module**

**TABLE 3-28  
SOP COMMON MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
SOP Class UID	(0008,0016)	1	Uniquely identifies the SOP Class. Always set to "1.2.840.10008.5.1.4.1.1.20"
SOP Instance UID	(0008,0018)	1	Uniquely identifies the SOP Instance. Internally generated.
Specific Character Set	(0008,0005)	1C	Character Set that expands or replaces the Basic Graphic Set. Defined Terms include for locally created images: ISO_IR 100 = Latin Alphabet No. 1 For images created from MWL, the value is copied from value provided in MWL, if not empty, otherwise ISO_IR 100 is used. Always included into image
Instance Number	(0020,0013)	3	See 3.4.6.1 for more specialization

**3.4.6.15 Private Image Module**

**TABLE 3-29  
PRIVATE IMAGE MODULE ATTRIBUTES**

Attribute Name	Tag	Private Creator ID	Attribute Description
Origin	(0009,xx12)	QUASAR_INTERNAL_USE	The origin of the image. "isworklist" if scheduled in MWL, "regular" if locally scheduled.
Image Type	(0009,xx1B)	QUASAR_INTERNAL_USE	Image type string as passed in the scan request
Stop Reason	(0009,xx1D)	QUASAR_INTERNAL_USE	Defines condition that image was installed to database

DIR DOC0660253 REV 15

Send List Mode flag	(0009,xx24)	QUASAR_INTERNAL_USE	Defined if list mode data shall be sent to list mode data destination station Taken from user input for Discovery NM/CT 670 CZT implementation. Send as 0 for all other implementations.
Imageset Objects Count	(0037,xx90)	QUASAR_INTERNAL_USE	For Multi-Frame images (NM) - equal to value of tag (0028,0008) - Number of Frames. For CT Images contains number of DICOM objects in Imageset
CT Isocenter Distance From Floor	(0037,xx92)	QUASAR_INTERNAL_USE	Distance of CT Isocenter from floor, in mm. Shared by all images created by same Hybrid scan. Used for Optima NM/CT 640 implementation. Default value: 1040 mm
Source Translator	(0013,xx11)	GEMS_GENIE_1	Internal code of product DICOM implementation. Enumerated Value = 11.
Dataset Name	(0011,xx12)	GEMS_GENIE_1	List of dataset names
Acquisition Parent UID	(0011,xx31)	GEMS_GENIE_1	Shared by all images created by same Hybrid scan. In Discovery NM/CT 670* and Optima NM/CT 640 implementations is also used for correct registration. For more details about NM/CT registration refer to Frame Of Reference, Section 3.4.4.1 Not sent for QC Images, may be sent otherwise.
Bed Position	(0027,xx11)	APEX_PRIVATE	Linear position of table in cm

3.4.6.16 Private NM Image Module

TABLE 3-30  
PRIVATE NM IMAGE MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Rate Vector	(0009,xx01)	QUASAR_INTERNAL_USE	Rate for each frame
Count Vector	(0009,xx02)	QUASAR_INTERNAL_USE	Counts accumulated for each frame

DIR DOC0660253 REV 15

Time Vector	(0009,xx03)	QUASAR_INTERNAL_USE	Time for each frame. For TOMO and GATED TOMO images with type of detector Motion = ACQ DURING STEP contains values of M-View time (M-Time) calculated as Raw Time normalized according to scan motion acquisition time (C-Time) vs M-Time
Camera Shape	(0009,xx08)	QUASAR_INTERNAL_USE	Camera Shape L-mode flag Possible Values : 1-L mode 0-Other
Whole Body Spots	(0009,xx10)	QUASAR_INTERNAL_USE	Defines if WB Spots are created. Sent in images created by WB Scans. Possible values: 0 –in WHOLE BODY image created by WB Scan 1 – SPOTS images created by WB Scan
Scan Administration Plan	(0009,xx25)	QUASAR_INTERNAL_USE	List of Administration Events from Scan Administration Plan
WB From	(0009,xx61)	QUASAR_INTERNAL_USE	Bed position in mm where WHOLE BODY image is started. Sent when Image Type (0008,0008) value3 contains WHOLE BODY, otherwise is not sent. Used in all configurations except Brivo NM 615.
WB To	(0009,xx62)	QUASAR_INTERNAL_USE	Bed position in mm where WHOLE BODY image is ended. Sent when Image Type (0008,0008) value3 contains WHOLE BODY, otherwise is not sent. Used in all configurations except Brivo NM 615.
WB Image Start YOffset	(0009,xx66)	QUASAR_INTERNAL_USE	YOffset in pixels from beginning of WHOLEBODY Image, referenced to WB From (0009,xx61). Sent when Image Type (0008,0008) value3 contains WHOLE BODY, otherwise is not sent.
WB Image End YOffset	(0009,xx67)	QUASAR_INTERNAL_USE	YOffset in pixels from beginning of WHOLEBODY Image, referenced to WB To (0009,xx62). Sent when Image Type (0008,0008) value3 contains WHOLE BODY, otherwise is not sent.

DIR DOC0660253 REV 15

Yoffset	(0037,xx71)	QUASAR_INTERNAL_USE	<p>CT to NM Y difference</p> <p>In 800 Series*, Discovery NM/CT 670* and Optima NM/CT 640 implementations is also used for correct registration.</p> <p>For more details about NM/CT registration refer to Section 3.4.4.1.</p> <p>In Discovery NM 630 and Brivo NM 615 implementation is always sent as 0.</p>
Image Position Corrections Flag	(0037,xx72)	QUASAR_INTERNAL_USE	<p>Defined if additional image position corrections shall be performed.</p> <p>Always sent as “false”.</p>
NM Bed Position	(0037,xx73)	QUASAR_INTERNAL_USE	<ol style="list-style-type: none"> <li>1) In 800 Series*,Discovery NM/CT 670* implementations for Hybrid NM/CT Scans contains distance from NM FOV center to the position on bed where landmark was set (mm).</li> <li>2) In Optima NM/CT 640 implementation contains distance from NM FOV center to registration origin.</li> <li>3) In 800 Series*, Discovery NM/CT 670* and Optima NM/CT 640 implementations also used for correct registration For more details about NM/CT registration refer to Section 3.4.4.1.)</li> <li>4) For other scan types in 800 Series*, Discovery NM/CT 670* and Optima NM/CT 640 implementations contains linear bed position (in mm).</li> <li>5) In Discovery NM 630 and Brivo NM 615 implementations contains linear bed position (in mm)</li> </ol>
Xoffset	(0037,xx78)	QUASAR_INTERNAL_USE	<p>CT to NM X difference</p> <p>In 800 Series*, Discovery NM/CT 670* and Optima NM/CT 640 implementations is also used for correct registration.</p> <p>For more details about NM/CT registration refer to Section 3.4.4.1.</p> <p>In Discovery NM 630 and Brivo NM 615 implementations is always sent as 0.</p>
Collimator SQ	(0037,xx10)	QUASAR_INTERNAL_USE	<p>Contains information of collimators parameters. May contain 0 or 1 item.</p>
>Hole Diameter	(0037,xx1B)	QUASAR_INTERNAL_USE	<p>Collimator hole diameter</p>
>Hole Length	(0037,xx30)	QUASAR_INTERNAL_USE	<p>Collimator hole length</p>



DIR DOC0660253 REV 15

>Collimator Thickness	(0037,xx40)	QUASAR_INTERNAL_USE	Collimator thickness
>Septal Thickness	(0037,xx50)	QUASAR_INTERNAL_USE	Collimator septal thickness
>Intrinsic Resolution	(0037,xx60)	QUASAR_INTERNAL_USE	Collimator intrinsic resolution
>Blurring Slope	(0037,xx70)	QUASAR_INTERNAL_USE	Collimator blurring slope
Radio Nuclide Name	(0011,xx0D)	GEMS_GENIE_1	Name of radionuclide used Contain name of Energy Session, selected for acquisition
Detector Number	(0011,xx15)	GEMS_GENIE_1	Detector number image was acquired by. Sent when image is acquired using one detector only. For Brivo NM 615 Implementation WB Posterior image this field is set to "2", even though it is obtained by Detector-1. This is for SW compatibility purpose.

3.4.6.17 Private Image Tomo Module

TABLE 3-31  
PRIVATE IMAGE TOMO MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Angle Vector	(0009,xx07)	QUASAR_INTERNAL_USE	Angle for each TOMO frame. For each frame is tells what is the angle of the detector
Raw Time Vector	(0009,xx1A)	QUASAR_INTERNAL_USE	Raw time vector
M-View Normalization Factors Vectors	(0009,xx63)	QUASAR_INTERNAL_USE	Vector of normalization factors per M-View (M-View is a projection acquired during detector motion). Calculated as M-View Time (0009, xx03) vs. C-Time (0009,xx63). Send only for images where Type of Detector Motion (0054,0202) is ACQ DURING STEP, otherwise is not sent
C-Time	(0009,xx64)	QUASAR_INTERNAL_USE	Scan motion acquisition time in msec. Calculated as average of acquisition during motion time(M-Time) of all steps except the first step. Send for images where Type of Detector Motion (0054,0202) is ACQ DURING STEP
S&S Time	(0009,xx65)	QUASAR_INTERNAL_USE	The duration of acquiring step without motion in msec. Calculated as average of time vector (0009, xx03). Send for images where Type of Detector Motion (0054,0202) is ACQ DURING STEP
Send Joined Flag	(0009,xx23)	QUASAR_INTERNAL_USE	Defines if image shall be joined on send. Set to 1 – for TOMO and GATED TOMO images. For other image types is not sent.

DIR DOC0660253 REV 15

3.4.6.18 Private Image Multi-Gated Module

TABLE 3-32  
PRIVATE IMAGE MULTI-GATED MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Triggers Modification Flag	(0033,xx30)	GEMS_GENIE_1	Triggers Modification Flag
Number of triggers	(0033,xx33)	GEMS_GENIE_1	Number of triggers
Trigger size	(0033,xx34)	GEMS_GENIE_1	Size of one Trigger data slot
Trigger Data size	(0033,xx35)	GEMS_GENIE_1	Size of Trigger Data size
Trigger Data	(0033,xx36)	GEMS_GENIE_1	Buffer with trigger data information
Starting Heart Rate	(0009,xx37)	GEMS_GENIE_1	Heart rate at start of acquisition.

3.4.6.19 Private Image GSPECT Module

TABLE 3-33  
PRIVATE IMAGE GSPECT MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Average RR Time Vector	(0009,xx15)	QUASAR_INTERNAL_USE	Average r-r time vector
Low Limit Vector	(0009,xx16)	QUASAR_INTERNAL_USE	Low window limit vector
High Limit Vector	(0009,xx17)	QUASAR_INTERNAL_USE	High window limit vector
Begin Index Vector	(0009,xx18)	QUASAR_INTERNAL_USE	Begin index vector: link to heart beat vector
End Index Vector	(0009,xx19)	QUASAR_INTERNAL_USE	End index vector: link to heart beat vector
Perfusion SOP Instance UID	(0009,xx47)	QUASAR_INTERNAL_USE	Summed image UID

3.5 STANDARD EXTENDED AND PRIVATE DATA ATTRIBUTES

The Product supports the Standard Extended and Private Attributes defined in the following sections in Standard Extended NM SOP Instances as Type 3 data elements.

3.5.1 Standard Extended Attributes

The Product supports the following attributes, not specified in the NM IOD, in SOP Instances as Type 3 data elements.

TABLE 3-34  
STANDARD EXTENDED ATTRIBUTES

Information Entity Name	Attribute Name	Tag	Use
Study	Allergies	(0010,2110)	Description of prior reaction to contrast agents, or other patient allergies or adverse reactions.

DIR DOC0660253 REV 15

	Pregnancy Status	(0010,21C0)	Describes pregnancy state of patient Enumerated Values: 0001 = not pregnant 0002 = possibly pregnant 0003 = definitely pregnant 0004 = unknown
	Requested Procedure Comments	(0040,1400)	User-defined Study notes
Series	Patient Position	(0018,5100)	Patient position descriptor relative to the Equipment.

3.5.2 Private Group QUASAR\_INTERNAL\_USE

TABLE 3-35  
PRIVATE GROUP QUASAR\_INTERNAL\_USE

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0009,00xx)	LO	1	QUASAR_INTERNAL_USE
Rate Vector	(0009,xx01)	UL	1-n	Rate for each frame
Count Vector	(0009,xx02)	UL	1-n	Counts accumulated for each frame
Time Vector	(0009,xx03)	UL	1-n	Time for each frame
Angle Vector	(0009,xx07)	UL	1-n	Angle for each TOMO frame.
Camera Shape	(0009,xx08)	US	1	Camera Shape
Whole Body Spots	(0009,xx10)	US	1	Defines if WB Spots are created. Sent in images created by WB Scans.
Origin	(0009,xx12)	LO	1	The origin of the image.
Sequence Type	(0009,xx13)	ST	1	Acquired Sequence Type
Sequence Name	(0009,xx14)	ST	1	Acquired Sequence Name
Average RR Time Vector	(0009,xx15)	UL	1-n	Average r-r time vector
Low Limit Vector	(0009,xx16)	UL	1-n	Low window limit vector
High Limit Vector	(0009,xx17)	UL	1-n	High window limit vector
Begin Index Vector	(0009,xx18)	UL	1-n	begin index vector: link to heart beat vector
End Index Vector	(0009,xx19)	UL	1-n	end index vector: link to heart beat vector
Raw Time Vector	(0009,xx1A)	UL	1-n	Raw time vector
Image Type	(0009,xx1B)	LO	1	Image type string as passed in the scan request
Stop Reason	(0009,xx1D)	US	1	Defines condition that image was installed to db
Auto-Processing Application	(0009,xx1E)	ST	1	Auto-Processing Application Description
Acquisitions Count	(0009,xx1F)	US	1	Number of acquisitions per study
Completed Acquisitions Count	(0009,xx20)	US	1	Number of completed acquisitions per study
Send Joined Flag	(0009,xx23)	US	1	Defines if image shall be joined on send. Sent for TOMO and GATED TOMO images only

DIR DOC0660253 REV 15

Send List Mode Flag	(0009,xx24)	US	1	Defines if list mode data of acquired scan is sent to list mode data destination station
Administrations Plan	(0009,xx25)	LO	1..n	List of scan planned administrations' names.
Patient Unique Key	(0009,xx39)	UI	1	Patient unique key
Protocol Scheduled Date	(0009,xx40)	DA	1	Protocol Scheduled Date
Protocol Scheduled Time	(0009,xx41)	TM	1	Protocol Scheduled Time
Acquisition flag	(0009,xx42)	LO	1	Used for indicating if the study is acquired
Matched protocol	(0009,xx43)	LO	1	The originally matched protocol vs. protocol name which is the protocol actually acquired
Private SPS ID	(0009,xx44)	SH	1	Defines protocols that were appended to the original MWL protocol. Contains SPS ID of original MWL Protocol.
Pre-Medication	(0009,xx45)	LO	1	Keeps the Pre-Medication as appears in the "To Do" list
Perfusion SOP Instance UID	(0009,xx47)	UI	1	Summed image UID .Sent for GATED TOMO images only
Anatomic Reference	(0009,xx48)	LO	1	Keeps the anatomic reference for the specific scan
Acquisition Mode	(0009,xx59)	IS	1	Acquisition Mode.
WB From	(0009,xx61)	DS	1	Bed position in mm where WHOLEBODY image is started. Sent for WHOLEBODY images only.
WB To	(0009,xx62)	DS	1	Bed position in mm where WHOLEBODY image is ended. Sent for WHOLEBODY images only.
M-View Normalization Factors Vectors	(0009,xx63)	UL	1..n	Vector of normalization factors per M-View Sent for TOMO and GATED TOMO images only.
C-Time	(0009,xx64)	UL	1	Scan motion acquisition time. Sent for TOMO and GATED TOMO images only.
S&S Time	(0009,xx65)	UL	1	The duration of acquiring step without motion. Sent for TOMO and GATED TOMO images only.
WB Image Start YOffset	(0009,xx66)	UL	1	YOffset in pixels from beginning of WHOLEBODY Image, referenced to WB From (0009,xx61). Sent for WHOLEBODY images only.
WB Image End YOffset	(0009,xx67)	UL	1	YOffset in pixels from beginning of WHOLEBODY Image, referenced to WB To (0009,xx62). Sent for WHOLEBODY images only.
Private Creator Identification	(0037,00xx)	LO	1	QUASAR_INTERNAL_USE

DIR DOC0660253 REV 15

Collimator SQ	(0037,xx10)	SQ	1	Contains information of collimators parameters.
Hole Diameter	(0037,xx1B)	LO	1	Collimator hole diameter
Hole Length	(0037,xx30)	LO	1	Collimator hole length
Collimator Thickness	(0037,xx40)	LO	1	Collimator thickness
Septal Thickness	(0037,xx50)	LO	1	Collimator Septal thickness
Intrinsic Resolution	(0037,xx60)	LO	1	Collimator intrinsic resolution
Blurring Slope	(0037,xx70)	LO	1	Collimator blurring slope
Yoffset	(0037,xx71)	FD	1	CT to NM Y difference
Image Position Corrections Flag	(0037,xx72)	SH	1	Defined if image position corrections are performed.
NM Bed Position	(0037,xx73)	FD	1	According to implementation and scan types contains either <ul style="list-style-type: none"> <li>- Distance from NM FOV center to the position on bed where landmark was set .</li> <li>or</li> <li>- Distance from NM FOV center to the registration origin</li> <li>or</li> <li>- Linear bed position</li> </ul>
Xoffset	(0037,xx78)	FD	1	CT to NM X difference
Imageset Objects Count	(0037,xx90)	IS	1	Number of DICOM objects in Imageset for CT/ number of frames for NM.
CT Isocenter Distance from Floor	(0037,xx92)	DS	1	Distance of CT Isocenter Distance Floor, in mm

3.5.3 Private Group GEMS\_GENIE\_1

TABLE 3-36  
PRIVATE GROUP GEMS\_GENIE\_1

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0009,00xx)	LO	1	GEMS_GENIE_1
Starting Heart Rate	(0009,xx37)	SL	1	Heart rate at start of acquisition.
Private Creator Identification	(0011,00xx)	LO	1	GEMS_GENIE_1
Radio Nuclide Name	(0011,xx0D)	LO	1	Name of radionuclide used.
Dataset Name	(0011,xx12)	LO	1-n	List of Dataset names.
Detector Number	(0011,xx15)	SL	1-n	Sent when image is acquired using one detector only.
Acquisition Parent UID	(0011,xx31)	LO	1-n	Shared by all images created by same scan
Private Creator Identification	(0013,00xx)	LO	1	GEMS_GENIE_1

DIR DOC0660253 REV 15

Source Translator	(0013,xx11)	SL	1	Internal code of product DICOM implementation. Enumerated Value = 11.
Private Creator Identification	(0033,00xx)	LO	1	GEMS_GENIE_1
Triggers Modification Flag	(0033,xx30)	SL	1	Triggers Modification Flag
Number of triggers	(0033,xx33)	SL	1	Number of triggers
Trigger size	(0033,xx34)	SL	1	Size of one Trigger data slot
Trigger Data size	(0033,xx35)	SL	1	Size of Trigger Data size
Trigger Data	(0033,xx36)	OB	1	Buffer with trigger data information

3.5.4 Private Group GEMS\_XELPRV\_01

TABLE 3-37  
PRIVATE GROUP GEMS\_XELPRV\_01

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0033,00xx)	LO	1	GEMS_XELPRV_01
Object Type	(0033,xx08)	CS	1	Object Type. Contains string "SERIES DATA "
Modified Flag	(0033,xx10)	SL	1	Default value = 0 (Not Modified)
Name	(0033,xx11)	LO	1	SDO Name
Database Object Unique ID	(0033,xx16)	LO	1	Database UID of SDO; contains value of SDO UID tag (0033,xx72) generated at time of object creation
Date	(0033,xx17)	SH	1	SDO Creation date
Time	(0033,xx18)	SH	1	SDO Creation time
SeriesDataFlags	(0033,xx19)	UL	1	SDO Flags. Default value = 0
ProtocolName	(0033,xx1A)	LO	1	Name of Protocol created SDO
RelevantDataUID	(0033,xx1B)	LO	1	UID(s) of SOP Instance(s) relative to SDO
BulkData	(0033,xx1C)	OB	1	SDO parameter(s) stored as binary buffer(s)
IntData	(0033,xx1D)	SL	1-n	List of SDO parameters stored as integers
Double Data	(0033,xx1E)	FD	1-n	List of SDO parameters stored as doubles
String Data	(0033,xx1F)	OB	1	List of SDO parameters stored as list of strings
BulkDataFormat	(0033,xx20)	OB	1	Format of bulk parameters; contains information about name and size of bulk buffers
IntDataFormat	(0033,xx21)	OB	1	Format of integer parameters; contains information about name and number of integers in list
DoubleDataFormat	(0033,xx22)	OB	1	Format of double parameters; contains information about name and number of doubles in list

DIR DOC0660253 REV 15

StringDataFormat	(0033,xx23)	OB	1	Format of string parameters; contains information about name and number of strings in list
Description	(0033,xx24)	LT	1	User or equipment generated SDO description
Series Data Sequence	(0033,xx70)	SQ	1	Sequence of item contains information about acquisition parameters. May contain from 1 to n Items. Each Items describes specific parameters set.
SDO Private SOP Class UID	(0033,xx71)	UI	1	SDO Private SOP Class UID- "1.2.840.113619.4.17"
SDO Instance UID	(0033,xx72)	UI	1	SDO Instance UID; Internally generated

3.5.5 Private Group APEX\_PRIVATE

TABLE 3-38  
PRIVATE GROUP APEX\_PRIVATE

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0027,00xx)	LO	1	APEX_PRIVATE
Bed Position	(0027,xx11)	DS	1	Linear position of table.

## 4. SECONDARY CAPTURE INFORMATION OBJECT IMPLEMENTATION

### 4.1 INTRODUCTION

**Note:** In this section term **NM Camera AE DICOM implementation** refers to **DICOM implementation of NM General Purpose 600/800 Series products except Discovery NM/CT 670\* and NM/CT 800\* in CT standalone mode if not specially mentioned.**

This section specifies the use of the DICOM SC Image IOD to represent the information included in SC Images produced by this implementation. Corresponding attributes are conveyed using the module construct.

Screen Save images, created on the **NM General Purpose 600/800 Series** systems, are sent as DICOM Secondary Capture images in the following cases:

- As reports of several Quality Control Operations (Daily QC, Weekly QC, etc.), referred as term “QC SC” further in this section.
- As summary reports, complimentary for CT Radiation Dose SR (in Optima NM/CT 640 Implementation only). Refer to Section 9 for more details. This type of SC object is referred as term “CT Dose SR SC” further in this section.
- As summary reports, complimentary Radiopharmaceutical Radiation Dose Structure Reports (except Brivo NM 615 implementation). Refer to Section 10 for more details. This type of SC object is referred as term “RRDR SC” further in this section.

Only single frame Secondary Capture Image IOD is supported.

### 4.2 NM CAMERA MAPPING OF DICOM ENTITIES

NM Camera AE maps DICOM Information Entities to local Information Entities in the product’s database and user interface.

**TABLE 4-1  
MAPPING OF DICOM ENTITIES TO NM CAMERA ENTITIES**

DICOM IE	NM Camera Entity
Patient	Patient
Study	Exam
Series	Series
Image	Image



4.3 IOD MODULE TABLE

The Secondary Capture Information Object Definition comprises the modules of the following table, plus Standard Extended and Private attributes. Standard Extended and Private attributes are described in Section 4.5.

TABLE 4-2  
SC IMAGE IOD MODULES

Entity Name	Module Name	Usage	Reference
Patient	Patient	Used for QC SC	4.4.1.1
		For CT Dose SR SC and RRDR SC (same description as for NM IOD)	3.4.1.1
	Clinical Trial Subject	Not Used	N/A
Study	General Study	Used for QC SC	4.4.2.1
		For CT Dose SR SC and RRDR SC (same description as for NM IOD)	3.4.2.1
	Patient Study	Used for CT Dose SR SC and RRDR SC (same description as for NM IOD)	3.4.2.2
	Standard Extended Study	Used for CT Dose SR SC and RRDR SC (same description as for NM IOD)	3.4.2.3
	Private Study	Used for CT Dose SR SC (same description as for NM IOD)	3.4.2.4
	Clinical Trial Study	Not Used	N/A
Series	General Series	Used	4.4.3.1
	Private SC Series	Used for CT Dose SR SC and RRDR SC	4.4.3.2
	Clinical Trial Series	Not Used	N/A
Equipment	General Equipment	Used (same description as for NM IOD)	3.4.5.1
	SC Equipment	Used	4.4.4.1
Image	General Image	Used	4.4.5.1
	Image Pixel	Used	4.4.5.2
	Device	Not Used	N/A
	Specimen	Not Used	N/A
	SC Image	Used	4.4.5.3
	Overlay Plane	Not Used	N/A
	Modality LUT	Not Used	N/A
	VOI LUT	Not Used	N/A
	SOP Common	Used	4.4.5.4
	Standard Extended Image	Used	4.4.5.5
	Private SC Image	Used	4.4.5.5

4.4 INFORMATION MODULE DEFINITIONS

Please refer to DICOM Standard Part 3 (Information Object Definitions) for a description of each of the entities and modules contained within the SC Information Object.

The following modules are included to convey Enumerated Values, Defined Terms, and Optional Attributes supported. Type 1 & Type 2 Attributes are also included for completeness and to define what values they may take and where these values are obtained from when generating the instance. It should be noted that they are the same ones as defined in the DICOM Standard Part 3 (Information Object Definitions). Also note that Attributes not present in tables are not supported.

4.4.1 Patient Entity Modules

4.4.1.1 Patient Module

TABLE 4-3  
PATIENT MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Patient's Name	(0010,0010)	2	For CT Dose SR SC and RRDR SC: Patient's full name. (*) Compound from user Last Name and First Name for locally scheduled protocols(**) (***) For QC SC: name of user that performed the QC protocol. For example, "Service^^^^ "
Patient ID	(0010,0020)	2	For CT Dose SR SC and RRDR SC: primary hospital identification number or code for the patient. (*) (**) (***) For QC SC: set to Camera Name from system configuration
Issuer of Patient ID	(0010,0021)	3	Not Used
Issuer of Patient ID Qualifiers Sequence	(0010,0024)	3	Not Used
Patient's Birth Date	(0010,0030)	2	For CT Dose SR SC and RRDR SC: Birth date of the patient. (*) (**) (***) For QC SC: execution date of the QC Test Sent as ZERO LENGTH value if value is not received either from MWL or from user input.
Patient's Sex	(0010,0040)	2	For CT Dose SR SC and RRDR SC: Sex of the named patient. (*) (**)(***) Enumerated Values: M = male F = female O = other Sent as ZERO LENGTH value if value is not received either from MWL or from user input. For QC SC: always sent as empty string

DIR DOC0660253 REV 15

Other Patient IDs	(0010,1000)	3	Not Used
Other Patient IDs Sequence	(0010,1002)	3	Not Used

**Note 1 :** (\*) - Attributes copied from the Worklist if the study source was actually copied from a Worklist query result.

**Note 2 :** (\*\*) - Attributes copied from the user input for Locally scheduled Protocols

**Note 3 :** (\*\*\*) - Cannot be modified by user if meaningful value is received from MWL

#### 4.4.2 Study Entity Modules

##### 4.4.2.1 General Study Module

**TABLE 4-4  
GENERAL STUDY MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Study Instance UID	(0020,000D)	1	Unique identifier for the Study. Internally generated by the system.
Study Date	(0008,0020)	2	Date the Study started.
Study Time	(0008,0030)	2	Time the Study started
Referring Physician's Name	(0008,0090)	2	Name of the patient's referring physician. Sent as empty string.
Study ID	(0020,0010)	2	Equipment generated Study identifier. Automatically assigned to default Value "QC-protocol"
Accession Number	(0008,0050)	2	Accession Number. Sent as an empty string
Study Description	(0008,1030)	3	Study Description. Automatically assigned to the description of the QC test performed. For example, "NM Daily QC", "Uniformity Test", "COR Test" etc.  Always sent as non-empty value

#### 4.4.3 Series Entity Modules

##### 4.4.3.1 General Series Module

**TABLE 4-5  
GENERAL SERIES MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Modality	(0008,0060)	1	SC Images created by NM General Purpose 600/800 Series products generally have this attribute set to the value found in the original image.  Defined Terms: <ul style="list-style-type: none"> <li>NM = Nuclear Medicine (QC SC and RRDR SC)</li> <li>CT= Computerized Tomography (CT Dose SR SC)</li> </ul>

Series Instance UID	(0020,000E)	1	Internally generated unique identifier of the Series.
Series Number	(0020,0011)	2	A number that identifies this Series. For QC SC : set as ZERO LENGTH value for For CT Dose SR SC :Set as random integer value For RRDR SC: Internally generated using default Base Value = 500 Calculated as <Base Value>+<N+1>, where <N> is number of RRDR SC objects already exists in Study
Protocol name	(0018,1030)	3	User-defined description of the conditions under which the Series was performed. For CT Dose SC: contains acquisition protocol Name For RRDR SC: always sent as “Exposure Monitoring”
Series Date	(0008,0021)	3	Date the Series started Always sent
Series Time	(0008,0031)	3	Time the Series started. Always sent
Series Description	(0008,103E)	3	Description of the Series. Always sent as non-empty value. For QC SC: defined by the name of QC report. For CT Dose SR SC : Composed from parent CT Series Name and Scan Type followed by suffix “Dose Secondary Capture” For RRDR SC: Composed from Administration Event name followed by suffix “ Secondary Capture”

4.4.3.2 Private SC Series

TABLE 4-6  
PRIVATE SC SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Origin	(0009,xx12)	QUASAR_INTERNAL_USE	The origin of the image. For CT Dose SR: Acquisition protocol name For RRDR SC: Always sent as “exposure monitoring”
Send List Mode flag	(0009, xx24)	QUASAR_INTERNAL_USE	Defined if list mode data shall be sent to list mode data destination station Always sent as 0.

DIR DOC0660253 REV 15

Sequence Type	(0009,xx13)	QUASAR_INTERNAL_USE	Type of sequence created the object Always send as "secondary capture"
Sequence Name	(0009,xx14)	QUASAR_INTERNAL_USE	Name of Sequence created the object. Always sent as non-empty value
Pre-Medication	(0009,xx45)	QUASAR_INTERNAL_USE	Keeps the Pre-Medication as appears in the "To-Do" list  Sent as empty string if not received from user input.

4.4.4 Equipment Entity Modules

4.4.4.1 SC Equipment Module

TABLE 4-7  
SC EQUIPMENT MODULE ATTRIBUTES

Attribute Name	Tag	Type	Use
Conversion Type	(0008,0064)	1	Specify Defined Terms used: SI = Scanned Image
Modality	(0008,0060)	3	SC Images created by NM General Purpose 600/800 Series products generally have this attribute set to the value found in the original image.  Always sent.  Defined Terms: <ul style="list-style-type: none"> <li>• NM = Nuclear Medicine (QC SC and RRDR SC)</li> <li>• CT= Computerized Tomography (CT Dose SR SC)</li> </ul>

4.4.5 Image Entity Modules

4.4.5.1 General Image Module

TABLE 4-8  
GENERAL IMAGE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Instance Number	(0020,0013)	2	A number that identifies this image. Sent as ZERO Length value
Content Date	(0008,0023)	2C	The date the SC Image pixel data creation started
Content Time	(0008,0033)	2C	The time the SC image pixel data creation started

Quality Control Image	(0028,0300)	3	Indicates that this image is a quality control image. Enumerated Values: YES - for QC SC NO - for CT Dose SC and RRDR SC
Derivation Description	(0008,2111)	3	A text description of how this image was derived. For QC SC: not sent For CT Dose SC: Always set to "X-Ray Radiation Dose Report" For RRDR SC: Always set to "Administration Secondary Capture"
Image Type	(0008,0008)	3	The following Enumerated Value of Value 1 is supported: <ul style="list-style-type: none"> <li>DERIVED - identifies a Derived Image</li> </ul> The following Enumerated Value of Value 2 is supported: <ul style="list-style-type: none"> <li>SECONDARY - identifies a Secondary Image</li> </ul>

4.4.5.2 Image Pixel Module

TABLE 4-9  
IMAGE PIXEL MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Samples per Pixel	(0028,0002)	1	Number of samples (planes) in this image. Possible values: 1, if Photometric Interpretation (0028,0004) equals "MONOCHROME2" 3, if Photometric Interpretation (0028,0004) equals "RGB"
Photometric Interpretation	(0028,0004)	1	Specifies the intended interpretation of the pixel data Defined Terms supported: <ul style="list-style-type: none"> <li>MONOCHROME2</li> <li>RGB (for QC SC)</li> </ul>
Rows	(0028,0010)	1	Number of rows in the image. Possible values: QC SC: 878, 934 CT Dose SC and RRDR SC: 300
Columns	(0028,0011)	1	Number of columns in the image Possible values: QC SC: 1004 CT Dose SC: 520 RRDR SC: 620

Bits Allocated	(0028,0100)	1	Number of bits allocated for each pixel sample. Each sample shall have the same number of bits allocated. Enumerated Values supported : QC SC: 8 CT Dose SC and RRDR SC: 16
Bits Stored	(0028,0101)	1	Number of bits stored for each pixel sample. Value equal to Bit Allocated (0028,0100)
High Bit	(0028,0102)	1	Most significant bit for pixel sample data. Value equal to Bit Stored (0028,0101) - 1
Pixel Representation	(0028,0103)	1	Data representation of the pixel samples. Each sample shall have the same pixel representation. Enumerated Values: 0000H = unsigned integer.
Pixel Data	(7FE0,0010)	1	A data stream of the pixel samples that comprise the Image.
Planar Configuration	(0028,0006)	1C	Sent for QC SC only Enumerated Values: 0000H = color-by-pixel

4.4.5.3 SC Image Module

TABLE 4-10  
SC IMAGE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Use
Pixel Spacing	(0028,0030)	1C	Not sent. Secondary Capture images created by product are not calibrated images, contain screen captures of QC result tables and Dose report information

4.4.5.4 SOP Common Module

TABLE 4-11  
SOP COMMON MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
SOP Class UID	(0008,0016)	1	Uniquely identifies the SOP Class. Always set to "1.2.840.10008.5.1.4.1.1.7"
SOP Instance UID	(0008,0018)	1	Uniquely identifies the SOP Instance. Internally generated.
Specific Character Set	(0008,0005)	1C	Character Set that expands or replaces the Basic Graphic Set. Defined Terms include for locally created images: ISO_IR 100 = Latin Alphabet No. 1 Always sent.
Instance Number	(0020,0013)	3	Sent as empty String

DIR DOC0660253 REV 15

4.4.5.5 Standard Extended Image Module

TABLE 4-12  
IMAGE STANDARD EXTENDED ATTRIBUTES

Attribute Name	Tag	Use
Counts Accumulated	(0018,0070)	Always sent as "1"

4.4.5.6 Private SC Image Module

TABLE 4-13  
PRIVATE SC MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Sequence Name	(0009,xx14)	"QUASAR_INTERNAL_USE"	Acquired Sequence Type
Patient Unique Key	(0009,xx39)	"QUASAR_INTERNAL_USE"	Patient unique key
Dataset Name	(0011,xx12)	"GEMS_GENIE_1"	List of Dataset names

4.5 STANDARD EXTENDED AND PRIVATE DATA ATTRIBUTES

The Product supports the Standard and Private Attributes defined in the following sections in Standard Extended SC SOP Instances as Type 3 data elements.

4.5.1 Standard Extended Attributes

The Product supports the following attributes, not specified in the SC IOD, in SOP Instances as Type 3 data elements.

TABLE 4-14  
STANDARD EXTENDED ATTRIBUTES

Information Entity Name	Attribute Name	Tag	Use
Image	Counts Accumulated	(0018,0070)	

4.5.2 Private Group QUASAR\_INTERNAL\_USE

Private Group QUASAR\_INTERNAL\_USE is modeled as part of the Image Information Entity.

TABLE 4-15  
PRIVATE GROUP QUASAR\_INTERNAL\_USE

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0009,00xx)	LO	1	QUASAR_INTERNAL_USE
Origin	(0009,xx12)	LO	1	The origin of the image.
Sequence Type	(0009,xx13)	ST	1	Acquired Sequence Type
Sequence Name	(0009,xx14)	ST	1	Acquired Sequence Name
Send List Mode Flag	(0009,xx24)	US	1	Defines if list mode data of acquired scan is sent to list mode data destination station
Patient Unique Key	(0009,xx39)	UI	1	Patient unique key



DIR DOC0660253 REV 15

Pre-Medication	(0009,xx45)	LO	1	Keeps the Pre-Medication as appears in the "To Do" list
----------------	-------------	----	---	---

**4.5.3 Private Group GEMS\_GENIE\_1**

Private Group GEMS\_GENIE\_1 is modeled as part of the Image Information Entity.

**TABLE 4-16  
PRIVATE GROUP GEMS\_GENIE\_1**

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0011,00xx)	LO	1	GEMS_GENIE_1
Dataset Name	(0011,xx12)	LO	1-n	List of Dataset names.

**4.6 STANDARD EXTENDED AND PRIVATE CONTEXT GROUPS**

NM Camera AE does not support any coded terminology.

## 5. MODALITY WORKLIST QUERY IMPLEMENTATION

### 5.1 INTRODUCTION

The Discovery NM/CT 670\* and NM/CT 800\* implementations (for NM standalone exams and NM/CT hybrid exams) , Discovery NM 630, NM 870, Brivo NM 615 and Optima NM/CT 640 implementations map DICOM Information Entities to local Information Entities in the product’s database and user interface.

In the Discovery NM/CT 670\*and NM/CT 800\* implementations for CT standalone exams, Worklist query is performed via the CT scanner console. For Modality Worklist Information Model Definition, refer to appropriate CT Scanner DICOM Conformance Statements (Section 1.6 )

**Note:** In this section term **NM Camera AE DICOM implementation** refers to **DICOM implementation of NM General Purpose 600/800 Series products except Discovery NM/CT 670\* and NM/CT 800\*in CT standalone mode.**

This section specifies the use of the DICOM Modality Worklist Information Model used to organize data and against which a Modality Worklist Query will be performed.

### 5.2 NM CAMERA MAPPING OF DICOM ENTITIES

The NM Camera AE maps DICOM Information Entities to local Information Entities in the product’s database and user interface.

**TABLE 5-1**  
**MAPPING OF DICOM ENTITIES TO NM CAMERA ENTITIES**

DICOM	NM Camera Entity
Scheduled Procedure Step	Protocol
Requested Procedure	Study
Imaging Service Request	Study
Visit	Study
Patient	Patient

Matching Requested Procedure Step to NM Camera protocol is done according to predefined configuration.

The configuration contains the following tags

- (0040,0007) – Scheduled Procedure Step Description
- (0032,1060) – Requested Procedure Description
- (0040,0008) – Scheduled Protocol Code Sequence - Code Meaning

The default configuration is (0040,0007) – Scheduled Procedure Step Description

If Scheduled Protocol Code Sequence - tag (0040,0008) - is selected for mapping, the protocol will be mapped according to the value of Code Meaning - tag (0008,0104).

**5.3 WORKLIST QUERY MODULE TABLE**

See DICOM PS 3.3 and PS 3.4 for a complete definition of the entities, modules, and attributes.

**TABLE 5-2  
MODALITY WORKLIST INFORMATION MODEL MODULES**

Entity Name	Module Name	Reference
Scheduled Procedure Step	SOP Common	5.4.1.1
	Scheduled Procedure Step	5.4.1.2
Requested Procedure	Requested Procedure	5.4.2.1
Imaging Service Request	Imaging Service Request	5.4.3.1
Visit	Visit Identification	5.4.4.1
	Visit Status	5.4.4.2
	Visit Relationship	5.4.4.3
	Visit Admission	5.4.4.4
Patient	Patient Relationship	5.4.5.1
	Patient Identification	5.4.5.2
	Patient Demographic	5.4.5.3
	Patient Medical	5.4.5.4

**5.4 WORKLIST QUERY MODULE DEFINITIONS**

Please refer to DICOM Standard PS 3.3. (Information Object Definitions) for a description of each of the query key attributes contained within the Modality Worklist Information Model.

Note that in all tables below information in “**Mapped into Instance/MPPS**” column is referenced to NM images and CT Images only if not explicitly specified. .

**5.4.1 Common Scheduled Procedure Step Entity Modules**

**5.4.1.1 SOP Common Module**

**TABLE 5-3  
SOP COMMON MODULE ATTRIBUTES**

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance/MPPS	Note
Specific Character Set	(0008,0005)	O	1C	Yes /Yes	See 5.4.1.1.1

**5.4.1.1.1 Specific Character Set**

The NM Scanner AE will use any Specific Character Set value returned in a Scheduled Procedure Step Identifier in the images created pursuant to that Scheduled Procedure Step. Text attributes, including Patient and Physician names, that include non-ASCII characters will be displayed as described in Section 2.7

For MPPS created from MWL request, the value is copied from value provided in MWL if not empty, otherwise ISO\_IR 100 is used.

5.4.1.2 Scheduled Procedure Step Module

TABLE 5-4  
SCHEDULED PROCEDURE STEP MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance / MPPS	Note
Scheduled Procedure Step Sequence	(0040,0100)	R	1	No / No	Only one Item is supported
>Scheduled Station AE Title	(0040,0001)	R	1	No / No	Single Value matching is used.
>Scheduled Procedure Step Start Date	(0040,0002)	R	1 *	Yes / No	Specified as range of date of the form: from Date - To Date. SPS Start Date is mapped into private attribute "Protocol Scheduled Date" in the image – Tag (0009,xx40), QUASAR_INTERNAL_USE Cannot be modified in UI if received from MWL
>Scheduled Procedure Step Start Time	(0040,0003)	R	1 *	Yes / No	Matching is not supported. SPS Start Time is mapped into private attribute "Protocol Scheduled Time" in the image - tag (0009,xx41), QUASAR_INTERNAL_USE Cannot be modified in UI if received from MWL.
>Modality	(0008,0060)	R	1 *	No / No	Single value matching is allowed. Possible Values: NM, CT
>Scheduled Performing Physician's Name	(0040,0006)	R	2	Yes / No	Wildcard matching is allowed by Last Name and First Name separately. User may enter matching values for Last Name and/or First Name separately in UI. The value sent in MWL request is created according to the following template : < Last Name>^< First Name> Mapped to Performing Physicians' Name tag (0008, 1050) for NM , CT and CT Dose SR IODs.
>Scheduled Procedure Step Description	(0040,0007)	O	1C *	Yes / Yes	Always included in the MWL request. May be used for Protocol Mapping.
>Scheduled Protocol Code Sequence	(0040,0008)	O	1C	Yes / Yes	Always included in the MWL request. MPPS: Copied from the first valid item of MWL request.

>>Code Value	(0008,0100)	O	1	Yes / Yes	Requested explicitly MPPS: Copied from the first valid item of MWL request
>>Coding Scheme Designator	(0008,0102)	O	1	Yes/ Yes	Requested explicitly MPPS: Copied from the first valid item of MWL request
>>Coding Scheme Version	(0008,0103)	O	3	Yes/Yes	Requested explicitly MPPS: Copied from the first valid item of MWL request
>>Code Meaning	(0008,0104)	O	3*	Yes/ Yes	Requested explicitly . Returned non-empty value may be used for Protocol Mapping If Protocol Mapping is configured to use this tag: <ul style="list-style-type: none"> <li>• Never displayed and stored in DB for mapped protocols.</li> <li>• If protocol is not mapped - displayed in "Scheduled Study" column in " To Do" List;</li> </ul> MPPS: Copied from the first valid item of MWL request
>Pre-Medication	(0040,0012)	O	2C*	Yes / No	Always included in the MWL request. Mapped into private attribute Pre-Medication – tag (0009,xx45) "QUASAR_INTERNAL_USE"
>Scheduled Procedure Step ID	(0040,0009)	O	1	Yes / Yes	Always included in the MWL request.
>Requested Contrast Agent	(0032,1070)	O	2C	No / No	Always included in the MWL request.
>Scheduled Procedure Step Status	(0040,0020)	O	3	No / No	Always included in the MWL request.

**Note:** \* in the *Expected Return Key Type* column indicates that this information is displayed on screen, if available

5.4.2 Common Requested Procedure Entity Modules

5.4.2.1 Requested Procedure Module

TABLE 5-5  
REQUESTED PROCEDURE MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance / MPPS	Note
----------------	-----	----------------------------	----------------------------	-----------------------------	------

Requested Procedure ID	(0040,1001)	O	1 *	Yes / Yes	Single value matching is allowed. Displayed in the "To Do" list.
Requested Procedure Description	(0032,1060)	O	1 *	Yes / Yes	Always included in the MWL request. Matching is not supported . May be used for Protocol Mapping. If Protocol Mapping is configured to use this tag: <ul style="list-style-type: none"> <li>• If protocol is not mapped - displayed in "Scheduled Study" column in "To Do" List.</li> <li>• Mapped to Study Description tag (0008,1030), if belongs to first Protocol in the Study</li> </ul>
Requested Procedure Code Sequence	(0032,1064)	O	1	Yes /Yes	Always included in the MWL request as ZERO Length Sequence Copied to Procedure Code Sequence (0008,1032) when stored to instance
>Code Value	(0008,0100)	O	1	Yes / Yes	Not Requested explicitly Copied to Procedure Code Sequence (0008,1032) when stored to instance
>Coding Scheme Designator	(0008,0102)	O	1	Yes /Yes	Not Requested explicitly Copied to Procedure Code Sequence (0008,1032) when stored to instance
>Coding Scheme Version	(0008,0103)	O	3	Yes / Yes	Not Requested explicitly Copied to Procedure Code Sequence (0008,1032) when stored to instance
>Code Meaning	(0008,0104)	O	3	Yes / Yes	Not Requested explicitly Copied to Procedure Code Sequence (0008,1032) when stored to instance
Study Instance UID	(0020,000D)	O	1	Yes / Yes	Always included in the MWL request.
Referenced Study Sequence	(0008,1110)	O	2	Yes / Yes	Always included in the MWL request as ZERO Length Sequence
>Referenced SOP Class UID	(0008,1150)	O	1C	Yes / Yes	Not Requested explicitly
>Referenced SOP Instance UID	(0008,1155)	O	1C	Yes / Yes	Not Requested explicitly

**Note:** \* in the *Expected Return Key Type* column indicates that this information is displayed on screen, if available

DIR DOC0660253 REV 15

5.4.3 Common Imaging Service Request Entity Modules

5.4.3.1 Imaging Service Request Module

TABLE 5-6  
IMAGING SERVICE REQUEST MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance / MPPS	Note
Accession Number	(0008,0050)	O	2 *	Yes / Yes	Single Value and Wildchar matching may be requested for this data element
Requesting Physician	(0032,1032)	O	2	No / No	Always included in the MWL request.
Referring Physician's Name	(0008,0090)	O	2 *	Yes / No	Always included in the MWL request Only First Name and Last Name are displayed on screen and stored in image; User can modify value arrived from MWL.
Placer Order Number / Imaging Service Request	(0040,2016)	O	3	No / No	Always included in the MWL request.
Filler Order Number / Imaging Service Request	(0040,2017)	O	3	No / No	Always included in the MWL request.

**Note:** \* in the *Expected Return Key Type* column indicates that this information is displayed on screen, if available

5.4.4 Common visit Entity Modules

5.4.4.1 Visit Identification

TABLE 5-7  
VISIT IDENTIFICATION MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance / MPPS	Note
Admission ID	(0038,0010)	O	2	No / No	Always included in the MWL request
Institution Name	(0008,0080)	O	3	No / No	Always included in the MWL request . Value from MWL response is not stored in the image – used value from System configuration instead.

5.4.4.2 Visit Status

TABLE 5-8  
VISIT STATUS MODULE ATTRIBUTES

DIR DOC0660253 REV 15

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance	Note
Current Patient Location	(0038,0300)	O	2	No / No	Always included in the MWL request

5.4.4.3 Visit Relationship

TABLE 5-9  
VISIT RELATIONSHIP MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance / MPPS	Note
Referenced Patient Sequence	(0008,1120)	O	2	No / Yes	

5.4.4.4 Visit Admission

TABLE 5-10  
VISIT ADMISSION MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance/ MPPS	Note
Referring Physician's Name	(0008,0090)	O	3	Yes / No	Always included in the MWL request Only First Name and Last Name are displayed on screen and stored in image; User can modify value arrived from MWL.

5.4.5 Common Patient Entity Modules

5.4.5.1 Patient Relationship

TABLE 5-11  
PATIENT RELATIONSHIP MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance / MPPS	Note
Referenced Visit Sequence	(0008,1125)	O	3	No / No	Always included in MWL request
Referenced Patient Alias Sequence	(0038,0004)	O	3	No / No	Always included in MWL request



5.4.5.2 Patient Identification

TABLE 5-12  
PATIENT IDENTIFICATION MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance / MPPS	Note
Patient's Name	(0010,0010)	R	1 *	Yes / Yes	Single Value or Wildchar matching is allowed Last name and First Name separately; Only First Name and Last Name are displayed on screen and stored into image. User cannot modify value received from MWL
Patient ID	(0010,0020)	R	1 *	Yes / Yes	Only Single Value matching is allowed for this data element. User cannot modify value received from MWL
Issuer of Patient ID	(0010,0021)	O	3	No / No	Always included in the MWL request
Other Patient IDs Sequence	(0010,1002)	O	3	No / No	Not Used
Issuer of Patient ID Qualifiers Sequence	(0010,0024)	O	3	No / No	Not Used
Other Patient IDs	(0010,1000)	O	3	No / No	Always included in the MWL request

**Note:** \* in the *Expected Return Key Type* column indicates that this information is displayed on screen, if available

5.4.5.3 Patient Demographic

TABLE 5-13  
PATIENT DEMOGRAPHIC MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance/ MPPS	Note
Patients Birth Date	(0010,0030)	O	2 *	Yes / Yes	Always included in the MWL request User cannot modify value received from MWL
Patient's Sex	(0010,0040)	O	2 *	Yes / Yes	Always included in the MWL request User cannot modify value received from MWL
Patient's Weight	(0010,1030)	O	2*	Yes / No	Always included in the MWL request User can modify value received from MWL

DIR DOC0660253 REV 15

Patient's Size	(0010,1020)	O	2*	Yes / No	Always included in the MWL request User can modify value received from MWL
Patient's Age	(0010,1010)	O	3	No / No	Always included in the MWL request Value calculated from Patient Birth Date is displayed on screen stored into image.
Confidentiality Constraint on Patient Data Description	(0040,3001)	O	3	No / No	Always included in the MWL request

**Note:** \* in the *Expected Return Key Type* column indicates that this information is displayed on screen, if available

5.4.5.4 Patient Medical

TABLE 5-14  
PATIENT MEDICAL MODULE ATTRIBUTES

Attribute Name	Tag	Expected Matching Key Type	Expected Returned Key Type	Mapped into Instance / MPPS	Note
Patient State	(0038,0500)	O	2	No / No	Always included in the MWL request
Pregnancy Status	(0010,21C0)	O	2*	Yes / No	Always included in the MWL request
Medical Alerts	(0010,2000)	O	2	No / No	Always included in the MWL request
Contrast Allergies	(0010,2110)	O	2*	Yes / No	Always included in the MWL request
Special Needs	(0038,0050)	O	2	No / No	Always included in the MWL request
Additional Patient History	(0010,21B0)	O	3	No / No	Always included in the MWL request

**Note:** \* in the *Expected Return Key Type* column indicates that this information is displayed on screen, if available

## 6. STORAGE COMMITMENT PUSH MODEL IMPLEMENTATION

### 6.1 STORAGE COMMITMENT PUSH MODEL INFORMATION OBJECT DEFINITION

**Note:** In this section term NM Camera AE DICOM implementation refers to DICOM implementation of NM General Purpose 600/800 Series products except Discovery NM/CT 670\* and NM/CT 800\*in CT standalone mode.

This section describes NM Camera AE Storage Commitment Push Model Implementation.

Please refer to DICOM Part 3 (Information Object Definitions) for a description of each of the attributes contained within the Storage Commitment Information Object.

The Storage Commitment Information Object is used both for N-ACTION Storage Commitment Requests by the SCU and N-EVENT-REPORT Storage Commitment Notifications by the SCP.

Note that requesting storage commitment for CT data is performed via the CT scanner console. Refer to appropriate CT Scanner DICOM Conformance Statements (Section 1.6).

#### 6.1.1 STORAGE COMMITMENT MODULE FOR N-ACTION

**TABLE 6-1  
STORAGE COMMITMENT MODULE FOR N-ACTION**

Attribute Name	Tag	SCU Use
Transaction UID	(0008,1195)	Internally generated.
Storage Media File-Set ID	(0088,0130)	Not used
Storage Media File-Set UID	(0088,0140)	Not used
Referenced SOP Sequence	(0008,1199)	May contain 1 or more items
>Referenced SOP Class UID	(0008,1150)	Storage SOP classes supported as SCU: 1.2.840.10008.5.1.4.1.1.20 Nuclear Medicine Image Storage SOP Class UID. 1.2.840.10008.5.1.4.1.1.7 - Secondary Capture Image Storage SOP Class UID
>Referenced SOP Instance UID	(0008,1155)	SOP Instance UID of the Image which Storage Commitment is required for.
>Storage Media File-Set ID	(0088,0130)	Not used
>Storage Media File-Set UID	(0088,0140)	Not used

6.1.2 STORAGE COMMITMENT MODULE FOR N-EVENT-REPORT

TABLE 6-2  
STORAGE COMMITMENT MODULE FOR N-EVENT-REPORT

Attribute Name	Tag	SCU Use
Transaction UID	(0008,1195)	Used to identify the N-ACTION Request which N-EVENT-REPORT is relevant to.
Retrieve AE Title	(0008,0054)	Not Used
Storage Media File-Set ID	(0088,0130)	Not used
Storage Media File-Set UID	(0088,0140)	Not used
Referenced SOP Sequence	(0008,1199)	Used to identify the images which storage commitment was successful and mark them as Archived.
>Referenced SOP Class UID	(0008,1150)	
>Referenced SOP Instance UID	(0008,1155)	
>Retrieve AE Title	(0008,0054)	Not Used
>Storage Media File-Set ID	(0088,0130)	Not used
>Storage Media File-Set UID	(0088,0140)	Not used
Failed SOP Sequence	(0008,1198)	Used to identify the images which storage commitment was failed to prevent marking them as Archived.
>Referenced SOP Class UID	(0008,1150)	
>Referenced SOP Instance UID	(0008,1155)	
>Failure Reason	(0008,1197)	See Section 6.1.2.1 for the list of processed values.

6.1.2.1 Processing of Failure Reason when received in a N-Event-Report

When receiving a N-Event-Report request with a Event Type ID equal to 2, meaning that Storage Commitment is complete, but failure exists, following is the set of value that this Storage Commitment SCU AE is able to process:

Failure Reason	Meaning	Application Behavior When Receiving Reason Code
0110H	Processing failure	Transfer failure is logged (*)
0112H	No such object instance	Transfer failure is logged (*)
0213H	Resource limitation	Transfer failure is logged (*)
0122H	Referenced SOP Class not supported	Transfer failure is logged (*)
0119H	Class / Instance conflict	Transfer failure is logged (*)
0131H	Duplicate transaction UID	Transfer failure is logged (*)
*	Other Failure Reason code values.	Transfer failure is logged (*)

**Note:** (\*) - In all failure reasons an Appropriate error message is logged in `/home/ctuser/neuvo/logfiles/nwscp.log`. If the error message contains failed sop instance uid's those will logged in the log file. The image/series/studies will not be marked as archived in the Data Management Panel.

## 7. MODALITY PERFORMED PROCEDURE STEP IMPLEMENTATION

### 7.1 INTRODUCTION

**Note:** In this section term **NM Camera AE DICOM implementation** refers to **DICOM implementation of NM General Purpose 600/800 Series products except Discovery NM/CT 670\* and NM/CT 800\* in CT standalone mode.**

This section specifies the use of the DICOM Modality Performed Procedure Step information to be communicated to the Hospital/Radiology information system.

This feature works in conjunction with DICOM Modality Worklist feature, if installed. However the conformance of this feature is independent of Modality Worklist feature. For information on conformance of Modality Worklist feature to DICOM standard please refer to the appropriate section in this document.

Note that Modality Performed Procedure Step SOP Class is supported by Discovery NM/CT 670 system in CT Standalone mode via BrightSpeed Elite 16 slices CT Scanner, Discovery NM/CT 670 ES system in CT Standalone mode via Optima CT540 8 Slices CT scanner, by Discovery NM/CT 670 Pro, Discovery 670 DR, Discovery NM/CT 670 CZT and NM/CT 870\* systems in CT Standalone mode via Optima CT540 16 slices CT scanner, by NM/CT 850 system in standalone mode by CT 850 Scanner and by NM/CT 860 system in standalone mode by CT 860 Scanner. Refer to appropriate CT Scanner DICOM Conformance Statements (Section 1.6).

### 7.2 RELATIONSHIP BETWEEN SCHEDULED AND PERFORMED PROCEDURE STEPS

The NM Camera AE supports the following cases:

- Simple Case - one-to-one relationship between Scheduled Procedure Step and PPS
- Append Case - one -to-multiple relationship between Scheduled Procedure Step and PPS
- Unscheduled Case (Acquisition without MWL Data) - zero-to-one relationship between Scheduled Procedure Step and PPS
- Abandoned Case – one-to-zero relationship between Scheduled Procedure Step and PPS

Group case ( multiple-to-one relationship) is not supported.

### 7.3 MODALITY PERFORMED PROCEDURE STEP MODULE TABLE

See DICOM PS 3.3 and PS 3.4 for a complete definition of the entities, modules, and attributes.

**TABLE 7-1  
MODALITY PERFORMED PROCEDURE STEP MODULES**

Module Name	Reference
SOP Common	7.4.1
Performed Procedure Step Relationship	7.4.2
Performed Procedure Step Information	7.4.3
Image Acquisition Results	7.4.4
Radiation Dose	Not Used
Billing and Material Management Codes	Not Used

7.4 MODALITY PERFORMED PROCEDURE STEP MODULE DEFINITIONS

Please refer to DICOM Standard PS 3.3. (Information Object Definitions) for a description of each of the attributes contained within the Modality Performed Procedure Step Information Object Definition.

7.4.1 SOP Common Module

TABLE 7-2  
SOP COMMON MODULE ATTRIBUTES

Attribute Name	Tag	Type for SCU N-CREATE	Type for SCU N-SET	Usage in MPPS Instance
Specific Character Set	(0008,0005)	1C	1C	Character Set that expands or replaces the Basic Graphic Set. Defined Terms used for MPPS of acquisition without MWL Data: ISO_IR 100 = Latin Alphabet No. 1 For MPPS created from MWL request, the value is copied from value provided in MWL if not empty, otherwise ISO_IR 100 is used. Always sent.

7.4.2 Performed Procedure Step Relationship Module

TABLE 7-3  
PERFORMED PROCEDURE STEP RELATIONSHIP MODULE ATTRIBUTES

Attribute Name	Tag	Type for SCU - N-CREATE	
		Acquisition without MWL Entry	Acquisition with MWL Entry
Scheduled Step Attributes Sequence	(0040,0270)	1, contains only 1 item	1, contains only 1 item
>Study Instance UID	(0020,000D)	1, Locally generated UID	1, copied from MWL
>Referenced Study Sequence	(0008,1110)	2, Always sent as empty SQ	2, copied from MWL
>>Referenced SOP Class UID	(0008,1150)	1, Not sent	1, copied from MWL
>>Referenced SOP Instance UID	(0008,1155)	1, Not sent	1, copied from MWL.
>Accession Number	(0008,0050)	2, copied from value, manually entered in UI or sent as empty if no input is provided	2, copied from MWL
>Placer Order Number/Imaging Service Request	(0040,2016)	3, Not sent	3, Not sent
>Filler Order Number/Imaging Service Request	(0040,2017)	3, Not sent	3, Not sent
>Requested Procedure ID	(0040,1001)	2, Sent as empty string	2, copied from MWL

>Requested Procedure Code Sequence	(0032,1064)	3, Not sent	3, filled from first valid item from MWL request or sent empty if no valid items exist in MWL.
>>Code Value	(0008,0100)	1, Not sent	1, filled from first valid item from MWL request, not sent if no valid items exist in MWL
>>Coding Scheme Designator	(0008,0102)	1, Not sent	1, filled from first valid item from MWL request, not sent if no valid items exist in MWL
>>Code Meaning	(0008,0104)	1, Not sent	1, filled from first valid item from MWL request, not sent if no valid items exist in MWL. If value is missing in MWL request, filled by Code Value. (0008,0100)
>Requested Procedure Description	(0032,1060)	2, Sent as empty string	2, copied from MWL request
>Scheduled Procedure Step ID	(0040,0009)	2, Sent as empty string	2, copied from MWL request
>Scheduled Procedure Step Description	(0040,0007)	2, Sent as empty string	2, copied from MWL request
>Scheduled Protocol Code Sequence	(0040,0008)	2, Sent as empty SQ	2, filled from first valid item from MWL request, contains only 1 item
>>Code Value	(0008,0100)	N/A	1, filled from first valid item from MWL request
>>Coding Scheme Designator	(0008,0102)	N/A	1, filled from first valid item from MWL request
>>Coding Scheme Version	(0008,0103)	N/A	3, filled from first valid item from MWL request,
>>Code Meaning	(0008,0104)	N/A	3, filled from first valid item from MWL request. If value is missing in MWL request, filled by Code Value (0008,0100).
Patient's Name	(0010,0010)	2, copied from value, manually entered in UI. Always sent as non-empty value	2, copied from MWL request. Always sent as non-empty value, because MWL requests with missing Patient's Name are rejected by Product.
Patient ID	(0010,0020)	2, copied from value, manually entered in UI. Always sent as non-empty value	2, copied from MWL request. Always sent as non-empty value, because MWL requests with missing Patient ID are rejected by Product.
Issuer of Patient ID	(0010,0021)	3, Not sent	3, Not sent
Patient's Birth Date	(0010,0030)	2, copied from value, manually entered in UI or sent as empty if no input is provided	2, copied from MWL request

DIR DOC0660253 REV 15

Patient's Sex	(0010,0040)	2, copied from value, manually entered in UI or sent as empty if no input is provided.	2, copied from MWL request
Referenced Patient Sequence	(0008,1120)	2, Always sent as empty SQ	2, copied from MWL request or sent as empty SQ if no valid items present in this SQ in MWL request.
>Referenced SOP Class UID	(0008,1150)	N/A	1, copied from MWL request or not sent if no valid items present in this SQ in MWL request.
>Referenced SOP Instance UID	(0008,1155)	N/A	1, copied from MWL request or not sent if no valid items present in this SQ in MWL request.
Admission ID	(0038,0010)	3, Not sent	3, Not sent
Issuer of Admission ID	(0038,0011)	3, Not sent	3, Not sent
Service Episode ID	(0038,0060)	3, Not sent	3, Not sent
Issuer of Service Episode ID	(0038,0061)	3, Not sent	3, Not sent
Service Episode Description	(0038,0062)	3, Not sent	3, Not sent

7.4.3 Performed Procedure Step Information Module

TABLE 7-4  
PERFORMED PROCEDURE STEP INFORMATION MODULE ATTRIBUTES

Attribute Name	Tag	Type for SCU N-CREATE	Type for SCU N-SET	Usage in MPPS Instance
Performed Procedure Step ID	(0040,0253)	1	-	Equipment generated identifier of the protocol carried out within this step. The PPS ID is unique within a study. For MWL scheduled protocols set with "W_" + <SPS ID>. For locally scheduled protocols set with "L_" + numbered id starting from 1 (L_1, L_2).
Performed Station AE Title	(0040,0241)	1	-	Local AE Title, taken from DICOM configuration.
Performed Station Name	(0040,0242)	2	-	User defined name identifying the machine that produced the composite instances "Camera name" value taken from system configuration is used Sent as empty string if it is not received from user input.
Performed Location	(0040,0243)	2	-	Address of the institution where the equipment that produced the composite instances is located. "Institution Address" value taken from system configuration is used Sent as empty string if it is not received from user input.



DIR DOC0660253 REV 15

Performed Procedure Step Start Date	(0040,0244)	1	-	The date that the protocol (SPS) acquisition actually started (doesn't matter if the protocol originated from MWL or was locally scheduled)
Performed Procedure Step Start Time	(0040,0245)	1	-	The time that the protocol (SPS) acquisition actually started (doesn't matter if the protocol originated from MWL or was locally scheduled).
Performed Procedure Step Status	(0040,0252)	1	3	"IN PROGRESS" when acquisition starts "COMPLETED" when operator presses "Protocol Completed" or when all the scans in the protocol are completed. "DISCONTINUED" when operator marks protocol as Discontinued.
Performed Procedure Step Description	(0040,0254)	2	3	The full path of the performed protocol name. E.g. Factory&Cardiology&One Day Always sent as non-empty value
Performed Procedure Type Description	(0040,0255)	2	3	Always sent as empty string
Procedure Code Sequence	(0008,1032)	2	3	N-CREATE: Copied without change from Modality Worklist Requested Procedure Code Sequence (0032,1064) for Acquisition with MWL Entry in Simple or Appended Cases if the protocol originally matched from MWL was not changed. Otherwise sent as empty SQ. N-SET: Copied from Images if presents, otherwise sent as empty SQ.
>Include 'Code Sequence Macro'				
Performed Procedure Step End Date	(0040,0250)	2	3	N-CREATE: sent as empty N-SET: Date when user chooses Complete or Discontinue action.
Performed Procedure Step End Time	(0040,0251)	2	3	N-CREATE: sent as empty N-SET: Time when user chooses Complete or Discontinue action.
Comments on the Performed Procedure Step	(0040,0280)	3	3	User-defined comments on the Performed Procedure Step. Sent as empty string if no user input provided.
Performed Procedure Step Discontinuation Reason Code Sequence	(0040,0281)	3	3	N-CREATE: sent as empty SQ N-SET: if final state is DISCONTINUED, contains only 1 item with reason selected by patient from UI (see Table 7-5 for list of supported reasons). Otherwise, sent as empty SQ.

**TABLE 7-5  
PERFORMED PROCEDURE DISCONTINUE REASONS**

<b>Code Value (0008,0100)</b>	<b>Coding Scheme Designator (0008,0102)</b>	<b>Code Meaning (0008,0104) (Reason description )</b>
110500	DCM	Doctor cancelled procedure
110511	DCM	Nursing unit cancel
110502	DCM	Incorrect procedure ordered
110509	DCM	Change of procedure for correct charging
110510	DCM	Duplicate order
110514	DCM	Incorrect worklist entry selected
110503	DCM	Patient allergic to media/contrast
110505	DCM	Patient refused to continue procedure
110515	DCM	Patient condition prevented continuing
110506	DCM	Patient taken for treatment or surgery
110507	DCM	Patient did not arrive
110508	DCM	Patient pregnant
110504	DCM	Patient died
110501	DCM	Equipment failure
110516	DCM	Equipment change
110513	DCM	Discontinued for unspecified reason

**7.4.4 Image Acquisition Results Module**

**TABLE 7-6  
IMAGE ACQUISITION RESULTS MODULE ATTRIBUTES**

<b>Attribute Name</b>	<b>Tag</b>	<b>Type for SCU N-CREATE</b>	<b>Type for SCU N-SET</b>	<b>Usage in MPPS Instance</b>
Modality	(0008,0060)	1	-	For acquisitions with MWL: copied from MWL request For unscheduled case: set to default value = NM
Study ID	(0020,0010)	2	-	User or equipment generated Study identifier. Copied from Study Name entry of Study Info UI. Always sent as non-empty value.
Performed Protocol Code Sequence	(0040,0260)	2	3	Sent as empty SQ (Assisted Protocol Settings option is not supported)

DIR DOC0660253 REV 15

Performed Series Sequence	(0040,0340)	2	3	N-CREATE: sent as empty SQ N-SET: always sent, may contain 1 or more items, one item for each series reported.
>Performing Physician's Name	(0008,1050)	2	2	Name of the physician(s) administering this Series. Copied from Protocol Info UI. Sent as empty string if no user input provided.
>Protocol Name	(0018,1030)	1	1	User-defined description of the conditions under which the Series was performed. The full path of the performed protocol name. E.g. Factory&Service&Static
>Operator's Name	(0008,1070)	2	2	Name(s) of the operator(s) supporting the Series Copied from user input. Sent as empty string if no user input provided.
>Series Instance UID	(0020,000E)	1	1	Copied from image created
>Series Description	(0008,103E)	2	2	Copied from Image created
>Retrieve AE Title	(0008,0054)	2	2	Always sent as empty string
> Archive Requested	(0040,A494)	3	3	Not Sent
>Referenced Image Sequence	(0008,1140)	2	2	May contain 1 or more items: one item for each image created within the series
>>Referenced SOP Class UID	(0008,1150)	1	1	May contain the following SOP Class UIDs: Nuclear Medicine Image Storage - 1.2.840.10008.5.1.4.1.1.20 For Discovery NM/CT 670*, NM/CT 800* and Optima NM/CT 640 implementations may also contain: CT Image Storage - 1.2.840.10008.5.1.4.1.1.2 Secondary Capture Image Storage – 1.2.840.10008.5.1.4.1.1.7
>>Referenced SOP Instance UID	(0008,1155)	1	1	SOP Instance UID of image object created by PPS
>Referenced Non-Image Composite SOP Instance Sequence	(0040,0220)	2	2	For Discovery NM 630 and Brivo NM 615 implementations : always sent as empty SQ. For Discovery NM/CT 670*,Optima NM/CT 640 and NM/CT 800* implementations: <ul style="list-style-type: none"> <li>For Hybrid acquisitions may contain only 1 item: one item for each Dose SR Dose SR object.</li> <li>For non-hybrid acquisitions send as empty SQ.</li> </ul>

DIR DOC0660253 REV 15

>>Referenced SOP Class UID	(0008,1150)	1	1	Possible values: 1.2.840.10008.5.1.4.1.1.88.67 – X-Ray Radiation Dose SR SOP Class (for Optima NM/CT 640 implementation ) 1.2.840.10008.5.1.4.1.1.88.22 – Enhanced SR SOP Class (For Discovery 670*,NM/CT 800* and Optima NM/CT 640 implementations)
>>Referenced SOP Instance UID	(0008,1155)	1	1	SOP Instance UID of Dose SR object created by PPS.

**7.5 STANDARD EXTENDED AND PRIVATE DATA ATTRIBUTES**

The Product does not support any of the Standard Extended and Private Attributes.

**7.6 STANDARD EXTENDED AND PRIVATE CONTEXT GROUPS**

The Product does not support any Standard Extended, Private, and Configurable Context Groups.

## 8. CT INFORMATION OBJECT IMPLEMENTATION

### 8.1 INTRODUCTION

This section specifies the use of the DICOM CT Image IOD to represent the information included in CT Images produced by Optima NM/CT 640 implementation. Corresponding attributes are conveyed using the module construct.

**Note:** In this section term NM Camera AE DICOM implementation refers to Optima NM/CT 640 DICOM implementation.

For description of CT Images created by Discovery NM/CT 670\*and NM/CT 870\*implementations, refer to CT Scanner DICOM Conformance Statement ([Reference A](#) in Section 1.6).

For description of CT Images created by NM/CT 850 and NM/CT 860 implementations, refer to CT Scanner DICOM Conformance Statement ([Reference B](#) in Section 1.6).

### 8.2 NM CAMERA MAPPING OF DICOM ENTITIES

NM Camera AE maps DICOM Information Entities to local Information Entities in the product’s database and user interface.

**TABLE 8-1**  
**MAPPING OF DICOM ENTITIES TO NM CAMERA ENTITIES**

DICOM IE	NM Camera Entity
Patient	Patient
Study	Study
Series	Series
Image	Image

### 8.3 IOD MODULE TABLE

The Computed Tomography Information Object Definition comprises the modules of the following table, plus Standard Extended and Private attributes. Standard Extended and Private attributes are described in Section 8.5.

**TABLE 8-2**  
**CT IMAGE IOD MODULES**

Entity Name	Module Name	Usage	Reference
Patient	Patient	Used (same description as for NM IOD)	3.4.1.1
	Clinical Trial Subject	Not Used	NA
Study	General Study	Used (same description as for NM IOD)	3.4.2.1

	Patient Study	Used (same description as for NM IOD)	3.4.2.2
	Standard Extended Study	Used (same description as for NM IOD)	3.4.2.3
	Private Study	Used (same description as for NM IOD)	3.4.2.4
	Clinical Trial Study	Not Used	NA
Series	General Series	Used (same description as for NM IOD)	3.4.3.1
	Private Series	Used (same description as for NM IOD)	3.4.3.4
	Private CT Series	Used	8.4.1.1
	Clinical Trial Series	Not Used	NA
Frame of Reference	Frame of Reference	Used (same description as for NM IOD)	3.4.4.1
Equipment	General Equipment	Used (same description as for NM IOD)	3.4.5.1
Image	General Image	Used	8.4.3.1
	Image Plane	Used	8.4.3.2
	Image Pixel	Used	8.4.3.3
	Contrast/Bolus	Used for acquisitions which use contrast.	8.4.3.4
	Device	Not Used	NA
	Standard Extended Image	Used	8.4.3.5
	Private Image	Used(same description as for NM IOD)	3.4.6.15
	CT Image	Used	8.4.3.6
	Private CT Image	Used	8.4.3.9
	Overlay Plane	Not Used	NA
	VOI LUT	Used	8.4.3.7
	SOP Common	Used	8.4.3.8

**8.4 INFORMATION MODULE DEFINITIONS**

Please refer to DICOM Part 3 (Information Object Definitions) for a description of each of the entities, modules, and attributes contained within the CT Information Object.

The following modules are included to convey Enumerated Values, Defined Terms, and Optional Attributes supported. Type 1 & Type 2 Attributes are also included for completeness and to define what values they may take and where these values are obtained from when generating the instance. It should be noted that they are the same ones as defined in the DICOM Standard Part 3 (Information Object Definitions). Also note that Attributes not present in tables are not supported.

DIR DOC0660253 REV 15

8.4.1 Series Entity Modules

8.4.1.1 Private CT Series Module

TABLE 8-3  
PRIVATE CT SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Imagesets in Acquisition	(0037,xx91)	QUASAR_INTERNAL_USE	Number of Imagesets created by specific scan and stored within the Series.

8.4.2 Frame Of Reference Entity Modules

8.4.2.1 Frame Of Reference Module

TABLE 8-4  
FRAME OF REFERENCE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Frame of Reference UID	(0020,0052)	1	See 3.4.4.1 for further explanation.
Position Reference Indicator	(0020,1040)	2	See 3.4.4.1 for further explanation.

8.4.3 Image Entity Modules

8.4.3.1 General Image Module

TABLE 8-5  
GENERAL IMAGE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Instance Number	(0020,0013)	2	A number that identifies this image.
Patient Orientation	(0020,0020)	2C	Not Sent
Content Date	(0008,0023)	2C	The date the image pixel data creation started.
Content Time	(0008,0033)	2C	The time the image pixel data creation started
Image Type	(0008,0008)	3	See 8.4.3.6.1
Acquisition Date	(0008,0022)	3	The date the acquisition of data that resulted in this image started
Acquisition Time	(0008,0032)	3	The time the acquisition of data that resulted in this image started
Images in Acquisition	(0020,1002)	3	Number of images expected to be created by scan. Always sent.
Quality Control Image	(0028,0300)	3	Indicates whether or not this image is a quality control or phantom image. Enumerated Values: YES NO Always sent.

8.4.3.2 Image Plane Module

TABLE 8-6  
IMAGE PLANE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Pixel Spacing	(0028,0030)	1	Physical distance in the patient between the center of each pixel, specified by a numeric pair - adjacent row spacing (delimiter) adjacent column spacing in mm. (*)
Image Orientation (Patient)	(0020,0037)	1	The direction cosines of the first row and the first column with respect to the patient.
Image Position (Patient)	(0020,0032)	1	The x, y, and z coordinates of the upper left hand corner (center of the first voxel transmitted) of the image, in mm
Slice Thickness	(0018,0050)	2	Nominal slice thickness, in mm (*). Possible values: 2.5, 5 or 10 mm Always sent.
Slice Location	(0020,1041)	3	Relative position of the image plane expressed in mm. Explains the relationship of this value to the table and gantry (i.e. define the implementation specific reference point), and to the values in Image Orientation (Patient) (0020,0037) and Image Position (Patient) (0020,0032). Z is relative to reference point (=Landmark), which for the time being is the zero coordinate on pallet. Currently is equal to Z value of Image Position. Always sent.

**Note 1 :** (\*) - Attributes copied from the user input

8.4.3.3 Image Pixel Module

TABLE 8-7  
IMAGE PIXEL MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Samples per Pixel	(0028,0002)	1	See 8.4.3.5 for CT Images
Photometric Interpretation	(0028,0004)	1	See 8.4.3.5 for CT Images
Rows	(0028,0010)	1	Number of rows in the image (*).
Columns	(0028,0011)	1	Number of columns in the image(*)
Bits Allocated	(0028,0100)	1	See 8.4.3.5 for CT Images
Bits Stored	(0028,0101)	1	See 8.4.3.5 for CT Images
High Bit	(0028,0102)	1	See 8.4.3.5 for CT Images



Pixel Representation	(0028,0103)	1	Data representation of the pixel samples. Each sample shall have the same pixel representation. Enumerated Values used: 0000H = unsigned integer. 0001H = 2's complement
Pixel Data	(7FE0, 0010)	1	A data stream of the pixel samples that comprise the Image.
Planar Configuration	(0028,0006)	1C	Not Used (number of Samples per Pixel is always 1)
Pixel Aspect Ratio	(0028,0034)	1C	Not Used

**Note 1 :** (\*) - Attributes copied from the user input

8.4.3.4 Contrast/Bolus Module

**TABLE 8-8  
CONTRAST/BOLUS MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Contrast/Bolus Agent	(0018,0010)	2	Contrast or bolus agent (*) Sent as empty string if no user input is provided.
Contrast/Bolus Route	(0018,1040)	3	Administration route of contrast agent (*) Sent as empty string if no user input is provided. Possible values: <ul style="list-style-type: none"> <li>• ORAL</li> <li>• IV</li> <li>• ORAL and IV</li> </ul>

**Note 1 :** (\*) - Attributes copied from the user input

8.4.3.5 Standard Extended Image Module

**TABLE 8-9  
STANDARD EXTENDED IMAGE MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Counts Accumulated	(0018,0070)	3	Total counts in the pixel data
Spacing Between Slices	(0018,0088)	3	Spacing between slices, in mm. The spacing is measured from the center-to-center of each slice. (*) Sent, if Image Type (0008,0008) Value 3 is AXIAL.
Table Traverse	(0018,1131)	3	Table Traverse (*)

**Note 1 :** (\*) - Attributes copied from the user input

8.4.3.6 CT Image Module

TABLE 8-10  
CT IMAGE MODULE ATTRIBUTES

Attribute Name	Tag	Type	Use
Image Type	(0008,0008)	1	See 8.4.3.6.1
Samples per Pixel	(0028,0002)	1	Shall be 1.
Photometric Interpretation	(0028,0004)	1	Always set to MONOCHROME2
Bits Allocated	(0028,0100)	1	Shall be 16.
Bits Stored	(0028,0101)	1	Shall be 16.
High Bit	(0028,0102)	1	Shall have only the Enumerated Value of one less than the value sent in Bits Stored.
Rescale Intercept	(0028,1052)	1	The value b in relationship between stored values (SV) and the output units. Output units = m*SV+b If Image Type (0008,0008) Value 1 is ORIGINAL and Value 3 is not LOCALIZER, output units shall be Hounsfield Units (HU). If Image Type (0008,0008) Value 3 is AXIAL , sent as -1000 If Image Type (0008,0008) Value 3 is LOCALIZER, sent as 0
Rescale Slope	(0028,1053)	1	m in the equation specified in Rescale Intercept (0028,1052). Always sent as 1
Rescale Type	(0028,1054)	1C	Specifies the output units of Rescale Slope (0028,1053) and Rescale Intercept (0028,1052). If Image Type (0008,0008) Value 3 is AXIAL , sent as HU If Image Type (0008,0008) Value 3 is LOCALIZER, sent as US
KVP	(0018,0060)	2	Peak kilo voltage output of the x-ray generator used. (*) Possible Values: 120,140
Acquisition Number	(0020,0012)	2	A number identifying the single continuous gathering of data over a period of time which resulted in this image.

Scan Options	(0018,0022)	3	Parameters of scanning sequence. Defined Terms: HELICAL SCAN AXIAL SCAN STATIC SCAN SCOUT SCAN Always sent as non-empty value.
Data Collection Diameter	(0018,0090)	3	The diameter in mm of the region over which data were collected. See 8.4.3.6.2 Sent If Image Type (0008,0008) Value 3 is AXIAL. Default value = 500
Data Collection Center (Patient)	(0018,9313)	3	The x, y, and z coordinates (in the patient coordinate system) in mm of the center of the region in which data were collected. See 8.4.3.6.2 Sent, if Image Type (0008,0008) Value 3 is AXIAL. The following values are used: x is 0, y is 0 z defines Table Coordinate of Center.
Reconstruction Diameter	(0018,1100)	3	Diameter in mm of the region from within which data were used in creating the reconstruction of the image. See 8.4.3.6.2 Sent, if Image Type (0008,0008) Value 3 is AXIAL. (*) Default Value is 500.
Reconstruction Target Center (Patient)	(0018,9318)	3	The x, y, and z coordinates (in the patient coordinate system) of the reconstruction center target point as used for reconstruction in mm. See 8.4.3.6.2 Sent, if Image Type (0008,0008) Value 3 is AXIAL. The following values are used: x defines A/P OFF CENTER, default value is 0 y defines R/L OFF CENTER, default value is 0 z defines Table Coordinate of Center, equal to the Z value of (0018,9313) - Data Collection Center (Patient) .

Distance Source to Detector	(0018,1110)	3	Distance in mm from source to detector center. See 8.4.3.6.2 Always sent Possible Value: 949 mm
Distance Source to Patient	(0018,1111)	3	Distance in mm from source to iso-center (center of field of view). See 8.4.3.6.2 Always sent Possible Value: 541 mm
Gantry/Detector Tilt	(0018,1120)	3	Nominal angle of tilt in degrees of the scanning gantry. Sent, if Image Type (0008,0008) Value 3 is AXIAL. The value used: 0.
Table Height	(0018,1130)	3	The distance in mm of the top of the patient table to the center of rotation;(*). Always sent.
Rotation Direction	(0018,1140)	3	Direction of rotation of the source when relevant, about nearest principal axis of equipment. Enumerated Values: CW = clockwise Sent, if Image Type (0008,0008) Value 3 is AXIAL.
Exposure Time	(0018,1150)	3	Time of x-ray exposure in msec. Always sent .
X-ray Tube Current	(0018,1151)	3	X-Ray Tube Current in mA.(*) Always sent. May vary per slice, if Image Type (0008,0008) Value 3 is AXIAL and Exposure Modulation Type (0018,9323) is AUTO_MA
Exposure	(0018,1152)	3	The exposure expressed in mAs, for example calculated from Exposure Time and X-Ray Tube Current. Always sent. Calculated as $X\text{-Ray Tube Current} * (\text{Exposure Time} / 1000)$ May vary per slice, if Image Type (0008,0008) Value 3 is AXIAL and Exposure Modulation Type is AUTO_MA

Filter Type	(0018,1160)	3	Label for the type of filter inserted into the x-ray beam. Always sent. Possible value : BOWTIE
Generator Power	(0018,1170)	3	Power in kW to the x-ray generator. Always sent. Calculated as : (KVP* Tube Current/1000)
Focal Spot	(0018,1190)	3	Size of the focal spot in mm. Always sent as pair : Focal Spot width \Focal Spot length Possible values: Focal Spot width: 0.6 mm Focal Spot length: 0.8 mm
Convolution Kernel	(0018,1210)	3	A label describing the convolution kernel or algorithm used to reconstruct the data. (*) Always sent, if Image Type (0008,0008) Value 3 is AXIAL. Possible values: <ul style="list-style-type: none"> <li>• Standard</li> <li>• Bone</li> <li>• Smooth</li> <li>• Detailed</li> </ul>
Revolution Time	(0018,9305)	3	The time in seconds of a complete revolution of the source around the gantry orbit. Sent, if Image Type (0008,0008) Value 3 is AXIAL.(*) Possible values: 1 or 2 sec
Single Collimation Width	(0018,9306)	3	The width of a single row of acquired data (in mm). Always sent. Possible value: 2.5 mm
Total Collimation Width	(0018,9307)	3	The width of the total collimation (in mm) over the area of active x-ray detection. Note: This will be equal the number of effective detector rows multiplied by single collimation width. Always sent. Possible value: 10 mm

Table Speed	(0018,9309)	3	<p>The distance in mm that the table moves in one second during the gathering of data that resulted in this image.</p> <p>Sent if Image Type (0008,0008) Value 3 is AXIAL and Scan Options (0018,0022) is HELICAL SCAN or Image Type (0008,0008) Value 3 is LOCALIZER</p> <p>If Image Type (0008,0008) Value 3 is LOCALIZER, possible value is 100 mm/sec</p>
Table Feed per Rotation	(0018,9310)	3	<p>Motion of the table (in mm) during a complete revolution of the source around the gantry orbit.</p> <p>Sent, if Image Type (0008,0008) Value 3 is AXIAL and Scan Options (0018,0022) is HELICAL SCAN.</p> <p>Calculated as :</p> <p>Spiral Pitch Factor*Total Collimator Width</p>
Spiral Pitch Factor	(0018,9311)	3	<p>Ratio of the Table Feed per Rotation (0018,9310) to the Total Collimation Width (0018,9307).</p> <p>Sent, if Image Type (0008,0008) Value 3 is AXIAL and Scan Options (0018,0022) is HELICAL SCAN. (*)</p> <p>Possible values: 0.75 1.25, 1.75</p>
Exposure Modulation Type	(0018,9323)	3	<p>A label describing the type of exposure modulation used for the purpose of limiting the dose.</p> <p>Defined Terms:</p> <p>NONE</p> <p>AUTO_MA</p> <p>Sent, if Image Type (0008,0008) Value 3 is AXIAL.(*)</p>
Estimated Dose Saving	(0018,9324)	3	<p>A percent value of dose saving due to the use of Exposure Modulation Type (0018,9323). A negative percent value of dose savings reflects an increase of exposure.</p> <p>Sent, if Image Type (0008,0008) Value 3 is AXIAL.</p> <p>Sent as 0 if Exposure Modulation Type is NONE</p>

CTDIvol	(0018,9345)	3	Computed Tomography Dose Index (CTDIvol), in mGy. It describes the average dose for this image for the selected CT conditions of operation. Sent, if Image Type (0008,0008) Value 3 is AXIAL. May vary per slice, if Exposure Modulation Type is AUTO_MA
CTDI Phantom Type Code Sequence	(0018,9346)	3	The type of phantom used for CTDI Measurement. Only a single Item shall be permitted in this Sequence. Sent, if Image Type (0008,0008) Value 3 is AXIAL.
>Include Code Sequence Macro	DCID (4052) Phantom Devices is supported The following values are used: (DCM, 113690, IEC Head Dosimetry Phantom) (DCM, 113691, IEC Body Dosimetry Phantom)		

**Note 1 :** (\*) - Attributes copied from the user input

**8.4.3.6.1 Image Type**

The following Enumerated Values of Value 1 are created:

- ORIGINAL identifies an Original Image
- DERIVED identifies a Derived Image

The following Enumerated Values of Value 2 are created:

- PRIMARY identifies a Primary Image

The following Defined Terms of Value 3 are created:

- AXIAL identifies a CT Axial Image
- LOCALIZER identifies a CT Localizer Image

**8.4.3.6.2 Relationships Between CT Geometric Attributes**

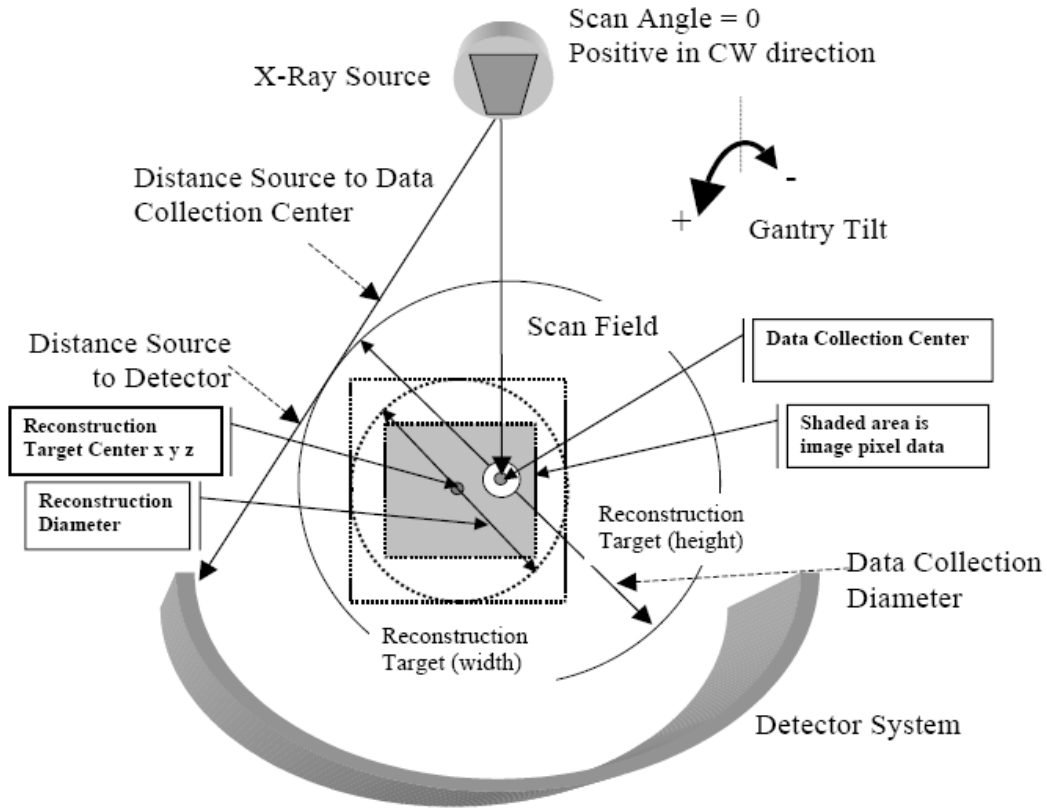


Figure C.8-19: Geometry of CT Acquisition System

8.4.3.7 VOI LUT module

TABLE 8-11  
VOI LUT MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
VOI LUT Sequence	(0028,3010)	1C	Not Used
Window Center	(0028,1050)	1C	Window Center for display. Only single value is present. (*) Sent as 40, if Image Type (0008,0008) Value 3 is AXIAL Sent as 440, if Image Type (0008,0008) Value 3 is LOCALIZER Always sent Derived from protocol defined preset values.



Window Width	(0028,1051)	1C	Window Width for display. Only single value is present. (*) Sent as 400, if Image Type (0008,0008) Value 3 is AXIAL Sent as 1100, if Image Type (0008,0008) Value 3 is LOCALIZER Always sent Derived from protocol defined preset values.
--------------	-------------	----	---

**Note 1 :** (\*) - Attributes copied from the user input

8.4.3.8 SOP Common Module

**TABLE 8-12**  
**SOP COMMON MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
SOP Class UID	(0008,0016)	1	Uniquely identifies the SOP Class. Always set to 1.2.840.10008.5.1.4.1.1.2
SOP Instance UID	(0008,0018)	1	Uniquely identifies the SOP Instance. Internally generated.
Specific Character Set	(0008,0005)	1C	Character Set that expands or replaces the Basic Graphic Set. Defined Terms include for locally created images: ISO_IR 100 = Latin Alphabet No. 1 For images created from MWL, the value is copied from value provided in MWL, if not empty, otherwise ISO_IR 100 is used. Always included into image
Instance Number	(0020,0013)	3	See 3.4.6.1 for more specialization

8.4.3.9 Private CT Image Module

**TABLE 8-13**  
**PRIVATE CT IMAGE MODULE ATTRIBUTES**

Attribute Name	Tag	Private Creator ID	Attribute Description
ImagesetSOPInstanceUID	(0009,xx46)	GEMS_GENIE_1	UID of CT Imageset. Shared for CT slices dedicated to be joined in one dataset on Xeleris database
Dataset UID	(0009,xx1E)	GEMS_GENIE_1	UID of Dataset object in Xeleris WS database. Always equal to value of (0009,xx46) ImagesetSOPInstanceUID

ACGroup	(0037,xx80)	QUASAR_INTERNAL_USE	Defines AC CT ImageSet(s) Contains Number of AC Group (1,2,etc) for AC ImageSet, for all other Imagesets contains -1 Sent, if Image Type (0008,0008) Value 3 is AXIAL.
ScanDirection	(0037,xx81)	QUASAR_INTERNAL_USE	Scan direction. Possible values: TO_HEAD TO_FEET
ScanRange	(0037,xx82)	QUASAR_INTERNAL_USE	Range exposed to X-ray, in mm Always sent
BreathingTime	(0037,xx83)	QUASAR_INTERNAL_USE	Breathing Time in seconds (*) Value -1 means that Breathing Time is not defined for the scan. Always sent
BreathholdTime	(0037,xx84)	QUASAR_INTERNAL_USE	Breathhold Time in seconds (*). Value -1 means that Breathhold Time is not defined for the scan Always Sent
BreathHoldFlag	(0037,xx85)	QUASAR_INTERNAL_USE	Defines, if Imageset is acquired during breathhold time (*). Possible Values: TRUE, FALSE Always sent
Acquisition Parent UID (CT)	(0031, xx02)	GEHC_HYBRID_01	Shared by all images created by same Hybrid scan. Equal to Acquisition Parent UID (0011,xx31) GEMS_GENIE_1 tag value if it exists. Also used for correct registration. For more details about NM/CT registration refer to Frame Of Reference, Section 3.4.4.1

Note 1 : (\*) - Attributes copied from the user input

**8.5 STANDARD EXTENDED AND PRIVATE DATA ATTRIBUTES**

The Product supports the Standard and Private Attributes defined in the following sections in Standard Extended CT SOP Instances as Type 3 data elements.

**8.5.1 Standard Attributes**

The Product supports the following standard attributes, not specified in the CT IOD, in SOP Instances as Type 3 data elements.

**TABLE 8-14  
STANDARD EXTENDED ATTRIBUTES**

Information Entity Name	Attribute Name	Tag	Use
Study	Allergies	(0010,2110)	Description of prior reaction to contrast agents, or other patient allergies or adverse reactions.
	Pregnancy Status	(0010,21C0)	Describes pregnancy state of patient Enumerated Values: 0001 = not pregnant 0002 = possibly pregnant 0003 = definitely pregnant 0004 = unknown
	Requested Procedure Comments	(0040,1400)	User-defined Study notes
Image	Counts Accumulated	(0018,0070)	Total counts in pixel data
	Spacing Between Slices	(0018,0088)	Spacing between slices, in mm. The spacing is measured from the center-to-center of each slice.
	Table Traverse	(0018,1131)	Table Traverse

**8.5.2 Private Group GEMS\_GENIE\_1**

**TABLE 8-15  
PRIVATE GROUP PRIVATE GROUP GEMS\_GENIE\_1**

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0009,00xx)	LO	1	GEMS_GENIE_1
Imageset SOP Instance UID	(0009,xx46)	UI	1	UID of CT ImageSet. Shared for CT slices dedicated to be joined in one dataset on Xeleris database
Dataset UID	(0009,xx1E)	UI	1	UID of Dataset object in Xeleris WS database.

DIR DOC0660253 REV 15

8.5.3 Private Group QUASAR\_INTERNAL\_USE

TABLE 8-16  
PRIVATE GROUP QUASAR INTERNAL USE

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0037,00xx)	LO	1	QUASAR_INTERNAL_USE
AC Group	(0037,xx80)	IS	1	Defines AC CT ImageSet(s)
Scan Direction	(0037,xx81)	CS	1	Scan direction.
Scan Range	(0037,xx82)	DS	1	Range exposed to X-ray, in mm
Breathing Time	(0037,xx83)	IS	1	Breathing Time in seconds
Breathhold Time	(0037,xx84)	IS	1	Breathhold Time in seconds.
BreathHoldFlag	(0037,xx85)	CS	1	Defines, if Imageset is acquired during breathhold time
Imagesets in Acquisition	(0037,xx91)	IS	1	Number of Imagesets created by specific scan and stored within the Series.

8.5.4 Private Group GEHC\_HYBRID\_01

TABLE 8-17  
PRIVATE GROUP GEHC HYBRID\_01

Attribute Name	Tag	VR	VM	Attribute Description and Use
Acquisition Parent UID (CT)	(0031, xx02)	UI	1	Shared by all images created by same Hybrid scan. Equal to Acquisition Parent UID (0011,xx31) GEMS_GENIE_1 tag value if it exists. Used for correct registration.

## 9. CT RADIATION DOSE STRUCTURED REPORT INFORMATION OBJECT IMPLEMENTATION

### 9.1 INTRODUCTION

This section specifies the use of the DICOM CT Radiation Dose SR IOD to represent results produced by Optima NM/CT 640 DICOM implementation. Corresponding attributes are conveyed using the module construct.

**Note:** In this section term NM Camera AE DICOM implementation refers to Optima NM/CT 640 DICOM implementation.

### 9.2 NM CAMERA MAPPING OF DICOM ENTITIES

The NM Camera AE maps DICOM Information Entities to local Information Entities in the product’s database and user interface.

**TABLE 9-1  
MAPPING OF DICOM ENTITIES TO NM CAMERA ENTITIES**

DICOM IE	NM Camera Entity
Patient	Patient
Study	Study
Series	Series
Document	Report

#### 9.2.1 IOD Module Table

The CT Radiation Dose Structured Report Information Object Definitions comprise the modules of the following tables plus Standard Extended and Private Attributes. Standard Extended and Private Attributes are described in Section 9.4

The contents of the SR Document Content are constrained by the supported template, as identified in Section 9.3.3.2.1.1. Standard, Standard Extended and Private Templates are further described in Section 9.6.

**TABLE 9-2  
STRUCTURE REPORT IOD MODULES**

Entity Name	Module Name	Usage	Reference
Patient	Patient	Used (same description as for NM IOD)	3.4.1.1
	Clinical Trial Subject	Not Used	N/A
Study	General Study	Used (same description as for NM IOD)	3.4.2.1
	Patient Study	Used (same description as for NM IOD)	3.4.2.2
	Standard Extended Study	Used (same description as for NM IOD)	3.4.2.3
	Private Study	Used (same description as for NM IOD)	3.4.2.4

	Clinical Trial Study	Not Used	N/A
Series	SR Document Series	Used	9.3.1.1
	Standard Extended Series	Used	9.3.1.2
	Private Series	Used	9.3.1.3
	Clinical Trial Series	Not Used	N/A
Frame Of Reference	Synchronization	Not Used	N/A
Equipment	General Equipment	Used (same description as for NM IOD)	3.4.5.1
	Enhanced General Equipment	Used	9.3.2.1
Document	SR Document General	Used	9.3.3.1
	SR Document Content	Used	9.3.3.2
	Standard Extended SR Document	Used	9.3.3.3
	Private SR Document	Used	9.3.3.4
	SOP Common	Used	9.3.3.5

**9.3 INFORMATION MODULE DEFINITIONS**

Please refer to DICOM Part 3 (Information Object Definitions) for a description of each of the entities, modules, and attributes contained within the SR Information Objects.

The following modules are included to convey Enumerated Values, Defined Terms, and Optional Attributes supported. Type 1 & Type 2 Attributes are also included for completeness and to define what values they may take and where these values are obtained from when generating the instance. It should be noted that they are the same ones as defined in the DICOM Standard Part 3 (Information Object Definitions). Also note that Attributes which are not presented in tables are not supported.

**9.3.1 Series Entity Modules**

**9.3.1.1 SR Document Series Module**

**TABLE 9-3  
SR DOCUMENT SERIES MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Use
Modality	(0008,0060)	1	Defined Terms used for data created on this system: Value = SR
Series Instance UID	(0020,000E)	1	Internally generated unique identifier of the Series.
Series Number	(0020,0011)	1	A number that identifies this Series. Internally generated Default Value = 1

Series Description	(0008,103E)	3	Description of the Series. Always sent. Default value CT Radiation Dose Report
Series Date	(0008,0021)	3	Date the Series started
Series Time	(0008,0031)	3	Time the Series started.
Referenced Performed Procedure Step Sequence	(0008,1111)	2	Uniquely identifies the Performed Procedure Step SOP Instance to which the Series is related. The sequence has exactly 1 item. Sequence is added to all image(s) created by system.
>Referenced SOP Class UID	(0008,1150)	1	Set to 1.2.840.10008.3.1.2.3.3
>Referenced SOP Instance UID	(0008,1155)	1	Uniquely identifies the referenced SOP Instance. Internally generated.

9.3.1.2 Standard Extended Series Module

TABLE 9-4  
STANDARD EXTENDED SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Performing Physicians' Name	(0008,1050)	3	Name of the physician(s) administering this Series. (*) (**) Sent as empty string if no user input provided.
Protocol Name	(0018,1030)	3	User-defined description of the conditions under which the Series was performed.  The full path of the performed protocol name. E.g. Factory&Service&Static Always sent with non-empty value.
Request Attributes Sequence	(0040,0275)	3	Sequence that contains attributes from the Imaging Service Request. The sequence has exactly 1 item. Used when Series as created as result of MWL request, not sent otherwise.
>Requested Procedure ID	(0040,1001)	3	Identifier that identifies the Requested Procedure in the Imaging Service Request. (*)
>Accession Number	(0008,0050)	3	A RIS generated number that identifies the order for the Study (*)  May be updated by user manually.
>Study Instance UID	(0020,000D)	3	Unique identifier for the Study. (*) (***)
>Referenced Study Sequence	(0008,1110)	3	Uniquely identifies the Study SOP Instances associated with this SOP Instance. The sequence has exactly 1 item. (*) (***)
>> <i>Include 'SOP Instance Reference Macro'</i>			
>Requested Procedure Description	(0032,1060)	3	Institution-generated administrative description or classification of Requested Procedure. (*) (***)
>Scheduled Procedure Step ID	(0040,0009)	3	Identifier that identifies the Scheduled Procedure Step. (*) (***)

>Scheduled Procedure Step Description	(0040,0007)	3	Institution-generated description or classification of the Scheduled Procedure Step to be performed. (*) (***)
>Scheduled Protocol Code Sequence	(0040,0008)	3	Sequence describing the Scheduled Protocol following a specific coding scheme. (*) (***) The sequence has exactly 1 item. If MWL request contains more than 1 item, only the first valid item is copied.
>>Include 'Code Sequence Macro'			
Performed Procedure Step ID	(0040,0253)	3	Equipment generated identifier of the protocol carried out within this step. The PPS ID is unique within a study. For MWL scheduled protocols set with W_ + <SPS ID>. For locally scheduled protocols set with L_ + numbered id starting from 1 (L_1, L_2). Sent as empty string for QC protocols. Always sent.
Performed Procedure Step Start Date	(0040,0244)	3	The date that the protocol (SPS) acquisition actually started (doesn't matter if the protocol originated from MWL or was locally scheduled). A locally scheduled protocol is an SPS that is created/added in the camera. Always sent.
Performed Procedure Step Start Time	(0040,0245)	3	The time that the protocol (SPS) acquisition actually started (doesn't matter if the protocol originated from MWL or was locally scheduled). A locally scheduled protocol is an SPS that is created/added in the camera. Always sent
Performed Procedure Step Description	(0040,0254)	3	The full path of the performed protocol name. E.g. Factory&Cardiology&One Day Always sent

**Note 1 :** (\*) – Attributes copied from the Worklist if the study source was actually copied from a Worklist query result (if available).

**Note 2 :** (\*\*) – Attributes copied from the user input for Locally scheduled Protocols

**Note 3 :** (\*\*\*) – Cannot be modified by user if received from MWL

9.3.1.3 Private Series Module

TABLE 9-5  
PRIVATE SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
----------------	-----	--------------------	-----------------------



Protocol Scheduled Date	(0009,xx40)	QUASAR_INTERNAL_USE	Protocol Scheduled Date (*). Taken from the SPS Start date of the first SPS in the study – Tag (0040,0002).  For QC scans – sent as ZERO Length Always sent
Protocol Scheduled Time	(0009,xx41)	QUASAR_INTERNAL_USE	Protocol Scheduled Time (*) Taken from the SPS Start time of the first SPS in the study – Tag (0040,0003).  For QC scans – sent as ZERO Length. Always sent
Matched protocol	(0009,xx43)	QUASAR_INTERNAL_USE	For Worklist items. The originally matched protocol vs. protocol name which is the protocol actually acquired.  Not sent for locally scheduled protocols.
Private SPS ID	(0009,xx44)	QUASAR_INTERNAL_USE	Defines protocols that were appended to the original MWL protocol. Contains SPS ID of original MWL Protocol.  Sent empty for images of locally scheduled protocols and appended to originally locally scheduled protocols.
Pre-Medication	(0009,xx45)	QUASAR_INTERNAL_USE	Keeps the Pre-Medication as appears in the To-Do list (*) (**)  Sent as empty string if not received from user input.

**Note 1 :** (\*) – Attributes copied from the Worklist if the study source was actually copied from a Worklist query result (if available).

**Note 2 :** (\*\*) – Attributes copied from the user input for Locally scheduled Protocols

9.3.2 Equipment Entity Modules

9.3.2.1 Enhanced General Equipment Module

**TABLE 9-6  
ENHANCED GENERAL EQUIPMENT MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
Manufacturer	(0008,0070)	1	Manufacturer of the equipment that produced the composite instances.  Default Value GE MEDICAL SYSTEMS
Manufacturer's Model Name	(0008,1090)	1	Manufacturer's model name of the equipment that produced the composite instances.  Set to Tandem_Optima_640 for Optima NM/CT 640 implementation
Device Serial Number	(0018,1000)	1	Manufacturer's serial number of the equipment that produced the composite instances.  Always taken from system configuration.

Software Versions	(0018,1020)	1	Manufacturer's designation of software version of the equipment that produced the composite instances  Contains Software/Hardware versions of current release e.g. 1.003.429.0\HARDWARE_VERSION_1
-------------------	-------------	---	--

9.3.3 Document Entity Modules

9.3.3.1 SR Document General Module

TABLE 9-7  
SR DOCUMENT GENERAL MODULE ATTRIBUTES

Attribute Name	Tag	Type	Use
Instance Number	(0020,0013)	1	A number that identifies this document. Default Value = 1
Completion Flag	(0040,A491)	1	Enumerated Value used: COMPLETE = Complete content.
Verification Flag	(0040,A493)	1	Enumerated Values: UNVERIFIED = Not attested to.
Completion Flag Description	(0040,A492)	3	Explanation of the value sent in Completion Flag (0040,A491).
Content Date	(0008,0023)	1	The date of the document creation.
Content Time	(0008,0033)	1	The time of the document creation.
Verifying Observer Sequence	(0040,A073)	1C	Not sent.
Referenced Request Sequence	(0040,A370)	1C	Not sent
Performed Procedure Code Sequence	(0040,A372)	2	A Sequence that conveys the codes of the performed procedures pertaining to this SOP Instance. Always contains 0 Items.
Current Requested Procedure Evidence Sequence	(0040,A375)	1C	Not Sent
Pertinent Other Evidence Sequence	(0040,A385)	1C	Not Sent
Referenced Instance Sequence	(0008,114A)	1C	Not Sent

9.3.3.2 SR Document Content Module

TABLE 9-8  
SR DOCUMENT CONTENT MODULE ATTRIBUTES

Attribute Name	Tag	Type	Use
Value Type	(0040,A040)	1	The type of the value encoded in this Content Item.  Defined Term used: CONTAINER

Concept Name Code Sequence	(0040,A043)	1C	Code describing the concept represented by this Content Item. Only a single Item is included in this sequence.
>Include 'Code Sequence Macro'			EV (113701, DCM,X-Ray Radiation Dose Report)
Continuity of Content	(0040,A050)	1	Enumerated Value used: SEPARATE
Content Template Sequence	(0040,A504)	1	Template that describes the content of this Content Item and its subsidiary Content Items. Only a single Item is included in this sequence
>Mapping Resource	(0008,0105)	1	Mapping Resource that defines the template. Defined Terms: DCMR = DICOM Content Mapping Resource
>Template Identifier	(0040,DB00)	1	10011 - CT Radiation Dose (9.3.3.2.1.1)
Observation DateTime	(0040,A032)	1	The date and time on which this Content Item was completed. Always sent.
Content Sequence	(0040,A730)	1	Provides the hierarchical structuring of the Content Tree by recursively nesting Content Items. A parent (or source) Content Item has an explicit relationship to each child (or target) Content Item, conveyed by the Relationship Type (0040,A010) attribute. One or more Items may be included in this sequence.
> Relationship Type	(0040,A010)	1	Defined Term used: CONTAINS HAS CONCEPT MOD HAS OBS CONTEXT
> Insert SR Document Content Module			Recursive inclusion to create document content tree. See section 9.3.3.2.1 for the list of supported templates.

9.3.3.2.1 SR Document Content Descriptions

9.3.3.2.1.1 Content Template

The product supports the following root Templates for SR SOP Instances created by the product.

TABLE 9-9  
SR ROOT TEMPLATES

SOP Class	Template ID	Template Name	Use
X-Ray Radiation Dose SR	10011	CT Radiation Dose	Create

DIR DOC0660253 REV 15

Refer to section 9.6 for a detailed description of the supported templates.

9.3.3.3 Standard Extended SR Document Module

TABLE 9-10  
STANDARD EXTENDED DOSE SR MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Acquisition Date	(0008,0022)	3	The date of the irradiation event that resulted in this Dose SR started
Acquisition Time	(0008,0032)	3	The time of the irradiation event that resulted in this Dose SR started
Quality Control Image	(0028,0300)	3	Indicates whether or not this Dose report is created by a quality control procedure or not. Enumerated Values: YES NO

9.3.3.4 Private SR Document Module

TABLE 9-11  
PRIVATE SR DOCUMENT MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Origin	(0009,xx12)	QUASAR_INTERNAL_USE	The origin of the image.”isworklist” if scheduled in MWL,” regular” if locally scheduled.
Image Type	(0009,xx1B)	QUASAR_INTERNAL_USE	Image type string as passed in the scan request Always sent as “ct dose report”
Dataset Name	(0011,xx12)	GEMS_GENIE_1	List of dataset names Always sent as “X-Ray Radiation Dose Report”
Document XML	(0041,xx01)	QUASAR_INTERNAL_USE	XML representation of SR document used for DICOM object creation

9.3.3.5 SOP Common Module

TABLE 9-12  
SOP COMMON MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
SOP Class UID	(0008,0016)	1	Uniquely identifies the SOP Class. Always set to 1.2.840.10008.5.1.4.1.1.88.67
SOP Instance UID	(0008,0018)	1	Uniquely identifies the SOP Instance. Internally generated.
Specific Character Set	(0008,0005)	1C	Character Set that expands or replaces the Basic Graphic Set. Defined Terms include for locally created images: ISO_IR 100 = Latin Alphabet No. 1 For images created from MWL, the value is copied from value provided in MWL, if not empty, otherwise ISO_IR 100 is used. Always included into image
Instance Number	(0020,0013)	3	See 9.3.3.1 for more specialization

DIR DOC0660253 REV 15

**9.4 STANDARD EXTENDED AND PRIVATE DATA ATTRIBUTES**

The Product supports the Standard and Private Attributes defined in the following sections in Standard Extended SR SOP Instances as Type 3 data elements.

**9.4.1 Standard Attributes**

The Product supports the following attributes, not specified in the SR IODs, in SOP Instances as Type 3 data elements.

**TABLE 9-13  
STANDARD EXTENDED ATTRIBUTES**

Information Entity Name	Attribute Name	Tag	Use
Study	Allergies	(0010,2110)	Description of prior reaction to contrast agents, or other patient allergies or adverse reactions.
	Pregnancy Status	(0010,21C0)	Describes pregnancy state of patient
	Requested Procedure Comments	(0040,1400)	User-defined Study notes
Series	Performing Physicians' Name	(0008,1050)	Name of the physician(s) administering this Series.
	Protocol Name	(0018,1030)	User-defined description of the conditions under which the Series was performed.
	Request Attributes Sequence	(0040,0275)	Sequence that contains attributes from the Imaging Service Request.
	Performed Procedure Step ID	(0040,0253)	Equipment generated identifier of the protocol carried out within this step.
	Performed Procedure Step Start Date	(0040,0244)	The date that the protocol (SPS) acquisition actually started (doesn't matter if the protocol originated from MWL or was locally scheduled).
	Performed Procedure Step Start Time	(0040,0245)	The time that the protocol (SPS) acquisition actually started (doesn't matter if the protocol originated from MWL or was locally scheduled).
	Performed Procedure Step Description	(0040,0254)	The full path of the performed protocol name.
SR Document	Acquisition Date	(0008,0022)	The date the acquisition of data that resulted in this image started
	Acquisition Time	(0008,0032)	The time the acquisition of data that resulted in this image started
	Counts Accumulated	(0018,0070)	Used as flag to display Dose SR in application Patient List
	Quality Control Image	(0028,0300)	Indicates whether or not this Dose report is created by a quality control procedure or not.

**9.4.2 Private Group QUASAR\_INTERNAL\_USE**

Private Group QUASAR\_INTERNAL\_USE is modeled as part of the Document Information Entity.

**TABLE 9-14  
PRIVATE GROUP QUASAR INTERNAL USE**

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0009,0010)	LO	1	QUASAR_INTERNAL_USE
Origin	(0009,xx12)	LO	1	The origin of the image.
Image Type	(0009,xx1B)	LO	1	Image type string as passed in the scan request
Protocol Scheduled Date	(0009,xx40)	DA	1	Protocol Scheduled Date
Protocol Scheduled Time	(0009,xx41)	TM	1	Protocol Scheduled Time
Matched protocol	(0009,xx43)	LO	1	The originally matched protocol vs. protocol name which is the protocol actually acquired
Private SPS ID	(0009,xx44)	SH	1	Defines protocols that were appended to the original MWL protocol. Contains SPS ID of original MWL Protocol.
Pre-Medication	(0009,xx45)	LO	1	Keeps the Pre-Medication as appears in the To Do list
Private Creator Identification	(0041,0010)	LO	1	QUASAR_INTERNAL_USE
Document XML	(0041,xx01)	UT	1	XML representation of SR document used for DICOM object creation

**9.4.3 Private Group GEMS\_GENIE\_1**

Private Group GEMS\_GENIE\_1 is modeled as part of the Document Information Entity.

**TABLE 9-15  
PRIVATE GROUP GEMS GENIE 1**

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0011,0010)	LO	1	GEMS_GENIE_1
Dataset Name	(0011,xx12)	LO	1-n	List of Dataset names.

**9.5 STANDARD EXTENDED AND PRIVATE CONTEXT GROUPS**

The Product does not support any Standard Extended, Private, and Configurable Context Groups

**9.6 STANDARD, STANDARD EXTENDED AND PRIVATE TEMPLATES**

The Product supports the Standard defined in the following sections. The product does not support any Standard Extended and Private Templates.

**9.6.1 Standard Templates**

The Product supports the following standard templates for SOP Instances created by this product.

**9.6.1.1 Template ID 10011 CT Radiation Dose**

This template defines a container (the root) with subsidiary content items, each of which corresponds to a single CT X-Ray irradiation event entry.

DIR DOC0660253 REV 15

There is a defined recording observer (the system or person responsible for recording the log, generally the system). Accumulated values shall be kept for a single CT scan, which is logical part of study workflow. Multiple CT Radiation Dose objects may be created for one Study.

**TID 10011**  
**CT RADIATION DOSE**

**Type: Extensible      Order: Significant**

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			CONTAINER	EV (113701, DCM, X-Ray Radiation Dose Report)	1	M		Always contains 1 item (continuous type SEPARATE)
2	>	HAS CONCEPT MOD	CODE	EV (121058, DCM, Procedure reported)	1	M		Possible values: EV (P5-08000,SRT, Computed Tomography X-Ray)
3	>>	HAS CONCEPT MOD	CODE	EV (G-C0E8, SRT, Has Intent)	1	M		DCID (3629) Procedure Intent Possible values: EV(R-408C3,SRT,Diagnostic Intent)
4	>		INCLUDE	DTID (1002) Observer Context	1-n	M		See 9.6.1.2
5	>	HAS OBS CONTEXT	DATETIME	EV (113809, DCM, Start of X-Ray Irradiation)	1	M		Start, Date Time of the first CT Irradiation Event of the accumulation
6	>	HAS OBS CONTEXT	DATETIME	EV (113810, DCM, End of X-Ray Irradiation)	1	M		End, Date Time of the last CT Irradiation Event of the accumulation
7	>	HAS OBS CONTEXT	CODE	EV (113705, DCM, Scope of Accumulation)	1	M		DCID (10000) Scope of Accumulation Possible values: EV (113015,DCM,Series)
8	>>	HAS PROPERTIES	UIDREF	EV (112002,DCM,Series Instance UID)	1	M		Series UID of CT Series created during acquisition.
9	>	CONTAINS	INCLUDE	DTID (10012) CT Accumulated Dose Data	1	M		See 9.6.1.4
10	>	CONTAINS	INCLUDE	DTID (10013) CT Irradiation Event Data	1-n	M		See 9.6.1.5 May include 1 or more templates.

DIR DOC0660253 REV 15

12	>	CONTAINS	CODE	EV (113854, DCM, Source of Dose Information)	1-n	M		The primary source of information from which this dose object was constructed. DCID (10021) Source of CT Dose Information Possible values: EV(113856, DCM, Automated Data Collection) EV(113867, DCM, Computed From Image Attributes) . This value is used in the following case: If CT server failed during acquisition, application do not receive specific event with parameters needed to complete Dose SR. In this case automatic recovery process is started, CT reconstruction server creates CT projections from CT raw data stored on disk before crash and CT Dose report parameters are calculated from CT images created and stored in database.
13	>	CONTAINS	INCLUDE	DTID (1020) Person Participant	1	U		Not Used

**9.6.1.2 TID 1002 Observer Context**

Describes the observer (device) that created the Content Items to which this context applies.

**TID 1002**

**OBSERVER CONTEXT**

**Type: Non-Extensible Order: Significant**

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1		HAS OBS CONTEXT	CODE	EV (121005,DCM, Observer Type)	1	M		DCID (270) Observer Type Possible values: (121007,DCM, Device)
2		HAS OBS CONTEXT	INCLUDE	DTID (1003) Person observer identifying attributes	1	MC	IFF Row 1 value = (121006,DCM, Person) or Row 1 is absent	Not Used
3		HAS OBS CONTEXT	INCLUDE	DTID (1004) Device observer identifying attributes	1	M		Always used. See 9.6.1.3

**9.6.1.3 TID 1004 Device Observer Identifying Attributes**

This template (derived from the DICOM General Equipment Module of PS3.3) contains identifying (and optionally descriptive) attributes of devices that are observers.



**TID 1004**  
**DEVICE OBSERVER IDENTIFYING ATTRIBUTES**  
**Type: Extensible                      Order: Significant**

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			UIDREF	EV (121012,DCM, Device Observer UID)	1	M		Implementation UID of Device
2			TEXT	EV (121013,DCM, Device Observer Name)	1	U		Value of Station Name (0008,1010) in General Equipment Module
3			TEXT	EV (121014,DCM, Device Observer Manufacturer)	1	U		Value of Manufacturer (0008,0070) in General Equipment Module
4			TEXT	EV (121015,DCM, Device Observer Model Name)	1	U		Value of Manufacturer's Model Name (0008,1090) in General Equipment Module
5			TEXT	EV (121016,DCM, Device Observer Serial Number)	1	U		Value of Device Serial Number (0018,1000) in General Equipment Module
6			TEXT	EV (121017,DCM, Device Observer Physical Location during observation)	1	U		Not Used
7			CODE	EV (113876, DCM, Device Role in Procedure)	1-n	U		The device performing the observations has other roles, e.g., as the irradiating device in a dose report, this may be recorded here, if not implicit.  BCID (7445) Device Participating Roles (113859,DCM, Irradiating Device); only one Code is defined

**9.6.1.4                      TID 10012                      CT Accumulated Dose Data**

This general template provides detailed information on CT X-Ray dose value accumulations over the scope of accumulation specified for the report (Single CT Scan identified with the single Irradiation Event).

**TID 10012**  
**CT ACCUMULATED DOSE DATA**  
**Type: Extensible                      Order: Significant**

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			CONTAINER	EV (113811, DCM, CT Accumulated Dose Data)	1	M		Contains 2 items
2	>	CONTAINS	NUM	EV (113812, DCM, Total Number of Irradiation Events)	1	M		Units = EV ({events} UCUM, events) May be 1 or more
3	>	CONTAINS	NUM	EV (113813, DCM, CT Dose Length Product Total)	1	M		Units = EV (mGy.cm, UCUM, mGy.cm)

DIR DOC0660253 REV 15

9.6.1.5 TID 10013 CT Irradiation Event Data

This template conveys the dose and equipment parameters of a single irradiation event.

**TID 10013**  
**CT IRRADIATION EVENT DATA**  
**Type: Extensible      Order: Significant**

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			CONTAINER	EV (113819, DCM, CT Acquisition)	1	M		Contains 5 items
2	>	CONTAINS	TEXT	EV (125203, DCM, Acquisition Protocol)	1	U		
3	>	CONTAINS	CODE	EV (123014 , DCM, Target Region)	1	M		DCID (4031) Common Anatomic Regions See 9.6.1.5.1 for possible values
4	>	CONTAINS	CODE	EV (113820, DCM, CT Acquisition Type)	1	M		DCID (10013) CT Acquisition Types Possible values : EV (P5-08001, SRT, Spiral Acquisition) EV (113804, DCM, Sequenced Acquisition) EV (113805, DCM, Constant Angle Acquisition)
6	>	CONTAINS	UIDREF	EV (113769, DCM, Irradiation Event UID)	1	M		Generated Internally.
7	>	CONTAINS	CONTAINER	EV (113822, DCM, CT Acquisition Parameters)	1	M		Contains 7 items
8	>>	CONTAINS	NUM	EV (113824, DCM, Exposure Time)	1	M		Units = EV (s, UCUM, s)
9	>>	CONTAINS	INCLUDE	DTID (10014) Scanning Length	1	M		TID 10014 (included 1 item from Template) See 9.6.1.6
10	>>	CONTAINS	NUM	EV (113826, DCM, Nominal Single Collimation Width)	1	M		Units = EV (mm, UCUM, mm) Possible value: 2.5

11	>>	CONTAINS	NUM	EV (113827, DCM, Nominal Total Collimation Width)	1	M		Units = EV (mm, UCUM, mm) Possible value: 10.0
12	>>	CONTAINS	NUM	EV (113828, DCM, Pitch Factor)	1	MC	IF row 4 equals (P5-08001, SRT, Spiral Acquisition) or equals (113804, DCM, Sequenced Acquisition)	Units = EV ( {ratio}, UCUM, ratio) Possible values: 0.75,1.25, 1.75
13	>>	CONTAINS	NUM	EV (113823, DCM, Number of X-Ray Sources)	1	M		Units = EV ( {X-Ray sources}, UCUM, X-Ray sources) Value = 1 always.
14	>>	CONTAINS	CONTAINER	EV (113831, DCM, CT X-Ray Source Parameters)	1-n	M		Only 1 container is added Contains 5 items
15	>>>	CONTAINS	TEXT	EV (113832, DCM, Identification of the X-Ray Source)	1	M		Always sent as Optima NM/CT 640 Tube
16	>>>	CONTAINS	NUM	EV (113733, DCM, KVP)	1	M		Units = EV (kV, UCUM, kV)
17	>>>	CONTAINS	NUM	EV (113833, DCM, Maximum X-Ray Tube Current)	1	M		Units = EV (mA, UCUM, mA)
18	>>>	CONTAINS	NUM	EV (113734, DCM, X-Ray Tube Current)	1	M		Units = EV (mA, UCUM, mA)
19	>>>	CONTAINS	NUM	EV (113834, DCM, Exposure Time per Rotation)	1	MC	IF row 4 does not equal (113805, DCM, Constant Angle Acquisition)	Units = EV (s, UCUM, s)
21	>	CONTAINS	CONTAINER	EV (113839, DCM, Effective Dose)	1	MC	IF row 4 does not equal (113805, DCM, Constant Angle Acquisition)	Contains 3 items + 2items (containers) from DTID (10015) CT Dose Check Details
22	>>	CONTAINS	NUM	EV (113830, DCM, Mean CTDIvol)	1	M		Units = EV (mGy, UCUM, mGy)

DIR DOC0660253 REV 15

23	>>	CONTAINS	CODE	EV (113835, DCM, CTDIw Phantom Type)	1	M		DCID (4052) Phantom Devices Possible values: EV(113690,DCM, IEC Head Dosimetry Phantom) EV(113691,DCM, IEC Body Dosimetry Phantom)
26	>>	CONTAINS	NUM	EV (113838, DCM, DLP)	1	M		Units = EV (mGy.cm, UCUM, mGy.cm)
30	>>	CONTAINS	INCLUDE	DTID (10015) CT Dose Check Details	1	M		See 9.6.1.7

9.6.1.5.1 CID 4031 Common Anatomic Regions

Context ID 4031  
Common Anatomic Regions

Type: Extensible Version: 20091021

Coding Scheme Designator (0008,0102)	Code Value (0008,0100)	Code Meaning (0008,0104)
SRT	T-D4000	Abdomen
SRT	R-FAB57	Abdomen and Pelvis
SRT	T-15750	Ankle joint
SRT	T-D8200	Arm
SRT	T-04000	Breast
SRT	T-11501	Cervical spine
SRT	T-D3000	Chest
SRT	R-FAB55	Chest and Abdomen
SRT	R-FAB56	Chest, Abdomen and Pelvis
SRT	T-12310	Clavicle
SRT	T-11BF0	Coccyx
SRT	T-D8300	Elbow
SRT	T-D0010	Entire body
SRT	T-D0300	Extremity
SRT	T-D9700	Foot
SRT	T-D8700	Hand
SRT	T-D1100	Head
SRT	T-32000	Heart
SRT	T-15710	Hip joint
SRT	T-D1213	Jaw region
SRT	T-D9200	Knee
SRT	T-D9400	Leg

DIR DOC0660253 REV 15

SRT	T-11502	Lumbar spine
SRT	T-D1600	Neck
SRT	R-FAB52	Neck and Chest
SRT	R-FAB53	Neck, Chest and Abdomen
SRT	R-FAB54	Neck, Chest, Abdomen and Pelvis
SRT	T-D6000	Pelvis
SRT	R-FAB58	Pelvis and lower extremities
SRT	T-D2220	Shoulder
SRT	T-11100	Skull
SRT	T-11AD0	Sacrum
SRT	T-11502	Thoracic spine
DCM	114010	Value Unknown

9.6.1.6 TID 10014 Scanning Length

TID 10014

Scanning Length

Type: Extensible Order: Significant

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Item Type	Condition	Value Set Constraint
1			NUM	EV (113825, DCM, Scanning Length)	1	M			Units = EV (mm, UCUM, mm)

9.6.1.7 TID 10015 CT Dose Check Details

This template records details related to the use of the NEMA Dose Check Standard (NEMA XR-25-2010).

TID 10015

CT DOSE CHECK DETAILS

Type: Extensible Order: Significant

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			CONTAINER	EV (113900, DCM, Dose Check Alert Details)	1	MC	Always sent, Dose Check Alerts are implemented by product	Contains up to 7 Items
2	>	CONTAINS	CODE	EV (113901, DCM, DLP Alert Value Configured)	1	M		DCID (230) Yes-No Possible values: EV(R-0038D, SRT, Yes) EV(R-00339, SRT, No) Value Yes is sent for Dose Reports created by clinical scans, No – for utilities (no patients are scanned)

**NM GENERAL PURPOSE 600/800 SERIES**

**GE Healthcare**

**DICOM CONFORMANCE STATEMENT**

DIR DOC0660253 REV 15

3	>	CONTAINS	CODE	EV (113902, DCM, CTDIvol Alert Value Configured)	1	M		DCID (230) Yes-No EV(R-0038D, SRT, Yes) EV(R-00339, SRT, No) Value Yes is sent for Dose Reports created by clinical scans, No – for utilities (no patients are scanned)
4	>	CONTAINS	NUM	EV (113903, DCM, DLP Alert Value)	1	MC	IFF value of Row 2 is (R-0038D, SRT, Yes)	Units = EV (mGy.cm, UCUM, mGy.cm)
5	>	CONTAINS	NUM	EV (113904, DCM, CTDIvol Alert Value)	1	MC	IFF value of Row 3 is (R-0038D, SRT, Yes)	Units = EV (mGy, UCUM, mGy)
6	>	CONTAINS	NUM	EV(113905, DCM, Accumulated DLP Forward Estimate)	1	MC	IF Accumulated DLP Forward Estimate (Row 6) exceeds DLP Alert Value (Row 4)	Units = EV (mGy.cm, UCUM, mGy.cm)
7	>	CONTAINS	NUM	EV (113906, DCM, Accumulated CTDIvol Forward Estimate)	1	MC	IF Accumulated CTDIvol Forward Estimate (Row 7) exceeds CTDIvol Alert Value (Row 5)	Units = EV (mGy, UCUM, mGy)
9	>	CONTAINS	INCLUDE	DTID (1020) Person Participant	1	MC	IF Accumulated DLP Forward Estimate (Row 6) exceeds DLP Alert Value (Row 4) or Accumulated CTDIvol Forward Estimate (Row 7) exceeds CTDIvol Alert Value (Row 5)	\$PersonProcedureRole = EV (113850, DCM, Irradiation Authorizing) See 9.6.1.7.1
10			CONTAINER	EV (113908, DCM, Dose Check Notification Details)	1	MC	Always sent, Dose Check Alerts are implemented by product	Contains up to 7 Items
11	>	CONTAINS	CODE	EV (113909, DCM, DLP Notification Value Configured)	1	M		DCID (230) Yes-No EV(R-0038D, SRT, Yes) EV(R-00339, SRT, No) Value Yes is sent for Dose Reports created by clinical scans, No – for utilities (no patients are scanned)

DIR DOC0660253 REV 15

12	>	CONTAINS	CODE	EV (113910, DCM, CTDIvol Notification Value Configured)	1	M		DCID (230) Yes-No EV(R-0038D, SRT, Yes) EV(R-00339, SRT, No) Value Yes is sent for Dose Reports created by clinical scans, No – for utilities (no patients are scanned)
13	>	CONTAINS	NUM	EV (113911, DCM, DLP Notification Value)	1	MC	IFF value of Row 11 is (R-0038D, SRT, Yes)	Units = EV (mGy.cm, UCUM, mGy.cm)
14	>	CONTAINS	NUM	EV (113912, DCM, CTDIvol Notification Value)	1	MC	IFF value of Row 12 is (R-0038D, SRT, Yes)	Units = EV (mGy, UCUM, mGy)
15	>	CONTAINS	NUM	EV (113913, DCM, DLP Forward Estimate)	1	MC	IF DLP Forward Estimate (Row 15) exceeds DLP Notification Value (Row 13)	Units = EV (mGy.cm, UCUM, mGy.cm)
16	>	CONTAINS	NUM	EV (113914, DCM, CTDIvol Forward Estimate)	1	MC	IF CTDIvol Forward Estimate (Row 16) exceeds CTDIvol Notification Value (Row 14)	Units = EV (mGy, UCUM, mGy)
18	>	CONTAINS	INCLUDE	DTID (1020) Person Participant	1	UC	IFF DLP Forward Estimate (Row 15) exceeds DLP Notification Value (Row 13) or CTDIvol Forward Estimate (Row 16) exceeds CTDIvol Notification Value (Row 14)	\$PersonProcedureRole = EV (113850, DCM, Irradiation Authorizing) See 9.6.1.7.1

**9.6.1.7.1 TID 1020 Person Participant**

This template describes a person participating in an activity as other than an observer or subject. E.g. for a dose report documenting an irradiating procedure, participants include the person administering the irradiation and the person authorizing the irradiation.

This Template is included with specific contextual parameters from a parent Template.

**TID 1020 Parameters**

Parameter Name	Parameter Usage
\$PersonProcedureRole	Coded term or Context Group for the Concept Name that describes the nature of the person’s participation in this procedure.

**TID 1020**  
**PERSON PARTICIPANT**  
**Type: Extensible Order: Significant**

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			PNAME	EV (113870,DCM, Person Name)	1	M		
2	>	HAS PROPERTIES	CODE	EV (113875,DCM, Person Role in Procedure)	1	M		\$PersonProcedureRole

**9.6.2 Standard Extended Templates**

The Product does not support any extensions or restrictions to standard templates for SOP Instances created by this product.

**9.6.3 Private Templates**

The Product does not support any private templates for SOP Instances created by this product.



## 10. RADIOPHARMACEUTICAL RADIATION DOSE STRUCTURED REPORT INFORMATION OBJECT IMPLEMENTATION

### 10.1 INTRODUCTION

This section specifies the use of the DICOM Radiopharmaceutical Radiation Dose SR IOD to represent results produced by all products except Brivo NM 615 DICOM implementation. Corresponding attributes are conveyed using the module construct.

**Note:** In this section term **NM Camera AE DICOM implementation** refers to **Discovery NM/CT 670\*, NM/CT 800\*, Discovery NM 630 and Optima NM/CT 640 DICOM implementations.**

### 10.2 NM CAMERA MAPPING OF DICOM ENTITIES

The NM Camera AE maps DICOM Information Entities to local Information Entities in the product’s database and user interface.

**TABLE 10-1  
MAPPING OF DICOM ENTITIES TO NM CAMERA ENTITIES**

DICOM IE	NM Camera Entity
Patient	Patient
Study	Study
Series	Series
Document	Report

#### 10.2.1 IOD Module Table

The Radiopharmaceutical Radiation Dose Structured Report Information Object Definitions comprise the modules of the following tables plus Standard Extended and Private Attributes. Standard Extended and Private Attributes are described in Section 10.4

The contents of the SR Document Content are constrained by the supported template, as identified in Section 10.3.3.2.1.1. Standard, Standard Extended and Private Templates are further described in Section 10.6.1.

**TABLE 10-2  
STRUCTURE REPORT IOD MODULES**

Entity Name	Module Name	Usage	Reference
Patient	Patient	Used (same description as for NM IOD)	3.4.1.1
	Clinical Trial Subject	Not Used	N/A
Study	General Study	Used (same description as for NM IOD)	3.4.2.1
	Patient Study	Used (same description as for NM IOD)	3.4.2.2
	Standard Extended Study	Used (same description as for NM IOD)	3.4.2.3

	Private Study	Used (same description as for NM IOD)	3.4.2.4
	Clinical Trial Study	Not Used	NA
Series	SR Document Series	Used	10.3.1.1
	Standard Extended Series	Used	10.3.1.2
	Private Series	Used	10.3.1.3
	Clinical Trial Series	Not Used	N/A
Frame Of Reference	Synchronization	Not Used	N/A
Equipment	General Equipment	Used (same description as for NM IOD)	3.4.5.1
	Enhanced General Equipment	Used	10.3.2.1
Document	SR Document General	Used	10.3.3.1
	SR Document Content	Used	10.3.3.2
	Standard Extended SR Document	Used	10.3.3.3
	Private SR Document	Used	10.3.3.4
	SOP Common	Used	10.3.3.5

**10.3 INFORMATION MODULE DEFINITIONS**

Please refer to DICOM Part 3 (Information Object Definitions) for a description of each of the entities, modules, and attributes contained within the SR Information Objects.

The following modules are included to convey Enumerated Values, Defined Terms, and Optional Attributes supported. Type 1 & Type 2 Attributes are also included for completeness and to define what values they may take and where these values are obtained from when generating the instance. It should be noted that they are the same ones as defined in the DICOM Standard Part 3 (Information Object Definitions). Also note that Attributes which are not presented in tables are not supported.

**10.3.1 Series Entity Modules**

**10.3.1.1 SR Document Series Module**

**TABLE 10-3  
SR DOCUMENT SERIES MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Use
Modality	(0008,0060)	1	Defined Terms used for data created on this system: Value = SR
Series Instance UID	(0020,000E)	1	Internally generated unique identifier of the Series.
Series Number	(0020,0011)	1	A number that identifies this Series. Internally generated using default base value = 600 Calculated as <Base Value>+<N+1>, where <N> is number of RRDR objects already exists in Study

Series Description	(0008,103E)	3	Description of the Series. Always sent. Composed from “ <Administration Event name> followed by suffix “ Dose Report”
Series Date	(0008,0021)	3	Date the Series is created
Series Time	(0008,0031)	3	Time the Series is created
Referenced Performed Procedure Step Sequence	(0008,1111)	2	Uniquely identifies the Administration to which the Series is related. The sequence has exactly 1 item. Sequence is added to all DICOM objects created by system.
>Referenced SOP Class UID	(0008,1150)	1	Set to 1.2.840.10008.3.1.2.3.3
>Referenced SOP Instance UID	(0008,1155)	1	Uniquely identifies the referenced Administration. Internally Generated.

10.3.1.2 Standard Extended Series Module

TABLE 10-4  
STANDARD EXTENDED SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Protocol Name	(0018,1030)	3	User-defined description of the conditions under which the Series was performed. Always sent as “Exposure Monitoring”

10.3.1.3 Private Series Module

TABLE 10-5  
PRIVATE SERIES MODULE ATTRIBUTES

Attribute Name	Tag	Private Creator ID	Attribute Description
Pre-Medication	(0009,xx45)	QUASAR_INTERNAL_USE	Keeps the Pre-Medication as appears in the “To-Do” list Sent as empty string if not received from user input.
Origin	(0009,xx12)	QUASAR_INTERNAL_USE	The origin of the image. Always sent as “exposure monitoring”
Send List Mode flag	(0009, xx24)	QUASAR_INTERNAL_USE	Defined if list mode data shall be sent to list mode data destination station Always sent as 0.
Sequence Type	(0009,xx13)	QUASAR_INTERNAL_USE	Type of sequence created the object Always send as “dose report”

10.3.2 Equipment Entity Modules

10.3.2.1 Enhanced General Equipment Module

TABLE 10-6  
ENHANCED GENERAL EQUIPMENT MODULE ATTRIBUTES

Attribute Name	Tag	Type	Attribute Description
Manufacturer	(0008,0070)	1	Manufacturer of the equipment that produced the composite instances. Default Value GE MEDICAL SYSTEMS
Manufacturer's Model Name	(0008,1090)	1	Manufacturer's model name of the equipment that produced the composite instances. Set to "Tandem_Discovery_670" for Discovery NM/CT 670 implementation Set to "Tandem_Discovery_670_ES" for Discovery NM/CT 670 ES implementation Set to "Tandem_Discovery_670_Pro" for Discovery NM/CT 670 Pro and Discovery 670 DR implementations Set to "Tandem_Apollo" for Discovery NM/CT 670 CZT implementation Set to "Tandem_Discovery_630" for Discovery NM 630 implementation Set to "Tandem_Optima_640" for Optima NM/CT 640 implementation Set to "Tandem_830" for NM 830 implementation Set to "Tandem_850" for NM/CT 850 implementation Set to "Tandem_860" for NM/CT 860 implementation Set to "Tandem_870_DR" for NM/CT 870DR implementation Set to "Tandem_870_CZT" for NM/CT 870CZT implementation Always sent
Device Serial Number	(0018,1000)	1	Manufacturer's serial number of the equipment that produced the composite instances. Always taken from system configuration.
Software Versions	(0018,1020)	1	Manufacturer's designation of software version of the equipment that produced the composite instances Contains Software/Hardware versions of current release e.g. 1.003.608.0\HARDWARE_VERSION_1

10.3.3 Document Entity Modules

10.3.3.1 SR Document General Module

TABLE 10-7  
SR DOCUMENT GENERAL MODULE ATTRIBUTES

Attribute Name	Tag	Type	Use
Instance Number	(0020,0013)	1	A number that identifies this document. Default Value = 1
Completion Flag	(0040,A491)	1	Enumerated Value used: PARTIAL = partial content

Verification Flag	(0040,A493)	1	Enumerated Values: UNVERIFIED = Not attested to.
Completion Flag Description	(0040,A492)	3	Explanation of the value sent in Completion Flag (0040, A491). Contains information about Administration Date, Time and Delivered Dose in Mbq and mci
Content Date	(0008,0023)	1	The date of the document creation.
Content Time	(0008,0033)	1	The time of the document creation.
Verifying Observer Sequence	(0040,A073)	1C	Not sent.
Referenced Request Sequence	(0040,A370)	1C	Not sent
Performed Procedure Code Sequence	(0040,A372)	2	A Sequence that conveys the codes of the performed procedures pertaining to this SOP Instance. Always contains 0 Items.
Current Requested Procedure Evidence Sequence	(0040,A375)	1C	Not Sent
Pertinent Other Evidence Sequence	(0040,A385)	1C	Not Sent
Referenced Instance Sequence	(0008,114A)	1C	Not Sent

10.3.3.2 SR Document Content Module

TABLE 10-8  
SR DOCUMENT CONTENT MODULE ATTRIBUTES

Attribute Name	Tag	Type	Use
Value Type	(0040,A040)	1	The type of the value encoded in this Content Item. Defined Term used: CONTAINER
Concept Name Code Sequence	(0040,A043)	1C	Code describing the concept represented by this Content Item. Only a single Item is included in this sequence.
>Include 'Code Sequence Macro'			EV (113500, DCM, Radiopharmaceutical Radiation Dose Report)
Continuity of Content	(0040,A050)	1	Enumerated Value used: SEPARATE
Content Template Sequence	(0040,A504)	1	Template that describes the content of this Content Item and its subsidiary Content Items. Only a single Item is included in this sequence
>Mapping Resource	(0008,0105)	1	Mapping Resource that defines the template. Defined Terms: DCMR = DICOM Content Mapping Resource
>Template Identifier	(0040, DB00)	1	10021 – Radiopharmaceutical Radiation Dose
Observation DateTime	(0040,A032)	1C	Observation date time is set to Pre-Injection Date Time if defined, otherwise set to Post-Injection Date Time if defined, otherwise set to Administration date time. Always sent.

Content Sequence	(0040,A730)	1	Provides the hierarchical structuring of the Content Tree by recursively nesting Content Items. A parent (or source) Content Item has an explicit relationship to each child (or target) Content Item, conveyed by the Relationship Type (0040,A010) attribute.  One or more Items may be included in this sequence.
> Relationship Type	(0040,A010)	1	Defined Term used: CONTAINS HAS CONCEPT MOD HAS OBS CONTEXT HAS PROPERTIES
> <i>Insert SR Document Content Module</i>			Recursive inclusion to create document content tree. See section 10.3.3.2.1 for the list of supported templates.

**10.3.3.2.1 SR Document Content Descriptions**

**10.3.3.2.1.1 Content Template**

The product supports the following root Templates for Radiopharmaceutical Radiation Dose SR SOP Instances created by the product.

SOP Class	Template ID	Template Name	Use
Radiopharmaceutical Radiation Dose SR	10021	Radiopharmaceutical Radiation Dose	Create

Refer to section 10.6 for a detailed description of the supported templates.

**10.3.3.3 Standard Extended SR Document Module**

**TABLE 10-9  
STANDARD EXTENDED DOSE SR MODULE ATTRIBUTES**

Attribute Name	Tag	Use
Counts Accumulated	(0018,0070)	Always sent as “1”

**10.3.3.4 Private SR Document Module**

**TABLE 10-10  
PRIVATE SR DOCUMENT MODULE ATTRIBUTES**

Attribute Name	Tag	Private Creator ID	Attribute Description
Dataset Name	(0011,xx12)	GEMS_GENIE_1	List of dataset names Always sent as “Radiopharmaceutical Radiation Dose Report”

**10.3.3.5 SOP Common Module**

**TABLE 10-11  
SOP COMMON MODULE ATTRIBUTES**

Attribute Name	Tag	Type	Attribute Description
----------------	-----	------	-----------------------

DIR DOC0660253 REV 15

SOP Class UID	(0008,0016)	1	Uniquely identifies the SOP Class. Always set to 1.2.840.10008.5.1.4.1.1.88.68
SOP Instance UID	(0008,0018)	1	Uniquely identifies the SOP Instance. Internally generated.
Specific Character Set	(0008,0005)	1C	Character Set that expands or replaces the Basic Graphic Set. Defined Terms include for locally created images: ISO_IR 100 = Latin Alphabet No. 1 For images created from MWL, the value is copied from value provided in MWL, if not empty, otherwise ISO_IR 100 is used. Always included into image.
Instance Number	(0020,0013)	3	See 10.3.3.1 for more specialization

**10.4 STANDARD EXTENDED AND PRIVATE DATA ATTRIBUTES**

The Product supports the Standard and Private Attributes defined in the following sections in Standard Extended SR SOP Instances as Type 3 data elements.

**10.4.1 Standard Attributes**

The Product supports the following attributes, not specified in the SR IODs, in SOP Instances as Type 3 data elements.

**TABLE 10-12  
STANDARD EXTENDED ATTRIBUTES**

Information Entity Name	Attribute Name	Tag	Use
Study	Allergies	(0010,2110)	Description of prior reaction to contrast agents, or other patient allergies or adverse reactions.
	Pregnancy Status	(0010,21C0)	Describes pregnancy state of patient
	Requested Procedure Comments	(0040,1400)	User-defined Study notes
Series	Protocol Name	(0018,1030)	Always sent as "Exposure Monitoring"
SR Document	Counts Accumulated	(0018,0070)	Used as flag to display Dose SR in application Patient List

**10.4.2 Private Group QUASAR\_INTERNAL\_USE**

Private Group QUASAR\_INTERNAL\_USE is modeled as part of the Document Information Entity.

**TABLE 10-13  
PRIVATE GROUP QUASAR INTERNAL USE**

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0009,0010)	LO	1	QUASAR_INTERNAL_USE
Origin	(0009,xx12)	LO	1	The origin of the image.
Sequence Type	(0009,xx13)	ST	1	Type of sequence created the object
Send List Mode flag	(0009, xx24)	US	1	Defined if list mode data shall be sent to list mode data destination station
Pre-Medication	(0009,xx45)	LO	1	Keeps the Pre-Medication as appears in the To Do list

DIR DOC0660253 REV 15

**10.4.3 Private Group GEMS\_GENIE\_1**

Private Group GEMS\_GENIE\_1 is modeled as part of the Document Information Entity.

**TABLE 10-14  
PRIVATE GROUP GEMS\_GENIE\_1**

Attribute Name	Tag	VR	VM	Attribute Description and Use
Private Creator Identification	(0011,0010)	LO	1	GEMS_GENIE_1
Dataset Name	(0011,xx12)	LO	1-n	List of Dataset names.

**10.5 STANDARD EXTENDED AND PRIVATE CONTEXT GROUPS**

The Product does not support any Standard Extended, Private, and Configurable Context Groups

**10.6 STANDARD, STANDARD EXTENDED AND PRIVATE TEMPLATES**

The Product supports the Standard defined in the following sections. The product does not support any Standard Extended and Private Templates.

**10.6.1 Standard Templates**

The Product supports the following standard templates for SOP Instances created by this product.

**10.6.1.1 TID 10021 Radiopharmaceutical Radiation Dose**

This template defines a container (the root) with subsidiary content items, each of which corresponds to a single Radiopharmaceutical administration event entry.

Multiple Radiopharmaceutical Radiation Dose objects may be created for one Study.

**TID 10021  
RADIOPHARMACEUTICAL RADIATION DOSE  
Type: Extensible Order: Significant**

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			CONTAINER	EV (113500, DCM, "Radiopharmaceutical Radiation Dose Report")	1	M		Root node
2	>	HAS CONCEPT MOD	CODE	EV (G-C2D0, SRT, "Associated Procedure")	1	M		DCID 3108 "NM/PET Procedures" Possible Values: see 10.6.1.1.1
3	>>	HAS CONCEPT MOD	CODE	EV (G-C0E8, SRT, "Has Intent")	1	M		DCID 3629 "Procedure Intent" Possible Values: EV (R-408C3, SRT, "Diagnostic Intent" )
4	>	CONTAINS	INCLUDE	DTID 10022 "Radiopharmaceutical Administration Event Data"	1	M		See 10.6.1.2
5	>	CONTAINS	INCLUDE	DTID 10024 "Radiopharmaceutical Administration Patient Characteristics"	1	U		See 10.6.1.3



6	>	CONTAINS	TEXT	EV (121106, DCM, "Comment")	1	U		
---	---	----------	------	-----------------------------	---	---	--	--

10.6.1.1.1 CID 3108 “NM/PET Procedures”

NM/PET PROCEDURES

Type: Extensible Version: 20160314

Coding Scheme Designator	Code Value	Code Meaning
SRT	P5-D30F8	Nuclear medicine cardiovascular study
SRT	P5-D6000	Radioisotope study of endocrine system
SRT	P5-D6500	Radioisotope study of hematopoietic system
SRT	P5-D5000	Radioisotope study of gastrointestinal system
SRT	P5-D0063	Radionuclide study for localization of inflammatory disease
SRT	P5-D1000	Radioisotope study of musculoskeletal system
SRT	P5-D90F8	Nuclear medicine diagnostic procedure on nervous system
SRT	P5-D0040	Radionuclide localization of tumor
SRT	P5-D2000	Radioisotope study of respiratory system
SRT	P5-D7000	Radioisotope study of genitourinary system

10.6.1.2 TID 10022 Radiopharmaceutical Administration Event Data

TID 10022

RADIOPHARMACEUTICAL ADMINISTRATION EVENT DATA

Type: Extensible Order: Significant

NL	Rel	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint	
1			CONTAINER	EV (113502, DCM, "Radiopharmaceutical Administration")	1	M		
2	>	CONTAINS	CODE	EV (F-61FDB, SRT, "Radiopharmaceutical agent")	1	M		DCID 25 “Radiopharmaceuticals” Possible values: see 10.6.1.2.1
3	>>	HAS PROPERTIES	CODE	EV (C-10072, SRT, "Radionuclide")	1	M		DCID 18 “Isotopes in Radiopharmaceuticals”
4	>>	HAS PROPERTIES	NUM	EV (R-42806, SRT, "Radionuclide Half Life")	1	M		UNITS = EV (s, UCUM, "seconds")
5	>	CONTAINS	UIDREF	EV (113503, DCM, "Radiopharmaceutical Administration Event UID")	1	M		
6	>	CONTAINS	DATETIME	EV (123003, DCM, "Radiopharmaceutical Start DateTime")	1	M		

DIR DOC0660253 REV 15

7	>	CONTAINS	NUM	EV (113507, DCM, "Administered activity")	1	M		UNITS = EV (MBq, UCUM, "MBq")
9	>	CONTAINS	NUM	EV (113508, DCM, "Pre-Administration Measured Activity")	1	U		UNITS = EV (MBq, UCUM, "MBq")
10	>	CONTAINS	NUM	EV (113509, DCM, "Post-Administration Measured Activity")	1	U		UNITS = EV (MBq, UCUM, "MBq")
11	>	CONTAINS	CODE	EV (G-C340, SRT, "Route of administration")	1	M		BCID 11 "Route of Administration" Possible values: see 10.6.1.2.2
21	>>	HAS PROPERTIES	CODE	EV (G-C581, SRT, "Site of")	1	MC	IF Row 20 equals (G-D101, SRT, "Intravenous route") or (G-D103, SRT, "Intramuscular route");	DCID 3746 "Percutaneous Entry Site" Possible values: see 10.6.1.2.3
22	>>>	HAS CONCEPT MOD	CODE	EV (G-C171, SRT, "Laterality")	1	MC	IF Row 21 has laterality	DCID 244 "Laterality" Possible values: see 10.6.1.2.4
23	>	HAS OBS CONTEXT	INCLUDE	DTID 1020 "Person Participant"	1-n	M		\$PersonProcedureRole = EV (113851, DCM, "Irradiation Administering")

10.6.1.2.1 CID 25 "Radiopharmaceuticals"

**RADIOPHARMACEUTICALS**  
Type: Extensible Version: 20160314

Coding Scheme Designator	Code Value	Code Meaning
SRT	C-B1041	Gallium <sup>67</sup> citrate
SRT	C-145AB	Indium <sup>111</sup> Capromab Pendetide
SRT	C-145AA	Indium <sup>111</sup> Pentetreotide
SRT	C-B1061	Indium <sup>111</sup> pentetate
SRT	C-B1135	Indium <sup>111</sup> oxyquinoline
SRT	C-B110E	Iodine <sup>123</sup> 3-Iodobenzylguanidine MIBG
SRT	C-B1108	Iofetamine I <sup>123</sup> hydrochloride
SRT	C-B1081	Sodium iodide I <sup>123</sup>
SRT	C-B1086	Sodium iodide I <sup>131</sup>
SRT	C-173A5	Krypton <sup>81m</sup>
SRT	C-B1206	Sodium pertechnetate Tc <sup>99m</sup>
SRT	C-B1223	Technetium Tc <sup>99m</sup> exametazine
SRT	C-B1203	Technetium Tc <sup>99m</sup> microaggregated albumin

DIR DOC0660253 REV 15

Coding Scheme Designator	Code Value	Code Meaning
SRT	C-B1225	Technetium Tc <sup>99m</sup> N-substituted iminodiacetate
SRT	C-163B0	Technetium Tc <sup>99m</sup> pentetate
SRT	C-163AB	Technetium Tc <sup>99m</sup> sestamibi
SRT	C-B1215	Technetium Tc <sup>99m</sup> pyro and polyphosphates
SRT	C-B1224	Technetium Tc <sup>99m</sup> tagged red cells
SRT	C-B1222	Technetium Tc <sup>99m</sup> sulfur colloid
SRT	C-163AD	Technetium Tc <sup>99m</sup> Tetrofosmin
SRT	C-163B8	Technetium <sup>99m</sup> Mercaptoacetyl triglycine MAG3
SRT	C-163BD	Technetium <sup>99m</sup> Dimercaptosuccinic Acid DMSA
SRT	C-B1231	Thallos chloride Tl <sup>201</sup>
SRT	C-B1218	Technetium Tc <sup>99m</sup> medronate
GEHC-NUC	10001	Iodine <sup>123</sup> DaTSCAN Ioflupane
GEHC-NUC	10003	Technetium <sup>99m</sup> Ethyl cysteinate dimer
GEHC-NUC	10004	Lutetium <sup>177</sup>
GEHC-NUC	10005	Xeneisol Xe <sup>133</sup>
GEHC-NUC	10006	Technetium <sup>99m</sup> TRODAT-1

10.6.1.2.2 CID 11 “Route of Administration”

**ROUTE OF ADMINISTRATION**  
 Type: Extensible Version: 20160314

Coding Scheme Designator	Code Value	Code Meaning
SRT	G-D101	Intravenous route
SRT	G-D102	Intra-arterial route
SRT	G-D103	Intramuscular route
SRT	G-D104	Subcutaneous route
SRT	G-D106	Intraperitoneal route
SRT	G-D107	Intramedullary route
SRT	G-D108	Intrathecal route
SRT	G-D109	Intra-articular route
NCIt	C38244	Intraepithelial route
SRT	G-D112	Topical route
SRT	G-D140	Oral route
NCIt	C38306	Transluminal route
SRT	G-D144	Intraluminal route
NCIt	C38213	Extraluminal route
SRT	R-40B32	By inhalation
SRT	G-D160	Per rectum
SRT	G-D164	Vaginal route

DIR DOC0660253 REV 15

Coding Scheme Designator	Code Value	Code Meaning
SRT	R-F2C86	Intraventricular route - cardiac
DCM	127070	Retro-orbital route
SRT	G-D172	Nasal route
SRT	G-D17D	Intradermal route
SRT	R-F2CD4	Intratumor route

10.6.1.2.3 CID 3746 “Percutaneous Entry Site”

**PERCUTANEOUS ENTRY SITE**  
 Type: Extensible Version: 20110609

Coding Scheme Designator	Code Value	Code Meaning
SRT	G-D067	Via femoral artery
SRT	G-D1E4	Via radial artery
SRT	G-D05F	Via brachial artery
SRT	G-D054	Via artery
SRT	G-D0C6	Via arm vein
SRT	G-D071	Via femoral vein
SRT	G-D052	Via vein

10.6.1.2.4 DCID 244 “Laterality”

**LATERALITY**  
 Type: Non-Extensible Version: 20030108

Coding Scheme Designator	Code Value	Code Meaning
SRT	G-A100	Right
SRT	G-A101	Left

10.6.1.2.5 TID 1020 Person Participant

This template describes a person participating in an activity as other than an observer or subject. E.g. for a dose report documenting an irradiating procedure, participants include the person administering the irradiation and the person authorizing the irradiation.

This Template is included with specific contextual parameters from a parent Template.

**TID 1020 Parameters**

Parameter Name	Parameter Usage
\$PersonProcedureRole	Coded term or Context Group for the Concept Name that describes the nature of the person’s participation in this procedure.

**TID 1020**  
**PERSON PARTICIPANT**  
 Type: Extensible Order: Significant

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			PNAME	EV (113870,DCM, Person Name)	1	M		
2	>	HAS PROPERTIES	CODE	EV (113875,DCM, Person Role in Procedure)	1	M		\$PersonProcedureRole

10.6.1.3 TID 10024 Radiopharmaceutical Administration Patient Characteristics

TID 10024

RADIOPHARMACEUTICAL ADMINISTRATION PATIENT CHARACTERISTICS

Type: Extensible Order: Significant

	NL	Rel	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			CONTAINER	EV (121118, DCM, "Patient Characteristics")				
2	>	CONTAINS	CODE	EV (121032, DCM, "Subject Sex")	1	U		DCID 7455 "Sex" Possible values: see 10.6.1.3.1
3	>	CONTAINS	NUM	EV (8302-2, LN, "Patient Height")	1	U		UNITS = EV (cm, UCUM, "cm")
4	>	CONTAINS	NUM	EV (29463-7, LN, "Patient Weight")	1	U		UNITS = EV (kg, UCUM, "kg").
5	>	CONTAINS	NUM	EV (F-01860, SRT, "Body Mass Index")	1	U		UNITS = EV (kg/m <sup>2</sup> , UCUM, "kg/m <sup>2</sup> ")
6	>>	INFERRED FROM	CODE	EV (121420, DCM, "Equation")	1	U	Sent only if Row 5 is sent (2 digits)	DT (122265, DCM, "BMI = Wt/Ht <sup>2</sup> ")

10.6.1.3.1 DCID 7455 "Sex"

SEX

Type: Non-Extensible

Version: 20040112

Coding Scheme Designator	Code Value	Code Meaning	Patient's Sex (0010,0040) Equivalent
DCM	M	Male	M
DCM	F	Female	F
DCM	121103	Undetermined sex	O

10.6.2 Standard Extended Templates

The Product does not support any extensions or restrictions to standard templates for SOP Instances created by this product.

DIR DOC0660253 REV 15

**10.6.3 Private Templates**

The Product does not support any private templates for SOP Instances created by this product.